


SIGVARIS
GROUP

Compression Solutions for Veno-Lymphatic Disorders



Inelastic Wraps,
Lymphedema Wear,
Bandages and
Therapy Supplies

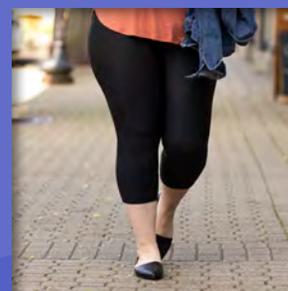

SIGVARIS
GROUP

Established in 1864, SIGVARIS GROUP is, and has always been, a family-owned company. Founded and headquartered in Switzerland, our passion for helping people feel their best is why we've been a global market leader for 158 years and counting.

Our products restore, assist, and enhance people's lives through brands like:

sigvaris

Medical Compression Socks, Hosiery,
and Solutions for Medically Complex Edema

SIGVARIS WELL BEING

Consumer Healthcare: Simplified Sizing
Compression Socks and Hosiery

We are renowned for making exceptional products that combine advanced technology and quality craftsmanship with diverse styles that fit our modern customers' lives.

When heritage, competence, and dedication come together, there's no limit to where we'll go.

Live your life. Be yourself.

For more information, please visit:
sigvaris.com



Live your life. Be yourself.

Medically Complex Edema

The SIGVARIS GROUP line of products for medically complex edema are for the treatment and management of conditions impacting the lymphatic and venous circulatory systems.

The product portfolio contains ready-to-wear garments for the active treatment and maintenance of lymphatic disorders in the upper and lower extremity and for patients managing chronic venous diseases.

Inelastic Wraps: Lower Extremity

Product	Page
Compreflex® Standard Calf & Foot	8
Compreflex Standard Calf	10
Compreflex Transition Calf	12
Compreflex Complete Calf	14
Compreflex Standard Thigh	16
Compreflex Standard Knee	18
Comprefit® Standard Calf & Foot	22
Comprefit Standard Calf	24
Comprefit Plus Calf & Foot	26
Comprefit Standard Thigh	28
Coolflex® Standard Calf	32
Compreboot Lite Foot	34
Compreboot Standard Foot	34
Compreboot Plus Foot	34
Coolflex Standard Foot	34
Compreflex Standard Knee	36
Compreknee Standard Knee (Closed Patella)	36

Lymphedema Wear: Upper Extremity

Product	Page
Neck and Mandible Pad	40
Comprevest	42
Chipsleeve® Standard Arm	44
Medaform® Standard Arm	46
Compreflex Standard Arm	48
Coolflex Standard Arm	50
560 Specialty Secure Armsleeve, Glove, and Gauntlet	52
Medahand®	54
Medaglove®	54
Dorsal Pocket Glove™	54

Lymphedema Wear: Lower Extremity

Product	Page
Compreshort® Short & Capri	58
Compreshape™ Short & Capri	60
Chipsleeve Standard Calf & Foot	64
Chipsleeve Full Leg	66
Medaform Standard Calf & Foot	70
Medaform Full Leg	72
Compreflex Reduce Calf	76
Compreflex Reduce Thigh	78
Compreflex Reduce Knee	80
Compreflex Reduce Foot	82

Table of Contents

Bandages & Therapy Supplies

Product	Page
Bandages, Stockinette, and Tape	86
Foam and Padding	87
Chip Pads	88
Liners	90
Strap Extenders	91

Medical Compression Stockings

Product	Page
230 Essential Cotton	132
860 Essential Opaque for Women	134
860 Essential Opaque for Men	136
550 Specialty Secure	138
602 Specialty Diabetic Compression Socks	140

Custom Garments

Product	Page
Chipsleeve Custom Arm	94
Chipvest	96
Compreshort Custom	98
Chipflow Calf & Foot	102
Chipflow Full Leg	104
Armassist®	106
Legassist® Calf & Foot	113
Legassist Calf & Foot with Shelf	116
Legassist Lobule Compression System	120
Legassist Full Leg	124
Medaboot® Custom Foot	128

Accessories

Product	Page
Gloves	148
Sim-Slide®	148
Magnide On / Off®	148
Cotton Underliner	149
Washing Solution	149
IT STAYS!®	149
Foot Slip	150
Sigvaris DOFF N' DONNER®	150
Sigvaris DOFF N' DONNER CONE	150
MELANY®	150





sigvaris

Inelastic Wraps

LOWER EXTREMITY

Therapeutic Compression for patients with Chronic Venous Disease, Healed or Active Venous Leg Ulcers & Other Complex Edemas

Adjustable Compression Devices are ideal for individuals with chronic venous disease and lymphatic disorders. Designed with inelastic materials, each wrap functions to increase the working pressure applied to the limb while safely maintaining a low resting pressure. With features that aid in application and support independent use, the Sigvaris line of inelastic wraps provides the compression each patient needs throughout their treatment journey.

Compreflex

Ideal for patients in need of an easy to adjust garment that allows for independent use.

Compression Level

20–50mmHg



Colors



Black

Beige



About Compreflex

The Compreflex range of products are low profile compression devices designed to deliver measurable compression to the upper and lower extremities.

Each Compreflex device is made of a soft conforming inelastic material, offering low to minimal stretch for proper fit and gradient application of compression.

The Compreflex garments are designed with key features that support the patient. Each garment includes easy-to-use, patented Accutab® Devices which ensure the accurate application of the prescribed dosage of therapeutic compression in 20–30mmHg, 30–40mmHg, or 40–50mmHg. They also contain an interior stretch panel for easy donning, a contoured calf for limb shaping, and rounded, opposing straps for optimal fit and durability.

Designed for patient independence

Features

Contoured design

- Garment features a curved design to provide a conforming fit

Easy to use

- Straps roll back and attach to garment to prevent tangling

Stretch panel

- The front stretch panel secures garment in place for easy application and high compliance

Ideal for

- Chronic Venous Disease
- Chronic Venous Insufficiency (CVI) including, but not limited to: ulcer due to venous stasis
- Lymphedema - Active and Maintenance
- Non-specific or post surgical edema

Inelastic Wraps

LOWER EXTREMITY

Easily and accurately set compression with Accutab

How to use


- Gently lay straps in place without stretching the fabric
- Attach an Accutab at the end of each strap
- Pull each strap to the prescribed dosage



Compression ranges are indicated by the lines on the Accutab: 20-30mmHg, 30-40mmHg, or 40-50mmHg.

How it works

Proven accuracy for resting pressure

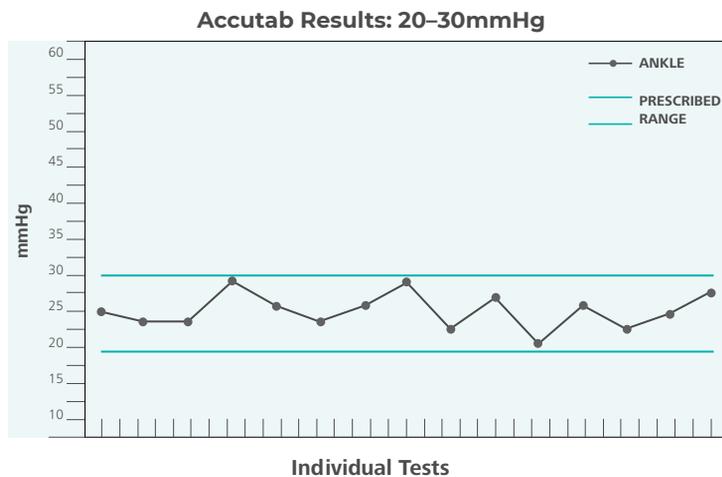


Accutab has been rigorously tested according to established guidelines in the International Compression Club Consensus documents.^{1,2} Testing was conducted using PicoPress, which is recognized for its high accuracy and reproducibility.³

Effective working pressure

- Inelastic compression reduces leg edema and improves venous blood flow^{4,5}
- Inelastic compression can increase arterial perfusion in patients with mild to moderate peripheral arterial disease^{4,5}

Resting Pressure

Donning
Compreflex

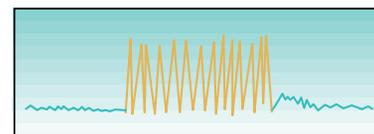

Resting Pressure

Laying Down
Sitting Still

Working pressure

Walking, stair climbing, sit to stand, toe raises or other therapeutic exercises to actuate the calf muscle

Safe, precise resting pressure

- Under the supervision of a healthcare provider, inelastic compression can be used safely for patients with mild to moderate peripheral arterial disease.^{4,5}
- Accutab delivers precise dosages for controlled initial resting pressure when donning a Sigvaris Compreflex garment

Compreflex Standard **CALF & FOOT**

Compression Level

20–50mmHg



Colors



Features

The Compreflex Standard Calf and Foot offers compression in the foot and calf. The inelastic calf wrap is paired with the Compreboot Standard Foot. The Compreboot is low profile and closes completely around the foot. When combined, the products create a complete compression system for the lower leg.

Ideal for

- Chronic Venous Disease
- Chronic Venous Insufficiency (CVI) including, but not limited to: ulcer due to venous stasis
- Lymphedema - Active and Maintenance
- Non-specific or post surgical edema

Includes

- Compreflex Standard Calf
- Compreboot Standard Foot
- 1 Pair of non-compressive Basic Liners
- Accutab Devices

Product Selection Tip

May be right for patients managing chronic venous disease with edema in the foot



Contoured design

Garment features a curved design to provide a conforming fit

Stretch panel

The front stretch panel secures garment in place for easy application and high compliance



Compreflex Standard

CALF & FOOT


Patient Name:	Contact Name:
Account Name:	Account #:
Account Phone #:	P.O. #:
Ship Name & Address:	Date:

Circumference

	Left	Right
C	_____	_____
A	_____	_____
I	_____	_____
J	_____	_____

G Length
 Regular: 30–33.9cm
 Tall: 34–42cm

G	_____
K	_____

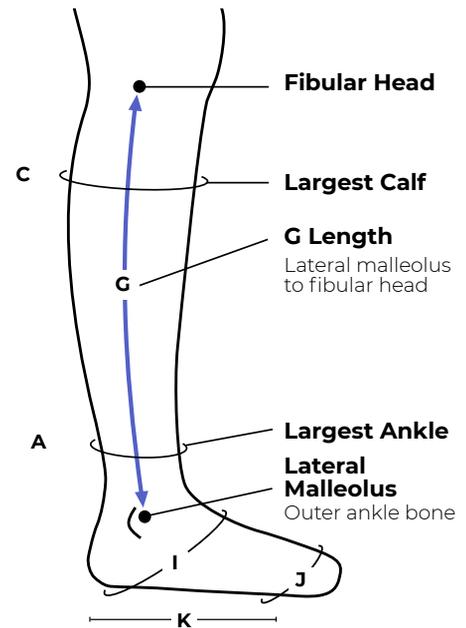

K: Measure **medial** length from heel to 1st metatarsal head

Product Information

Product includes one Compreflex Standard Calf wrap, one Compreboot Standard Foot wrap, and 1 pair of Basic Liners. **Note:** If Compreboot is needed in a different size, select the Compreflex Standard Calf with a separate foot wrap. The Tall Compreflex Standard ships with a regular foot wrap.

<input type="checkbox"/> Right Calf	Size:	Item #:
Color: <input type="checkbox"/> Black <input type="checkbox"/> Beige	Length:	Quantity:
<input type="checkbox"/> Left Calf	Size:	Item #:
Color: <input type="checkbox"/> Black <input type="checkbox"/> Beige	Length:	Quantity:

Companion Products (sold separately)

<input type="checkbox"/> Compreflex Strap Extender	Item #:
Color: <input type="checkbox"/> Black <input type="checkbox"/> Beige	Quantity:

Note: Liner and additional Compreboot options available, each sold separately.

Compreflex Standard Calf Size Chart (A–C)

	Small	Medium	Large	X Large	XX Large
C	29–39cm	34–44cm	40–50cm	46–56cm	54–64cm
A	19–26cm	21–30cm	26–36cm	31–41cm	36–46cm
Black					
Regular	1401-BKR	1402-BKR	1403-BKR	1404-BKR	1405-BKR
Tall	1401-BKT	1402-BKT	1403-BKT	1404-BKT	1405-BKT
Beige					
Regular	1411-BKR	1412-BKR	1413-BKR	1414-BKR	1415-BKR
Tall	1411-BKT	1412-BKT	1413-BKT	1414-BKT	1415-BKT

Compreboot Standard Foot Size Chart (I–K)

	Small		Medium / Large		X Large / XX Large	
	Regular	Long	Regular	Long	Regular	Long
I	28–36cm	28–36cm	≤39cm	≤39cm	≤44cm	≤44cm
J	22–26cm	22–26cm	≤30cm	≤30cm	≤33cm	≤33cm
K	14–18cm	19–23cm	18–20cm	23–25cm	20–22cm	25–27cm
Black	1801-BTR	1801-BTL	1803-BTR	1803-BTL	1805-BTR	1805-BTL
Beige	1811-BTR	1811-BTL	1813-BTR	1813-BTL	1815-BTR	1815-BTL

Compreflex Strap Extender Chart

		1 Strap	10 Unit Pack
Black	Regular	1401-SE1	1401-SEB
Beige	Tall	1411-SE1	1411-SEB



Compreflex Standard CALF

Compression Level

20–50mmHg



Colors



Black

Beige

Features

The Compreflex Standard Calf includes the inelastic calf wrap and a pair of non-compressive Basic Liners. The Compreflex Standard Calf is ideal for patients with chronic edema without foot involvement.

Ideal for

- Chronic Venous Disease
- Chronic Venous Insufficiency (CVI) including, but not limited to: ulcer due to venous stasis
- Lymphedema - Active & Maintenance
- Non-specific or post surgical edema

Includes

- Compreflex Standard Calf
- 1 Pair of non-compressive Basic Liners
- Accutab Devices

Product Selection Tip

May be right for patients managing chronic venous disease with no foot involvement or those needing a different foot piece



Contoured design

Garment features a curved design to provide a conforming fit

Stretch panel

The front stretch panel secures garment in place for easy application and high compliance



Compreflex Standard CALF

Patient Name:	Contact Name:
Account Name:	Account #:
Account Phone #:	P.O. #:
Ship Name & Address:	Date:

Circumference

Left Right


C _____

A _____

G Length

Regular: 30–33.9cm
Tall: 34–42cm

G _____

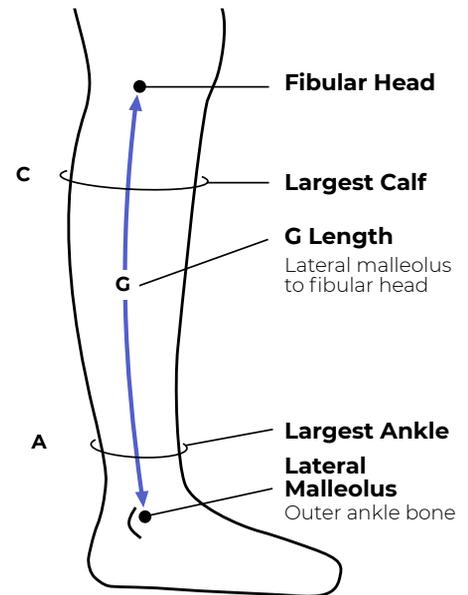

Product Information

Product includes one Compreflex Standard Calf wrap and one pair of non-compressive Basic Liners.

Right Calf Size: Item #:
Color: Black Beige Length: Quantity:

Left Calf Size: Item #:
Color: Black Beige Length: Quantity:

Companion Products (sold separately)

Compreflex Strap Extender Item #:
Color: Black Beige Quantity:

Note: Liner and Compreboot options available, each sold separately.

Compreflex Standard Calf Size Chart (A–C)

	Small	Medium	Large	X Large	XX Large
C	29–39cm	34–44cm	40–50cm	46–56cm	54–64cm
A	19–26cm	21–30cm	26–36cm	31–41cm	36–46cm
Black					
Regular	1401-NF-BKR	1402-NF-BKR	1403-NF-BKR	1404-NF-BKR	1405-NF-BKR
Tall	1401-NF-BKT	1402-NF-BKT	1403-NF-BKT	1404-NF-BKT	1405-NF-BKT
Beige					
Regular	1411-NF-BKR	1412-NF-BKR	1413-NF-BKR	1414-NF-BKR	1415-NF-BKR
Tall	1411-NF-BKT	1412-NF-BKT	1413-NF-BKT	1414-NF-BKT	1415-NF-BKT

Compreflex Strap Extender Chart

		1 Strap	10 Unit Pack
Black	Regular	1401-SE1	1401-SEB
Beige	Tall	1411-SE1	1411-SEB



Compreflex Transition

CALF

Compression Level

20–50mmHg



Colors



Black

Beige

Features

The Compreflex Transition Calf includes the inelastic calf wrap and a pair of Transition Liners. The Transition Liner offers 10–15mmHg of compression in the foot and ankle. The Compreflex Transition Calf is ideal for patients with active venous ulcers and is billable using the A6545-AW code (must be billed while wound is active).

Ideal for

- Chronic Venous Disease
- Chronic Venous Insufficiency (CVI) including, but not limited to: ulcer due to venous stasis
- Lymphedema - Active & Maintenance
- Non-specific or post surgical edema

Includes

- Compreflex Standard Calf
- 1 Pair of Transition Liners
- Accutab Devices

Product Selection Tip

May be right for patients with venous leg ulcers



**Billable under
Medicare code
A6545-AW***

*Code must be billed while wound is active.

Contoured design

Garment features a curved design to provide a conforming fit

Stretch panel

The front stretch panel secures garment in place for easy application and high compliance




Compreflex Transition

CALF

Patient Name:	Contact Name:
Account Name:	Account #:
Account Phone #:	P.O. #:
Ship Name & Address:	Date:

Circumference

	Left	Right
C	_____	_____
A	_____	_____
G Length	Regular: 30–33.9cm Tall: 34–42cm	
G	_____	_____

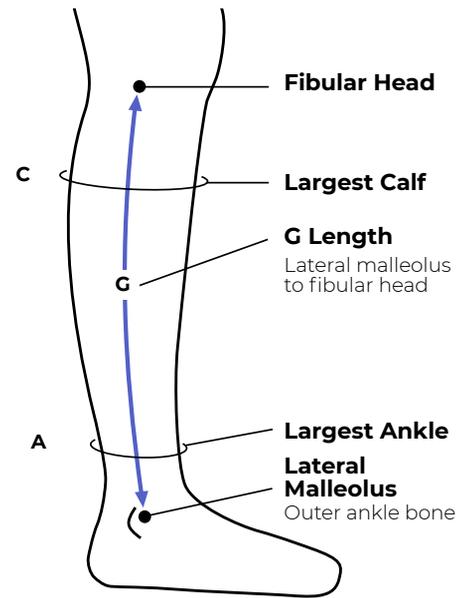

Product Information

Product includes one Compreflex Transition Calf wrap and one pair of Transition Liners. Liner included matches the garment size. Example: The medium garment includes a medium liner.

<input type="checkbox"/> Right Calf	Size:	Item #:
Color: <input type="checkbox"/> Black <input type="checkbox"/> Beige	Length:	Quantity:
<input type="checkbox"/> Left Calf	Size:	Item #:
Color: <input type="checkbox"/> Black <input type="checkbox"/> Beige	Length:	Quantity:

Companion Products (sold separately)

<input type="checkbox"/> Compreflex Strap Extender	Item #:
Color: <input type="checkbox"/> Black <input type="checkbox"/> Beige	Quantity:

Note: Liner and Compreboot options available, each sold separately.

Compreflex Transition Calf Size Chart (A–C)

	Small	Medium	Large	X Large	XX Large
C	29–39cm	34–44cm	40–50cm	46–56cm	54–64cm
A	19–26cm	21–30cm	26–36cm	31–41cm	36–46cm
Black					
Regular	1401-UC-BKR	1402-UC-BKR	1403-UC-BKR	1404-UC-BKR	1405-UC-BKR
Tall	1401-UC-BKT	1402-UC-BKT	1403-UC-BKT	1404-UC-BKT	1405-UC-BKT
Beige					
Regular	1411-UC-BKR	1412-UC-BKR	1413-UC-BKR	1414-UC-BKR	1415-UC-BKR
Tall	1411-UC-BKT	1412-UC-BKT	1413-UC-BKT	1414-UC-BKT	1415-UC-BKT

Compreflex Strap Extender Chart

		1 Strap	10 Unit Pack
Black	Regular	1401-SE1	1401-SEB
Beige	Tall	1411-SE1	1411-SEB



Compreflex Complete CALF

Compression Level

20–50mmHg



Colors



Black

Beige

Features

The Compreflex Complete Calf includes the inelastic calf wrap and a pair of Complete Liners. The Complete Liner offers 10–15mmHg of compression in the foot, ankle, and calf. The Compreflex Complete Calf is ideal for patients with active venous ulcers and is billable using the A6545-AW code (must be billed while wound is active).

Ideal for

- Chronic Venous Disease
- Chronic Venous Insufficiency (CVI) including, but not limited to: ulcer due to venous stasis
- Lymphedema - Active & Maintenance
- Non-specific or post surgical edema

Includes

- Compreflex Standard Calf
- 1 Pair of Complete Liners
- Accutab Devices

Product Selection Tip

May be right for patients with venous leg ulcers under medical provider's orders



Contoured design

Garment features a curved design to provide a conforming fit

Stretch panel

The front stretch panel secures garment in place for easy application and high compliance




Compreflex Complete CALF

Patient Name:	Contact Name:
Account Name:	Account #:
Account Phone #:	P.O. #:
Ship Name & Address:	Date:

Circumference

	Left	Right
C	_____	_____
A	_____	_____
G Length Regular: 30–33.9cm Tall: 34–42cm		
G	_____	_____

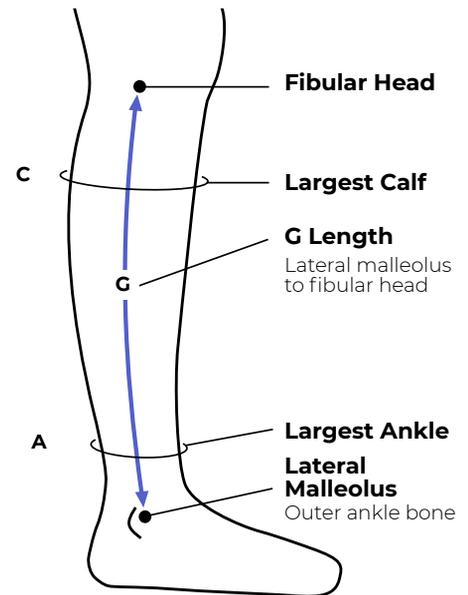

Product Information

Product includes one Compreflex Complete Calf wrap and one pair of Complete Liners. Liner included matches the garment size. Example: The medium garment includes a medium liner.

Right Calf Size: _____ Item #: _____

Color: Black Beige Length: _____ Quantity: _____

Left Calf Size: _____ Item #: _____

Color: Black Beige Length: _____ Quantity: _____

Companion Products (sold separately)

Compreflex Strap Extender Item #: _____

Color: Black Beige Quantity: _____

Note: Liner and Compreboot options available, each sold separately.

Compreflex Complete Calf Size Chart (A-C)

	Small	Medium	Large	X Large	XX Large
C	29–39cm	34–44cm	40–50cm	46–56cm	54–64cm
A	19–26cm	21–30cm	26–36cm	31–41cm	36–46cm
Black					
Regular	1401-FC-BKR	1402-FC-BKR	1403-FC-BKR	1404-FC-BKR	1405-FC-BKR
Tall	1401-FC-BKT	1402-FC-BKT	1403-FC-BKT	1404-FC-BKT	1405-FC-BKT
Beige					
Regular	1411-FC-BKR	1412-FC-BKR	1413-FC-BKR	1414-FC-BKR	1415-FC-BKR
Tall	1411-FC-BKT	1412-FC-BKT	1413-FC-BKT	1414-FC-BKT	1415-FC-BKT

Compreflex Strap Extender Chart

		1 Strap	10 Unit Pack
Black	Regular	1401-SE1	1401-SEB
Beige	Tall	1411-SE1	1411-SEB



Compreflex Standard

THIGH

Compression Level

20–50mmHg



Colors



Black

Beige

Features

The Compreflex Standard Thigh is an inelastic compression device for the thigh. The product includes a Belt and Hip Attachment to support the product throughout daily use.

Pair the Compreflex Standard Thigh with the Compreflex Standard Knee (sold separately) to create a complete compression system for the upper leg.

Ideal for

- Chronic Venous Disease
- Chronic Venous Insufficiency (CVI) including, but not limited to: ulcer due to venous stasis
- Lymphedema - Active & Maintenance
- Non-specific or post surgical edema

Includes

- Compreflex Standard Thigh
- 1 Pair of non-compressive Cotton Liners
- Belt and Hip Attachment
- Accutab Devices

Product Selection Tip

May be right for patients with edema in the upper leg or as a component of a full leg system




Compreflex Standard THIGH

Patient Name:	Contact Name:
Account Name:	Account #:
Account Phone #:	P.O. #:
Ship Name & Address:	Date:

Circumference

Left Right

F	_____	_____
E	_____	_____
D	_____	_____
C1	_____	_____

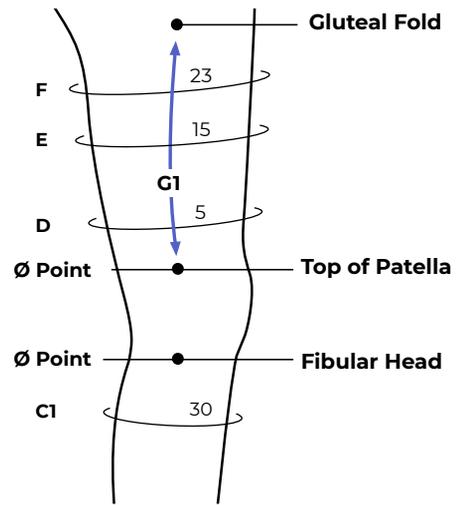

G1 Length

Regular: 20–24.9cm
Tall: 25–30cm

G1	_____	_____
----	-------	-------

Product Information

Product includes one Compreflex Standard Thigh wrap, one pair of Cotton Liners, and one Belt and Hip attachment.

Right Thigh Size: _____ Item #: _____

Color: Black Beige Length: _____ Quantity: _____

Left Thigh Size: _____ Item #: _____

Color: Black Beige Length: _____ Quantity: _____

Companion Products (sold separately)

Compreknee Standard Color: Black Beige

Size: _____ Item #: _____ Quantity: _____

Size: _____ Item #: _____ Quantity: _____

Compreflex Standard Knee Color: Black Beige

Size: _____ Item #: _____ Quantity: _____

Size: _____ Item #: _____ Quantity: _____

Compreflex Standard Thigh Size Chart (D-F)

	Small	Medium	Large	X Large
F	48–58cm	56–66cm	64–74cm	74–84cm
E	43–53cm	51–61cm	58–68cm	68–78cm
D	38–48cm	46–56cm	53–63cm	63–73cm

		Black			
Left	Regular	1401-TCR-L	1402-TCR-L	1403-TCR-L	1404-TCR-L
	Tall	1401-TCT-L	1402-TCT-L	1403-TCT-L	1404-TCT-L
Right	Regular	1401-TCR-R	1402-TCR-R	1403-TCR-R	1404-TCR-R
	Tall	1401-TCT-R	1402-TCT-R	1403-TCT-R	1404-TCT-R
		Beige			
Left	Regular	1411-TCR-L	1412-TCR-L	1413-TCR-L	1414-TCR-L
	Tall	1411-TCT-L	1412-TCT-L	1413-TCT-L	1414-TCT-L
Right	Regular	1411-TCR-R	1412-TCR-R	1413-TCR-R	1414-TCR-R
	Tall	1411-TCT-R	1412-TCT-R	1413-TCT-R	1414-TCT-R

Compreknee Size Chart (C1-D)

	Small	Medium	Large / X Large
D	38–48cm	46–56cm	53–73cm
C1	29–39cm	34–44cm	39–65cm

Standard (Closed Patella)

Color	Item #	Item #	Item #
Black	1101-KP	1102-KP	1103-KP
Beige	1111-KP	1112-KP	1113-KP

Compreflex Standard Knee Size Chart (C1-D)

	Small	Medium	Large	X Large
D	38–48cm	46–56cm	53–63cm	58–73cm
C1	29–39cm	34–44cm	39–49cm	44–55cm
Black	1401-KP	1402-KP	1403-KP	1404-KP
Beige	1411-KP	1412-KP	1413-KP	1414-KP

us_orders@sigvaris.com

P: 800-322-7744 • F: 800-481-5488



For a PDF of the product ordering form, scan the QR code or visit

sigvaris.group/mceforms

Compreflex Standard

KNEE

Compression Level

20–50mmHg



Colors



Black

Beige

Features

The Compreflex Standard Knee is an inelastic compression device for the knee. The product includes an additional strap known as the Butterfly Strap used for added compression to the popliteal or patella regions.

Pair the Compreflex Standard Knee with the Compreflex Standard Thigh and the Compreflex Standard Calf and Foot for total leg coverage.

Ideal for

- Chronic Venous Disease
- Chronic Venous Insufficiency (CVI) including, but not limited to: ulcer due to venous stasis
- Lymphedema - Active & Maintenance
- Non-specific or post surgical edema

Includes

- Compreflex Standard Knee
- Butterfly Strap
- Accutab Devices

Product Selection Tip

May be right for patients needing additional compression around the popliteal or patellar regions



Compreflex Standard

KNEE


Patient Name:	Contact Name:
Account Name:	Account #:
Account Phone #:	P.O. #:
Ship Name & Address:	Date:

Circumference

Left Right

D _____

C1 _____

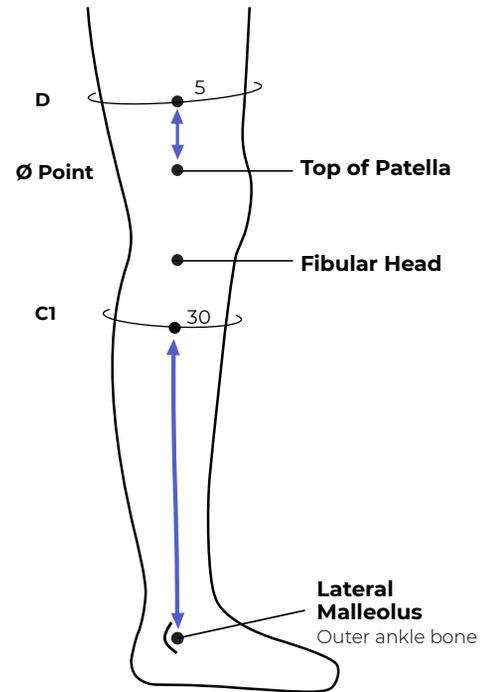

Product Information

<input type="checkbox"/> Compreflex Standard Knee	<input type="checkbox"/> Right	Item #:
Color: <input type="checkbox"/> Black <input type="checkbox"/> Beige	Size:	Quantity:
<input type="checkbox"/> Compreflex Standard Knee	<input type="checkbox"/> Left	Item #:
Color: <input type="checkbox"/> Black <input type="checkbox"/> Beige	Size:	Quantity:

Compreflex Standard Knee Size Chart (C1-D)

	Small	Medium	Large	X Large
D	38-48cm	46-56cm	53-63cm	58-73cm
C1	29-39cm	34-44cm	39-49cm	44-55cm
Black	1401-KP	1402-KP	1403-KP	1404-KP
Beige	1411-KP	1412-KP	1413-KP	1414-KP



Comprefit

Ideal for patients with chronic venous disease or those needing compression therapy following a vascular procedure.

Compression Level

20–50mmHg



Colors



Black



Beige

Strong containment, flexible fit

The Comprefit is a non-stretch, high-containing, inelastic compression garment. Individuals who have undergone a vascular procedure or those that experience extreme fluctuations in limb fluid volume will benefit from the firm containment and the flexibility the slotted closures this garment provides.

Each Comprefit garment features an overlapping band system. This band system mimics the short-stretch bandaging technique; increasing tissue pressure to combat fluctuations in limb fluid volume. The back spine of the garment firmly holds against the back of the limb and works as a stable but comfortable support base.

The Comprefit brand of products includes the Comprefit Plus which pairs the properties of the Comprefit Standard Calf with a foam liner to address the needs of patients with fibrotic or woody limbs.

Each Comprefit garment includes the non-compressive Basic Liner for the protection of the skin.



Ideal for

- Chronic Venous Disease
- Chronic Venous Insufficiency (CVI) including, but not limited to: ulcer due to venous stasis
- Lymphedema
- Non-specific or post surgical edema

Inelastic Wraps

LOWER EXTREMITY



Optimal fit

The slotted band system provides flexibility and conforming fit

Easy to use

Strap ends close to prevent tangling

Inelastic compression

Provides active, gradient compression for the reduction of edema

Overlapping band system

Offers a high level of containment

Comprefit Standard **CALF & FOOT**

Compression Level

20–50mmHg



Colors



Black

Beige

Features

The Comprefit Standard Calf and Foot offers compression in the foot, ankle, and calf. The inelastic calf wrap is paired with the Compreboot and the non-compressive Basic Liner. The Comprefit Standard Calf & Foot is ideal for patients needing a high level of containment.

Ideal for

- Chronic Venous Disease
- Chronic Venous Insufficiency (CVI) including, but not limited to: ulcer due to venous stasis
- Lymphedema
- Non-specific or post surgical edema

Includes

- Comprefit Standard Calf
- Compreboot Standard Foot
- 1 Pair of non-compressive Basic Liners

Product Selection Tip

May be right for patients with extreme fluctuations in limb fluid volume and those requiring foot compression



Comprefit Standard

CALF & FOOT

Patient Name:	Contact Name:
Account Name:	Account #:
Account Phone #:	P.O. #:
Ship Name & Address:	Date:

Product Information

Product includes one Comprefit Standard Calf wrap, one Compreboot Standard Foot wrap, and 1 pair of Basic Liners. The Compreboot Standard Foot wrap included matches the size of the Comprefit Standard Calf wrap selected. **Note:** If Compreboot needed is a different size, select the Comprefit Standard Calf with separate foot piece.

<input type="checkbox"/> Right Calf	Size:	Item #:
Color: <input type="checkbox"/> Black <input type="checkbox"/> Beige	Length:	Quantity:
<input type="checkbox"/> Left Calf	Size:	Item #:
Color: <input type="checkbox"/> Black <input type="checkbox"/> Beige	Length:	Quantity:

Companion Products (sold separately)

<input type="checkbox"/> Comprefit Strap Extender	Item #:
Color: <input type="checkbox"/> Black <input type="checkbox"/> Beige	Quantity:

Note: Liner and additional Compreboot options available, each sold separately.

Compreboot Standard Foot Size Chart (I-K)

	Small		Medium / Large		X Large / XX Large	
	Regular	Long	Regular	Long	Regular	Long
I	28-36cm	28-36cm	≤39cm	≤39cm	≤44cm	≤44cm
J	22-26cm	22-26cm	≤30cm	≤30cm	≤33cm	≤33cm
K	14-18cm	19-23cm	18-20cm	23-25cm	20-22cm	25-27cm
Black	1801-BTR	1801-BTL	1803-BTR	1803-BTL	1805-BTR	1805-BTL
Beige	1811-BTR	1811-BTL	1813-BTR	1813-BTL	1815-BTR	1815-BTL

Comprefit Strap Extender Chart

	1 Strap	10 Straps
Black	1100-SE1	1100-SEB
Beige	1110-SE1	1110-SEB


Circumference

	Left	Right
C1	_____	_____
C	_____	_____
B	_____	_____
A	_____	_____
I	_____	_____
J	_____	_____

G Length

Regular: 30-35.9cm
Tall: 36-42cm

G	_____	_____
K	_____	_____

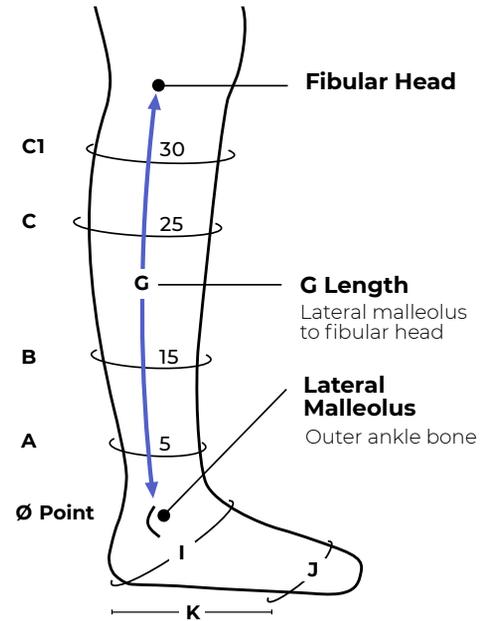

K: Measure **medial** length from heel to 1st metatarsal head

Comprefit Standard Calf Regular Size Chart (A-C1)

Regular		Small	Medium	Large	X Large	XX Large
	C1	29-39cm	34-44cm	39-49cm	48-58cm	55-65cm
C	29-39cm	34-44cm	39-49cm	48-58cm	55-65cm	
B	24-34cm	29-39cm	33-43cm	41-51cm	44-55cm	
A	20-29cm	21-30cm	25-36cm	32-42cm	33-43cm	
Black		1101-BKR	1102-BKR	1103-BKR	1104-BKR	1105-BKR
Beige		1111-BKR	1112-BKR	1113-BKR	1114-BKR	1115-BKR

Comprefit Standard Calf Tall Size Chart (A-C1)

Tall		Small Tall	Medium Tall	Large Tall	X Large Tall	XX Large Tall
	C1	29-39cm	34-44cm	39-49cm	48-58cm	55-65cm
C	26-36cm	31-41cm	35-45cm	44-54cm	50-60cm	
B	21-31cm	25-35cm	30-40cm	36-46cm	40-50cm	
A	20-29cm	21-30cm	25-36cm	32-42cm	33-43cm	
Black		1101-BKT	1102-BKT	1103-BKT	1104-BKT	1105-BKT
Beige		1111-BKT	1112-BKT	1113-BKT	1114-BKT	1115-BKT



Comprefit Standard **CALF**

Compression Level

20–50mmHg



Colors



Features

The Comprefit Standard Calf includes the inelastic calf wrap and a pair of non-compressive Basic Liners. The Comprefit Standard Calf is ideal for patients who require a high level of containment due to fluctuation in limb fluid volume but have no edema in the foot.

Ideal for

- Chronic Venous Disease
- Chronic Venous Insufficiency (CVI) including, but not limited to: ulcer due to venous stasis
- Lymphedema
- Non-specific or post surgical edema

Includes

- Comprefit Standard Calf
- 1 Pair of non-compressive Basic Liners

Product Selection Tip

May be right for patients with extreme fluctuations in limb fluid volume with exclusion of edema in the foot



Comprefit Standard


CALF

Patient Name:	Contact Name:
Account Name:	Account #:
Account Phone #:	P.O. #:
Ship Name & Address:	Date:

Circumference

Left Right

C1	_____	_____
C	_____	_____
B	_____	_____
A	_____	_____
G	_____	_____

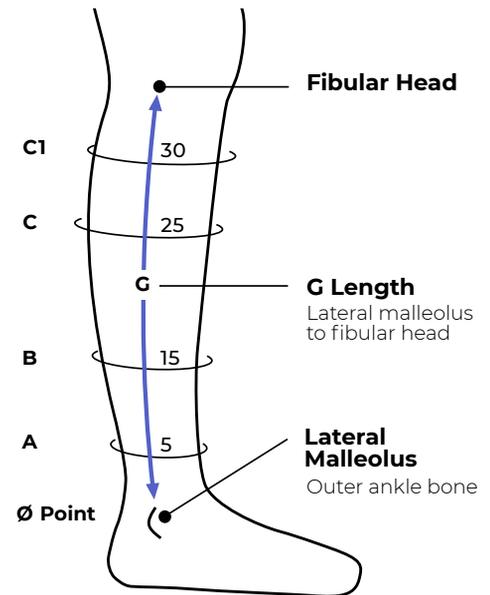

G Length
Regular: 30–35.9cm
Tall: 36–42cm

Product Information

Product includes one Comprefit Standard Calf and one pair of non-compressive Basic Liners.

Right Calf Size: Item #:

Color: Black Beige Length: Quantity:

Left Calf Size: Item #:

Color: Black Beige Length: Quantity:

Companion Products (sold separately)

Comprefit Strap Extender Size: Item #:

Color: Black Beige Quantity:

Note: Liner and Compreboot options available, each sold separately.

Comprefit Strap Extender Chart

	1 Strap	10 Straps
Black	1100-SE1	1100-SEB
Beige	1110-SE1	1110-SEB

Comprefit Standard Calf Regular Size Chart (A–C1)

	Small	Medium	Large	X Large	XX Large	
Regular	C1	29–39cm	34–44cm	39–49cm	48–58cm	55–65cm
	C	29–39cm	34–44cm	39–49cm	48–58cm	55–65cm
	B	24–34cm	29–39cm	33–43cm	41–51cm	44–55cm
	A	20–29cm	21–30cm	25–36cm	32–42cm	33–43cm
	Black	1101-NF-BKR	1102-NF-BKR	1103-NF-BKR	1104-NF-BKR	1105-NF-BKR
Beige	1111-NF-BKR	1112-NF-BKR	1113-NF-BKR	1114-NF-BKR	1115-NF-BKR	

Comprefit Standard Calf Tall Size Chart (A–C1)

	Small Tall	Medium Tall	Large Tall	X Large Tall	XX Large Tall	
Tall	C1	29–39cm	34–44cm	39–49cm	48–58cm	55–65cm
	C	26–36cm	31–41cm	35–45cm	44–54cm	50–60cm
	B	21–31cm	25–35cm	30–40cm	36–46cm	40–50cm
	A	20–29cm	21–30cm	25–36cm	32–42cm	33–43cm
	Black	1101-NF-BKT	1102-NF-BKT	1103-NF-BKT	1104-NF-BKT	1105-NF-BKT
Beige	1111-NF-BKT	1112-NF-BKT	1113-NF-BKT	1114-NF-BKT	1115-NF-BKT	



Comprefit Plus

CALF & FOOT

Compression Level

20–50mmHg



Color



Black

Features

The Comprefit Plus Calf & Foot pairs the high containing properties of the Comprefit Standard Calf with the stimulating properties of a foam liner. Combined, the calf and foot compression system works to breakdown fibrotic tissue while actively directing lymphatic fluid to the paths optimal for drainage.

Ideal for

- Chronic Venous Disease
- Chronic Venous Insufficiency (CVI) including, but not limited to: ulcer due to venous stasis
- Lymphedema
- Non-specific or post surgical edema

Includes

- Comprefit Standard Calf
- Compreboot Standard Foot
- Open-Cell Flat Foam Liner
- 1 Pair of non-compressive Basic Liners

Product Selection Tip

May be right for patients with fibrotic limbs and those that experience fluctuations in limb fluid volume



Comprefit Plus

CALF & FOOT

Patient Name:	Contact Name:
Account Name:	Account #:
Account Phone #:	P.O. #:
Ship Name & Address:	Date:

Product Information

Product includes one Comprefit Standard Calf wrap, one flat foam liner, one Compreboot Standard Foot wrap, and 1 pair of Basic Liners. The Compreboot Standard Foot wrap included matches in size of the Comprefit Standard Calf wrap selected.

<input type="checkbox"/> Right Calf	Size:	Item #:
Color: Black	Length:	Quantity:
<input type="checkbox"/> Left Calf	Size:	Item #:
Color: Black	Length:	Quantity:

Companion Products (sold separately)

<input type="checkbox"/> Comprefit Strap Extender	Item #:
Color: <input type="checkbox"/> Black	Quantity:

Compreboot Standard Foot Size Chart (I-K)

	Small		Medium / Large		X Large / XX Large	
	Regular	Long	Regular	Long	Regular	Long
I	28-36cm	28-36cm	≤39cm	≤39cm	≤44cm	≤44cm
J	22-26cm	22-26cm	≤30cm	≤30cm	≤33cm	≤33cm
K	14-18cm	19-23cm	18-20cm	23-25cm	20-22cm	25-27cm
Black	1801-BTR	1801-BTL	1803-BTR	1803-BTL	1805-BTR	1805-BTL

Comprefit Strap Extender Chart


	1 Strap	10 Straps
Black	1100-SE1	1100-SEB

Circumference

	Left	Right
C1	_____	_____
C	_____	_____
B	_____	_____
A	_____	_____
I	_____	_____
J	_____	_____

G Length
Regular: 30-35.9cm
Tall: 36-42cm

G	_____
K	_____

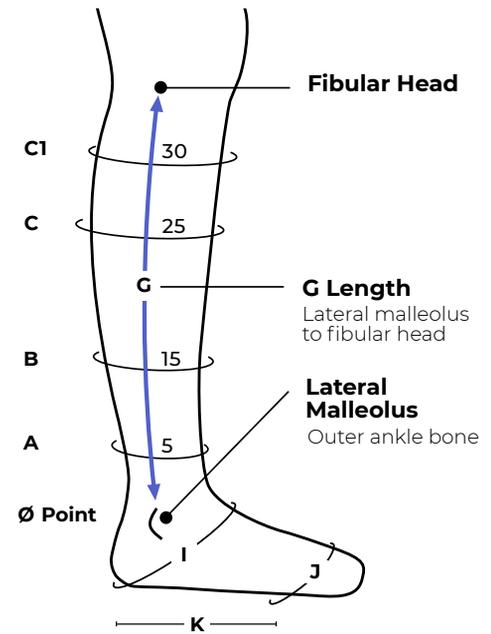

K: Measure **medial** length from heel to 1st metatarsal head

Comprefit Plus Calf Regular Size Chart (A-C1)

		Small	Medium	Large	X Large	XX Large
		Regular	C1	24-34cm	29-39cm	34-44cm
	C	24-34cm	29-39cm	34-44cm	43-53cm	50-60cm
	B	19-29cm	24-34cm	28-38cm	36-46cm	39-50cm
	A	15-24cm	16-25cm	20-31cm	27-37cm	28-38cm
Black		1101-BPR	1102-BPR	1103-BPR	1104-BPR	1105-BPR

Comprefit Plus Calf Tall Size Chart (A-C1)

		Small Tall	Medium Tall	Large Tall	X Large Tall	XX Large Tall
		Tall	C1	24-34cm	29-39cm	34-44cm
	C	21-31cm	26-36cm	30-40cm	39-49cm	45-55cm
	B	16-26cm	20-30cm	25-35cm	31-41cm	35-45cm
	A	15-24cm	16-25cm	20-31cm	27-37cm	28-38cm
Black		1101-BPT	1102-BPT	1103-BPT	1104-BPT	1105-BPT



Comprefit Standard **THIGH**

Compression Level

20–50mmHg



Colors



Black

Beige

Features

Comprefit Standard Thigh provides compression to the thigh. The product includes a Belt and Hip Attachment to secure the product in place for all-day wear.

Pair the Comprefit Standard Thigh with a knee garment (sold separately) for complete upper leg compression.

Ideal for

- Chronic Venous Disease
- Chronic Venous Insufficiency (CVI) including, but not limited to: ulcer due to venous stasis
- Lymphedema
- Non-specific or post surgical edema

Includes

- Comprefit Standard Thigh
- Belt and Hip Attachment
- 1 Pair of non-compressive Cotton Liners

Product Selection Tip

May be right for patients with chronic edema in the upper leg




Comprefit Standard THIGH

Patient Name:	Contact Name:
Account Name:	Account #:
Account Phone #:	P.O. #:
Ship Name & Address:	Date:

Circumference

Left Right

F	_____	_____
E	_____	_____
D	_____	_____
Ø Point	_____	_____
Ø Point	_____	_____
C1	_____	_____

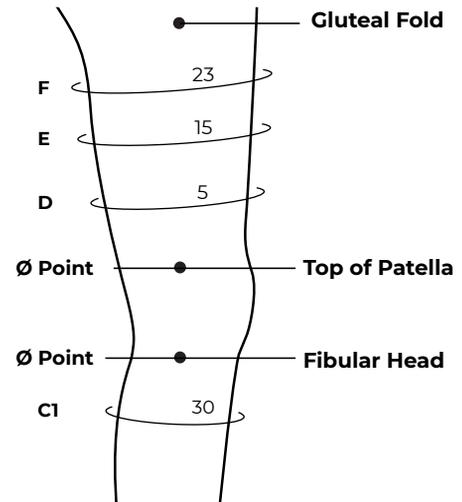

Product Information

Product includes one Comprefit Standard Thigh wrap, one pair of Cotton Liners, and one Belt and Hip attachment.

Right Thigh Size: Item #:

Color: Black Beige Quantity:

Left Thigh Size: Item #:

Color: Black Beige Quantity:

Companion Products (sold separately)

Compreknee Standard Color: Black Beige

Size: Item #: Quantity:

Size: Item #: Quantity:

Compreflex Standard Knee Color: Black Beige

Size: Item #: Quantity:

Size: Item #: Quantity:

Comprefit Standard Thigh Size Chart (D-F)

	Small	Medium	Large	X Large
F	48-58cm	56-66cm	64-74cm	74-84cm
E	43-53cm	51-61cm	58-68cm	68-78cm
D	38-48cm	46-56cm	53-63cm	63-73cm
Black	1101-TC	1102-TC	1103-TC	1104-TC
Beige	1111-TC	1112-TC	1113-TC	1114-TC

Compreflex Standard Knee Size Chart (C1-D)

	Small	Medium	Large	X Large
D	38-48cm	46-56cm	53-63cm	58-73cm
C1	29-39cm	34-44cm	39-49cm	44-55cm
Black	1401-KP	1402-KP	1403-KP	1404-KP
Beige	1411-KP	1412-KP	1413-KP	1414-KP

Compreknee Size Chart (C1-D)

	Small	Medium	Large / X Large
D	38-48cm	46-56cm	53-73cm
C1	29-39cm	34-44cm	39-65cm
Standard (Closed Patella)			
Black	1101-KP	1102-KP	1103-KP
Beige	1111-KP	1112-KP	1113-KP



Coolflex

Ideal for patients needing a lightweight compression garment with fragile or broken skin.

Compression Level

20–50mmHg



Color



Black


Lightweight, breathable compression wear

Coolflex is a lightweight compression device, optimal for individuals with fragile or broken skin. Made of a non-stretch mesh material, the Coolflex transfers heat and wicks away moisture, leaving the skin cool and dry.

Ideal for both venous and lymphatic limb disorders, the garments are available in standard and max sizing and can accommodate limbs that feature disproportionately small ankles with a larger calf.

Each Coolflex calf garment includes the unique, EZ-Close™ pressure system. The EZ-Close pressure system is positioned laterally, allowing the patient to apply the product without bending forward. The bands of the closure system are made to an optimal length, making them easy to grip and close using minimal hand strength and little dexterity.

The features of the closure system, along with the built-in guidelines for application, promote independence and encourage improved patient outcomes.

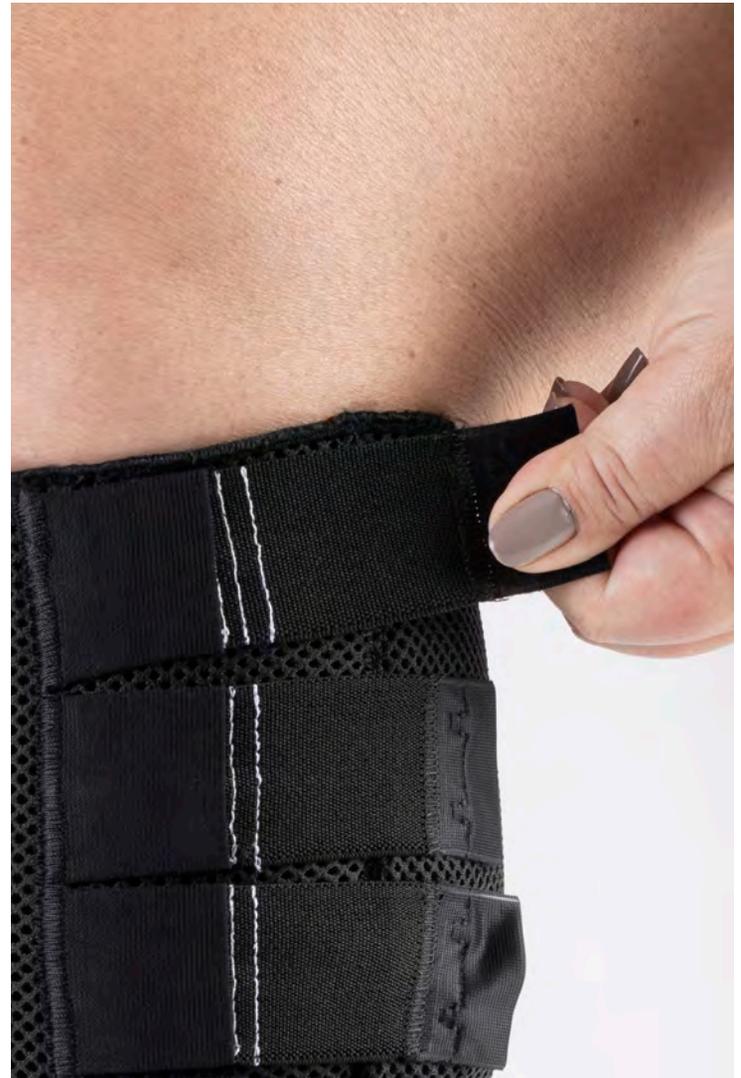

Ideal for

- Chronic Venous Disease
- Chronic Venous Insufficiency (CVI) including, but not limited to: ulcer due to venous stasis
- Fragile, sensitive skin
- Lymphedema
- Lipedema

This product is patent pending.

Inelastic Wraps

LOWER EXTREMITY



Most Breathable Garment

The Spacer Fabric is designed to be highly breathable for cool-wearing comfort

One-Handed Donning

EZ-Close Pressure System of eight tabs makes one-handed, lateral donning simple, even for patients with dexterity issues

Innovative Design

Durable, metal hook closures for maximum strength and durability



Accurate compression is easy with the Coolflex

The integrated EZ-Close Pressure System has visual indicators for the accurate delivery of 20–30mmHg, 30–40mmHg, or 40–50mmHg on each tab

us_orders@sigvaris.com

P: 800-322-7744 • F: 800-481-5488

Coolflex Standard **CALF**

Compression Level

20–50mmHg



Color



Black

Features

The Coolflex Standard Calf includes the inelastic calf wrap and a pair of non-compressive Basic Liners.

Ideal for

- Chronic Venous Disease
- Chronic Venous Insufficiency (CVI) including, but not limited to: ulcer due to venous stasis
- Fragile, sensitive skin
- Lymphedema
- Lipedema

Includes

- Coolflex Standard Calf
- 1 Pair of non-compressive Basic Liners

Product Selection Tip

May be right for patients with active venous leg ulcers and recurring ulcers or those with a small ankle and a large calf



Coolflex Standard

CALF


Patient Name:	Contact Name:
Account Name:	Account #:
Account Phone #:	P.O. #:
Ship Name & Address:	Date:

Circumference

	Left	Right
C1	_____	_____
C	_____	_____
B	_____	_____
A	_____	_____

G Length
 Regular: 30-35.9cm
 Tall: 36-42cm

G	_____
---	-------

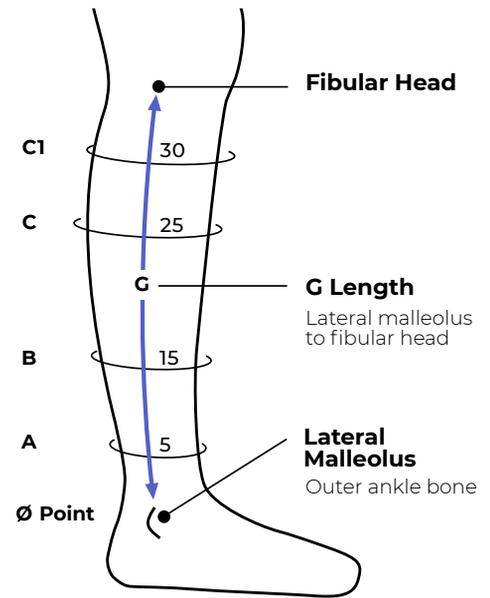

Product Information

Product includes one Coolflex Standard Calf and one pair of non-compressive Basic Liners.

<input type="checkbox"/> Right Calf	Size:	Item #:
Color: Black	Length:	Quantity:
<input type="checkbox"/> Left Calf	Size:	Item #:
Color: Black	Length:	Quantity:

Note: Liner and foot compression options available, each sold separately.

Coolflex Standard Calf Size Chart (A-C1) (Available in Black Only)

	Small	Small Max	Medium / Large	Medium / Large Max	X Large / XX Large
C1	28-39cm	35-50cm	34-50cm	45-64cm	45-64cm
C	28-39cm	35-50cm	34-50cm	45-64cm	45-64cm
B	23-33cm	27-38cm	29-43cm	35-49cm	37-54cm
A	16-26cm	16-26cm	23-33cm	23-33cm	30-42cm

		Black				
Left	Regular	2101-NF-BKR-L	2101-NF-BMR-L	2103-NF-BKR-L	2103-NF-BMR-L	2105-NF-BKR-L
	Tall	2101-NF-BKT-L	2101-NF-BMT-L	2103-NF-BKT-L	2103-NF-BMT-L	2105-NF-BKT-L
Right	Regular	2101-NF-BKR-R	2101-NF-BMR-R	2103-NF-BKR-R	2103-NF-BMR-R	2105-NF-BKR-R
	Tall	2101-NF-BKT-R	2101-NF-BMT-R	2103-NF-BKT-R	2103-NF-BMT-R	2105-NF-BKT-R



Compreboot Lite Foot

Variable compression



Features

- Ultra-thin profile that easily fits into most footwear
- Highly adjustable design may be trimmed to fit
- Provides mild therapeutic compression

Color



Black



Compreboot Standard Foot

Variable compression



Features

- Four-way stretch for firm compression around ankle
- Ideal for mild to moderate chronic edema
- Low stretch material conforms to foot shape

Colors



Black



Beige



Compreboot Plus Foot

Variable compression



Features

- Chip foam lining creates high / low pressure ridges
- Removable gripper patches provide traction
- Provides moderate therapeutic compression

Color



Black



Coolflex Standard Foot

Variable compression



Features

- Low profile design allows boot to fit in most footwear*
- Spacer Fabric is highly breathable for cool-wearing comfort
- Durable metal clasps for closure

Color



Black



Foot compression garments are pictured with the Basic Liner. Basic Liner is not included.

* Data on file

Foot Compression Wear

Patient Name:	Contact Name:
Account Name:	Account #:
Account Phone #:	P.O. #:
Ship Name & Address:	Date:

Circumference


Left Right

I _____

J _____

Length

K _____

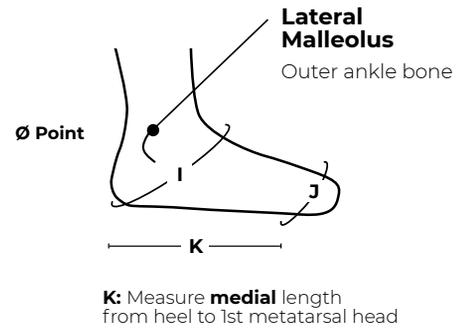

Product Information

Compreboot Lite Foot Right
Color: Black Size: Item #: Quantity:

Compreboot Lite Foot Left
Color: Black Size: Item #: Quantity:

Compreboot Standard Foot Right
Color: Black Beige Size: Item #: Quantity:

Compreboot Standard Foot Left
Color: Black Beige Size: Item #: Quantity:

Compreboot Plus Foot Right
Color: Black Size: Item #: Quantity:

Compreboot Plus Foot Left
Color: Black Size: Item #: Quantity:

Coolflex Standard Foot Right
Color: Black Size: Item #: Quantity:

Coolflex Standard Foot Left
Color: Black Size: Item #: Quantity:

Measuring Instructions

Foot Size:

- I: Measure circumference of ankle bend and heel
- J: Measure circumference across metatarsal heads
- K: Measure length from heel to 1st metatarsal head

Compreboot Lite Foot Size Chart (I-K) (Available in Black Only)

	Small / Medium / Large		X Large / XX Large	
	Regular	Long	Regular	Long
I	28-36cm	28-36cm	≤44cm	≤44cm
J	22-26cm	22-26cm	≤33cm	≤33cm
K	14-18cm	19-23cm	20-22cm	25-27cm
Black	1803-BLR	1803-BLL	1805-BLR	1805-BLL

Compreboot Standard Foot Size Chart (I-K)

	Small		Medium / Large		X Large / XX Large	
	Regular	Long	Regular	Long	Regular	Long
I	28-36cm	28-36cm	≤39cm	≤39cm	≤44cm	≤44cm
J	22-26cm	22-26cm	≤30cm	≤30cm	≤33cm	≤33cm
K	14-18cm	19-23cm	18-20cm	23-25cm	20-22cm	25-27cm
Black	1801-BTR	1801-BTL	1803-BTR	1803-BTL	1805-BTR	1805-BTL
Beige	1811-BTR	1811-BTL	1813-BTR	1813-BTL	1815-BTR	1815-BTL

Coolflex Standard Foot Size Chart (I-K) (Available in Black Only)

	Small		Medium / Large		X Large / XX Large	
	Regular	Long	Regular	Long	Regular	Long
I	28-38cm	28-38cm	30-43cm	30-43cm	34-49cm	34-49cm
J	22-26cm	22-26cm	25-30cm	25-30cm	28-33cm	28-33cm
K	14-18cm	19-21cm	18-20cm	21-23cm	20-22cm	23-25cm

Left

Black	2101-BTR-L	2101-BTL-L	2103-BTR-L	2103-BTL-L	2105-BTR-L	2105-BTL-L
--------------	------------	------------	------------	------------	------------	------------

Right

Black	2101-BTR-R	2101-BTL-R	2103-BTR-R	2103-BTL-R	2105-BTR-R	2105-BTL-R
--------------	------------	------------	------------	------------	------------	------------

Compreboot Plus Foot Size Chart (I-K) (Available in Black Only)

	Small		Medium / Large		X Large / XX Large	
	Regular	Long	Regular	Long	Regular	Long
I	28-36cm	28-34cm	≤39cm	≤39cm	≤44cm	≤44cm
J	22-26cm	22-25cm	≤30cm	≤28cm	≤33cm	≤33cm
K	14-18cm	19-23cm	18-20cm	23-25cm	20-22cm	25-27cm
Black	1801-BPR	1801-BPL	1803-BPR	1803-BPL	1805-BPR	1805-BPL

us_orders@sigvaris.com

P: 800-322-7744 • F: 800-481-5488



For a PDF of the product ordering form, scan the QR code or visit

sigvaris.group/mceforms

Compreflex Standard Knee

20-50mmHg



Features

- Provides accurate, inelastic compression levels of 20-50mmHg
- Includes inner donning sleeve for easy application
- Product includes a Butterfly Strap for difficult to compress areas

Colors



Black



Beige



Compreknee Standard Knee

Variable compression



Features

- Garment covers the knee, providing full coverage
- Straps provide a full-range of adjustment to create a comfortable fit
- Conforming material for a firm secure fit

Colors



Black



Beige



Knee Compression Wear


Patient Name:	Contact Name:
Account Name:	Account #:
Account Phone #:	P.O. #:
Ship Name & Address:	Date:

Circumference

Left Right

D _____

C1 _____

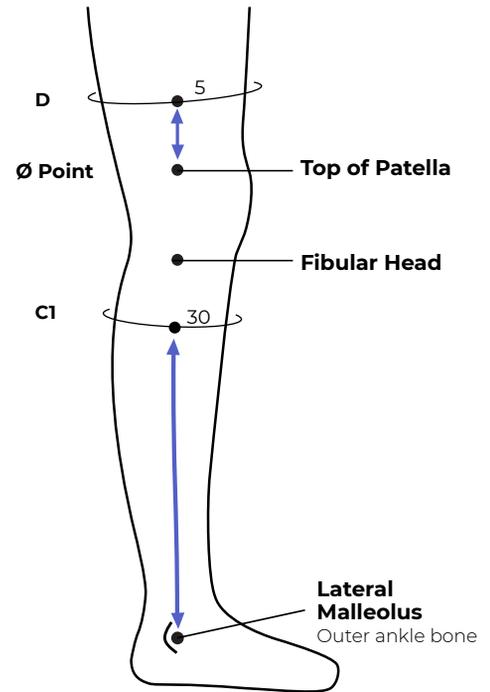

Product Information

Compreknee Standard Right Item #:
 Color: Black Beige Size: Quantity:

Compreknee Standard Left Item #:
 Color: Black Beige Size: Quantity:

Compreflex Standard Knee Right Item #:
 Color: Black Beige Size: Quantity:

Compreflex Standard Knee Left Item #:
 Color: Black Beige Size: Quantity:

Compreflex Standard Knee Size Chart (C1-D)


	Small	Medium	Large	X Large
D	38-48cm	46-56cm	53-63cm	58-73cm
C1	29-39cm	34-44cm	39-49cm	44-55cm
Black	1401-KP	1402-KP	1403-KP	1404-KP
Beige	1411-KP	1412-KP	1413-KP	1414-KP

Compreknee Size Chart (C1-D)

	Small	Medium	Large / X Large
D	38-48cm	46-56cm	53-73cm
C1	29-39cm	34-44cm	39-65cm
Standard (Closed Patella)			
Black	1101-KP	1102-KP	1103-KP
Beige	1111-KP	1112-KP	1113-KP





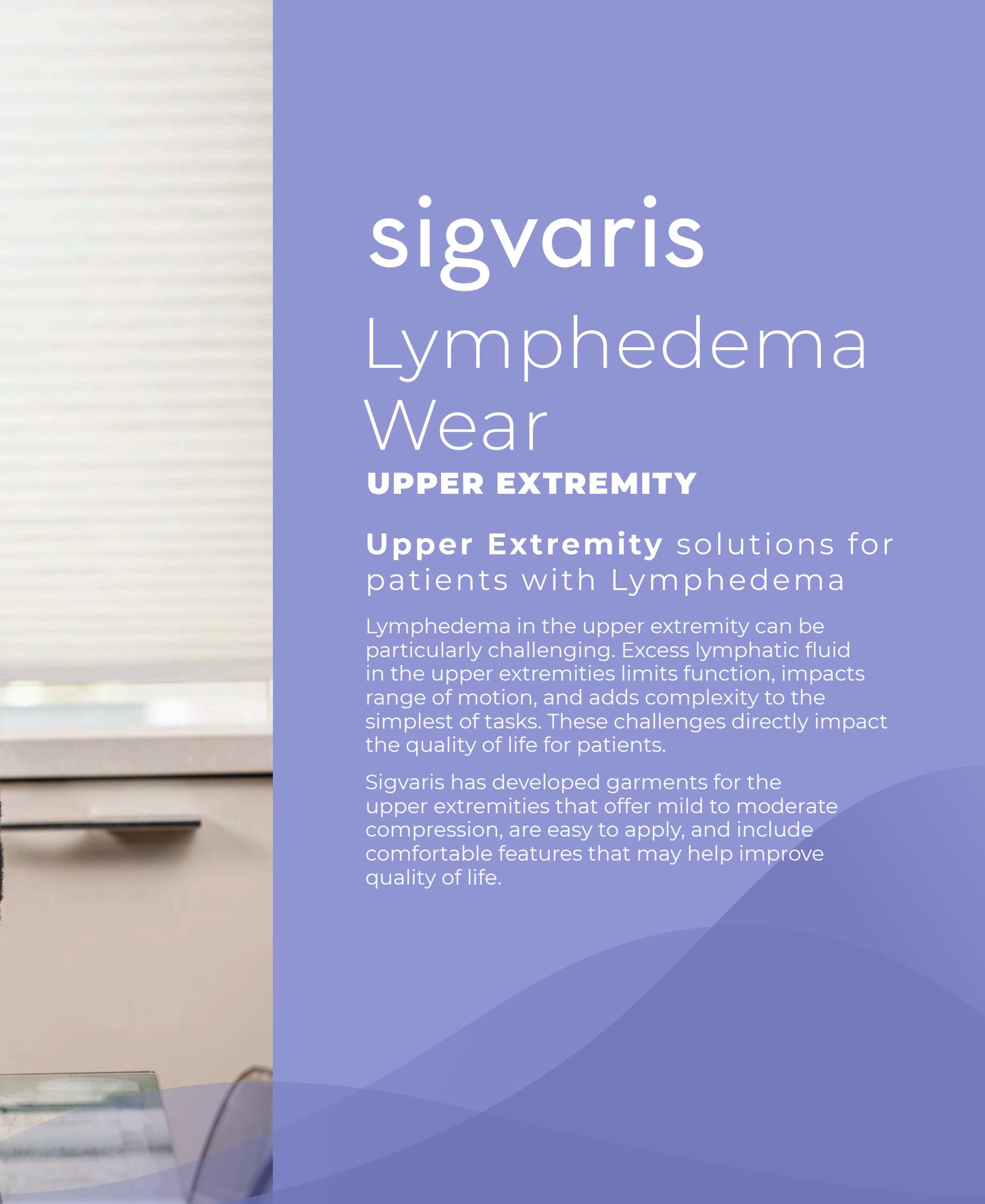

sigvaris

Lymphedema Wear

UPPER EXTREMITY

Upper Extremity solutions for patients with Lymphedema

Lymphedema in the upper extremity can be particularly challenging. Excess lymphatic fluid in the upper extremities limits function, impacts range of motion, and adds complexity to the simplest of tasks. These challenges directly impact the quality of life for patients.

Sigvaris has developed garments for the upper extremities that offer mild to moderate compression, are easy to apply, and include comfortable features that may help improve quality of life.

Neck and Mandible PAD

Compression Level

8-10mmHg



Color



Black

Features

The Neck and Mandible Pad is made from soft materials that conform to the mandible and neck. The Neck and Mandible Pad provides mild compression and uses easily adjustable straps that secure the garment in place.

Ideal for

- Lymphedema
- Post-surgical or non-specific edema

Includes

- Neck and Mandible Pad

Product Selection Tip

May be right for patients needing compression following head and neck surgery or radiation therapy



Directional channels

Channel formations facilitate efficient lymph drainage

Foam chips

Open-cell foam chips massage superficial lymphatics


Neck and Mandible PAD

Patient Name:	Contact Name:
Account Name:	Account #:
Account Phone #:	P.O. #:
Ship Name & Address:	Date:

Circumference

Head _____

Neck _____

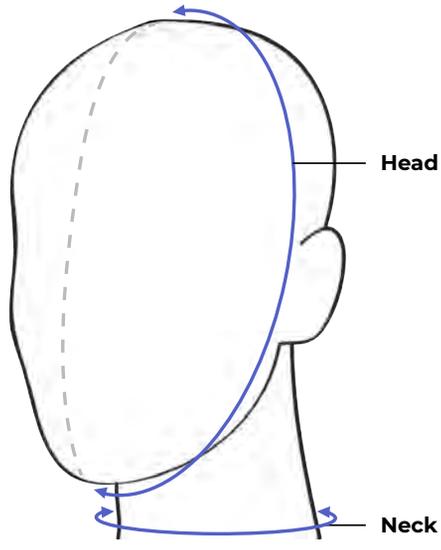

Product Information

Product includes one Neck and Mandible Pad.

<input type="checkbox"/> Neck & Mandible	Color: Black
Size:	Quantity: Item #:

Measuring Instructions

Head:
Measure the circumference by encircling the head; beginning at the top, moving below the chin; and finishing at the top of the head.

Neck:
Measure circumference at largest part of neck.

Neck & Mandible Pad Sizing Chart (Available in Black Only)

	Small	Medium	Large	X Large
Head	46-56cm	52-62cm	56-66cm	58-68cm
Neck	28-33cm	33-38cm	38-43cm	43-58cm
Black	2231-NM	2232-NM	2233-NM	2234-NM


Comprevest

Compression Level

15–25mmHg



Colors

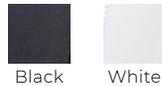

Black

White

Features

The Comprevest offers mild to moderate compression to the chest and back. Designed with features for an improved fit and overall comfort; the Comprevest includes wide straps and adjustable hook closures. The Comprevest is offered in sizes Small to XX-Large with cup sizes A/B–DD.

For added stimulation of the effected area, pair the Comprevest with the Chip Pad Bra. With chip foam sewn into channels, the Chip Pad Bra massages and breaks down fibrotic tissue caused by edema or radiation, while the channeling directs the fluid for drainage. The Chip Pad Bra is designed to be inserted into the interior side pocket of the Comprevest.

Ideal for


- Post-surgical procedures including mastectomy, lumpectomy or radical reconstruction
- Lymphedema

Includes

- Comprevest

Product Selection Tip

May be helpful for post-mastectomy patients

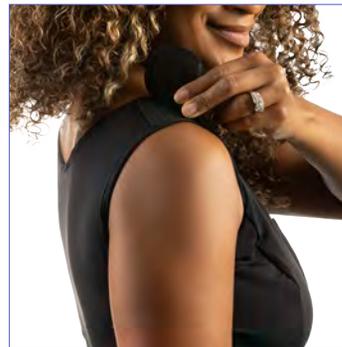

- Anterior hook adjustment for optimal fit
- Straps are easily adjusted
- Interior pockets hold Chip Pad Bra or prosthetic insert

Chip Pads

Chip Pad bra is easily inserted into the interior side pocket of the Comprevest



Chip Pad Bra Half

Cup: 20cmX20cm
Extension: 15cm
Total Length: 35cm
Width: 14cm



Chip Pad Bra Full

Cup: 20cmX20cm
Extension: 30cm
Total Length: 50cm
Width: 14cm


Comprevest

Patient Name:	Contact Name:
Account Name:	Account #:
Account Phone #:	P.O. #:
Ship Name & Address:	Date:

Step One

Measure the circumference in inches at the fullest part of bust.


Bust Size

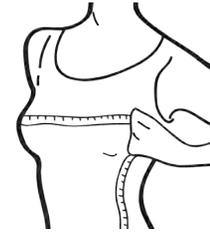

Step Two

Measure the circumference in inches directly under the bust.

Band Size

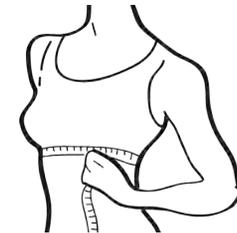

Product Information

Product includes one Comprevest.

Comprevest Cup Size: Item #:
 Color: Black White Band Size: Quantity:

Companion Products (sold separately)

Chip Pad Bra Half Item #:
 Color: Black White Quantity:
 Chip Pad Bra Full Item #:
 Color: Black White Quantity:

Step Three

Determine the cup by subtracting the band from the bust.

Cup Size

2" or less = **A/B**
 2" - 4" = **C/D**
 4" - 6" = **DD/DDD**

Caution: Elastic band contains natural rubber latex.

Chip Pad Bra Size Chart

Size	Single
Black	
Half	2230-BRH
Full	2230-BRF
White	
Half	2280-BRH
Full	2280-BRF

Comprevest Size Chart

	Small	Medium	Large	X Large	XX Large	
Band Size	28-30in	32-34in	36-38in	40-42in	44-46in	
Cup Size						
Black	A/B	2201-BRR-A	2202-BRR-A	2203-BRR-A	2204-BRR-A	2205-BRR-A
	C/D	2201-BRR-C	2202-BRR-C	2203-BRR-C	2204-BRR-C	2205-BRR-C
	DD/DDD	2201-BRR-D	2202-BRR-D	2203-BRR-D	2204-BRR-D	2205-BRR-D
White	A/B	2281-BRR-A	2282-BRR-A	2283-BRR-A	2284-BRR-A	2285-BRR-A
	C/D	2281-BRR-C	2282-BRR-C	2283-BRR-C	2284-BRR-C	2285-BRR-C
	DD/DDD	2281-BRR-D	2282-BRR-D	2283-BRR-D	2284-BRR-D	2285-BRR-D



Chipsleeve Standard

ARM

Compression Level

15–25mmHg



Color



Black

Features

Chipsleeve garments are made for the upper and lower extremity. Filled with chipped foam and designed with directional channels, each garment works to efficiently direct lymph fluid and promote drainage. Patients with fibrotic limbs will enjoy the comfortable fit and massaging effect of the chip foam filled garment.

The Chipsleeve Standard Arm provides mild to moderate gradient compression from the hand to the axilla. The product includes an Oversleeve, which is a high-stretch cover that improves the product fit, increasing the level of compression. Use the garment as a liner for bandaging or as an independent compression solution.

Ideal for

- Lymphedema
- Fibrotic tissue

Includes

- Chipsleeve Arm
- Oversleeve (Black)
- 1 Pair of non-compressive Cotton Liners

Product Selection Tip

May be right for patients requiring a nighttime garment for the management of lymphedema




Chipsleeve Standard **ARM**

Patient Name:	Contact Name:
Account Name:	Account #:
Account Phone #:	P.O. #:
Ship Name & Address:	Date:

Circumference

	Left	Right
C	_____	_____
B1	_____	_____
B	_____	_____
A	_____	_____

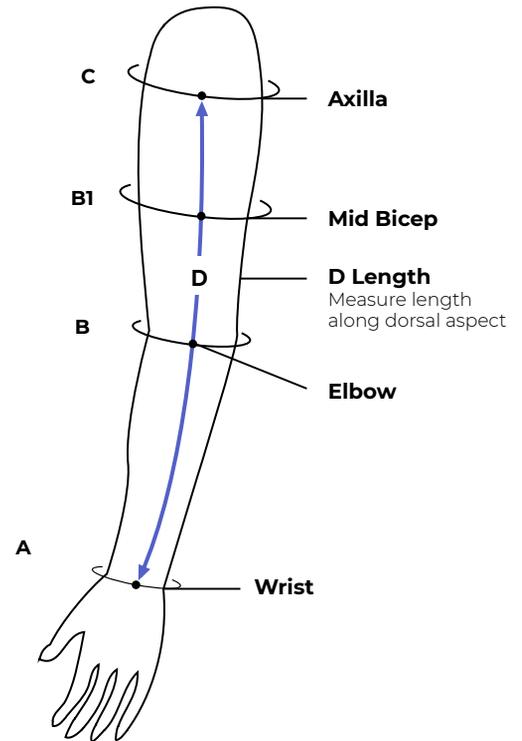

D Length
 Short: 38–42.9cm
 Regular: 43–47.9cm
 Long: 48–53cm

Measure length along dorsal aspect

D	_____	_____
---	-------	-------

Product Information

Product includes one Chipsleeve Standard Arm, one pair of Cotton Liners, and one black Oversleeve.

<input type="checkbox"/> Right Arm	Size:	Item #:
Color: Black	Length:	Quantity:
<input type="checkbox"/> Left Arm	Size:	Item #:
Color: Black	Length:	Quantity:

Oversleeve Arm Size Chart

Color	Small	Medium	Large	X Large	XX Large
Black	1801-ARR	1802-ARR	1803-ARR	1804-ARR	1805-ARR

Chipsleeve Standard Arm Size Chart (A–C) (Available in Black Only)

	Small	Medium	Large	X Large	XX Large
C	23–32cm	28–37cm	33–43cm	39–49cm	44–54cm
B1	22–30cm	26–34cm	30–39cm	35–44cm	39–49cm
B	20–27cm	24–31cm	28–35cm	32–39cm	36–43cm
A	14–17cm	15–18cm	16–19cm	18–20cm	19–22cm

Left

Short	2631-ARS-L	2632-ARS-L	2633-ARS-L	2634-ARS-L	2635-ARS-L
Regular	2631-ARR-L	2632-ARR-L	2633-ARR-L	2634-ARR-L	2635-ARR-L
Long	2631-ARL-L	2632-ARL-L	2633-ARL-L	2634-ARL-L	2635-ARL-L

Right

Short	2631-ARS-R	2632-ARS-R	2633-ARS-R	2634-ARS-R	2635-ARS-R
Regular	2631-ARR-R	2632-ARR-R	2633-ARR-R	2634-ARR-R	2635-ARR-R
Long	2631-ARL-R	2632-ARL-R	2633-ARL-R	2634-ARL-R	2635-ARL-R

Note: This product has a manufacturing and shipping lead time of up to 10 days.



Medaform Standard

ARM

Compression Level

30–40mmHg



Color



Blue &
Black

Features

The Medaform Standard Arm provides firm gradient compression from the hand to the axilla. The Medaform Standard Arm includes the Medahand. The Medahand provides comfortable moderate compression for the hand.

Ideal for

- Lymphedema
- Fibrotic tissue

Includes

- Medaform Standard Arm (upgraded Hand options are available at an additional cost)
- Medahand
- 1 Pair of non-compressive Cotton Liners
- Donning Loop
- Unistrap

Product Selection Tip

May be right for patients with moderate to severe lymphedema



Medaform Standard

ARM

Patient Name:	Contact Name:
Account Name:	Account #:
Account Phone #:	P.O. #:
Ship Name & Address:	Date:

Product Information

Product includes one Medaform, one Medahand, one Unistrap, one Donning Loop, and one pair of Cotton Liners. **Note:** Unless size is specified, the Medahand included will match the size of the Medaform selected.

<input type="checkbox"/> Right Arm	Size:	Item #:
Color: Blue & Black	Length:	Quantity
<input type="checkbox"/> Right Hand	Size:	Item #:
Color: Black	Length:	Quantity

<input type="checkbox"/> Left Arm	Size:	Item #:
Color: Blue & Black	Length:	Quantity
<input type="checkbox"/> Left Hand	Size:	Item #:
Color: Black	Length:	Quantity

Companion Products (sold separately)

<input type="checkbox"/> Medaglove	Size:	Item #:
<input type="checkbox"/> Right <input type="checkbox"/> Left	Length:	Quantity
<input type="checkbox"/> Dorsal Pocket Glove	Size:	Item #:
<input type="checkbox"/> Right <input type="checkbox"/> Left	Length:	Quantity

Hand Size Chart (A– E) (Available in Black Only)

	Small	Medium	Large	X Large
E	15–19cm	19–22cm	22–26cm	26–29cm
A	14–20cm	16–22cm	18–24cm	20–26cm

Black	Medahand			
Right	1201-HDR-R	1202-HDR-R	1203-HDR-R	1204-HDR-R
Left	1201-HDR-L	1202-HDR-L	1203-HDR-L	1204-HDR-L

Black	Medaglove			
Right	1201-GLR-R	1202-GLR-R	1203-GLR-R	1204-GLR-R
Left	1201-GLR-L	1202-GLR-L	1203-GLR-L	1204-GLR-L

Black	Dorsal Pocket Glove			
Right	1201-DGR-R	1202-DGR-R	1203-DGR-R	1204-DGR-R
Left	1201-DGR-L	1202-DGR-L	1203-DGR-L	1204-DGR-L

us_orders@sigvaris.com

P: 800-322-7744 • F: 800-481-5488


Circumference

	Left	Right
C	_____	_____
B1	_____	_____
B	_____	_____
A	_____	_____
E	_____	_____

D Length
 Short: 38–42.9cm
 Regular: 43–47.9cm
 Long: 48–53cm

Measure length along dorsal aspect

D	_____	_____
---	-------	-------

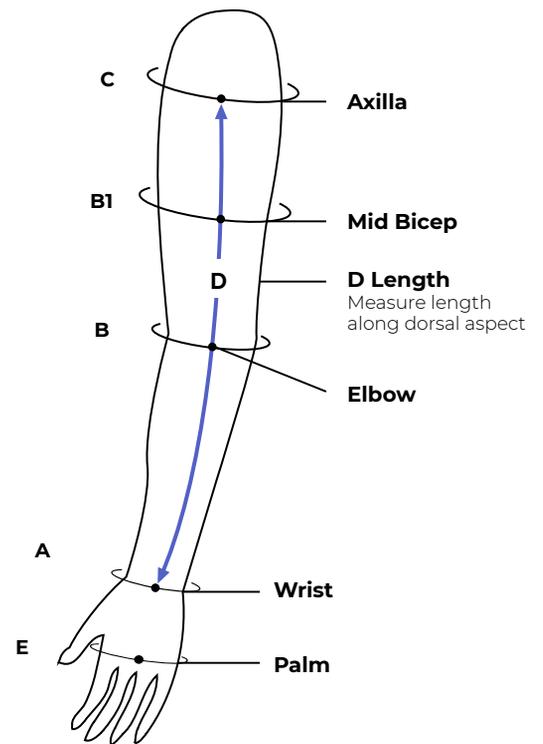

Medaform Arm Size Chart (A– E) (Available in Blue & Black Only)

	Small	Medium	Large	X Large
C	23–32cm	28–37cm	33–43cm	39–49cm
B1	22–30cm	26–34cm	30–39cm	35–44cm
B	20–27cm	24–31cm	28–35cm	32–39cm
A	14–17cm	15–18cm	16–19cm	18–20cm

		Blue & Black			
Left	Short	1201-ARS-L	1202-ARS-L	1203-ARS-L	1204-ARS-L
	Regular	1201-ARR-L	1202-ARR-L	1203-ARR-L	1204-ARR-L
	Long	1201-ARL-L	1202-ARL-L	1203-ARL-L	1204-ARL-L
Right	Short	1201-ARS-R	1202-ARS-R	1203-ARS-R	1204-ARS-R
	Regular	1201-ARR-R	1202-ARR-R	1203-ARR-R	1204-ARR-R
	Long	1201-ARL-R	1202-ARL-R	1203-ARL-R	1204-ARL-R

Note: This product has a manufacturing and shipping lead time of up to 10 days.



For a PDF of the product ordering form, scan the QR code or visit sigvaris.group/mceforms

Compreflex Standard

ARM

Compression Level

20–30mmHg



Colors



Black

Beige

Features

The Compreflex Standard Arm is an inelastic compression device for the arm. The device features an attached inner sleeve for easy application and the patented 90° Flexi-Zone™ for improved movement at the elbow. The Compreflex Standard Arm may be paired with the Dorsal Pocket Glove or Secure Gauntlet for compression in the hand.

Ideal for


- Lymphedema
- Post-surgical edema
- Non-specific edema

Includes

- Compreflex Standard Arm
- 1 Pair of non-compressive Cotton Liners
- 1 Butterfly Strap

Product Selection Tip

May be right for patients with lymphedema of the arm

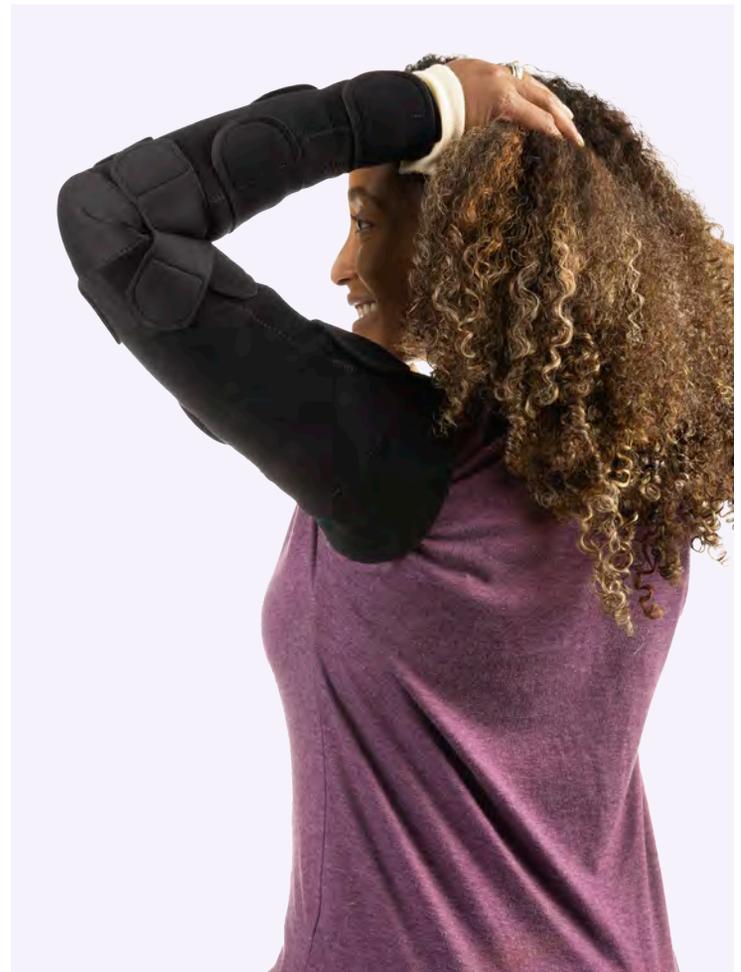

Note: Hand compression wear sold separately



Dorsal Pocket Glove



Secure Gauntlet


Compreflex Standard

ARM

Patient Name:	Contact Name:
Account Name:	Account #:
Account Phone #:	P.O. #:
Ship Name & Address:	Date:

Circumference

	Left	Right
C	_____	_____
B1	_____	_____
B	_____	_____
A	_____	_____
E	_____	_____

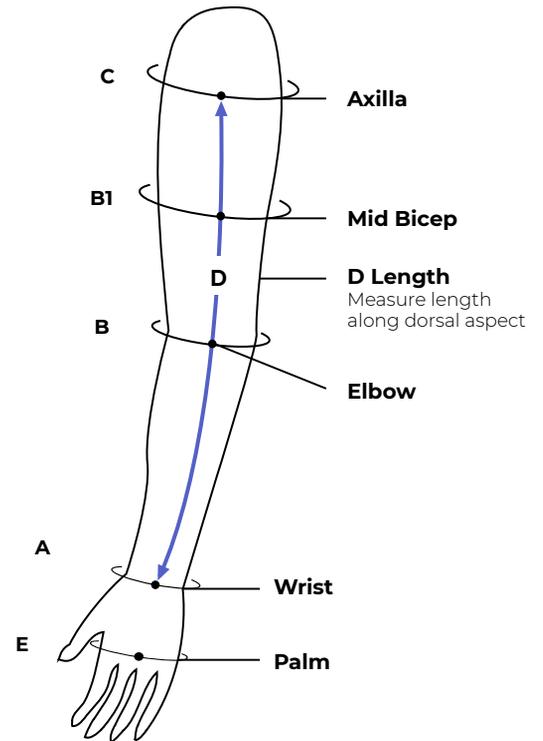

D Length

Short: 38–42.9cm
 Regular: 43–47.9cm
 Long: 48–53cm

Measure length along dorsal aspect

D	_____	_____
---	-------	-------

Product Information

Product includes one Compreflex Standard Arm and one pair of Cotton Liners.

<input type="checkbox"/> Right Arm	Size:	Item #:
Color: <input type="checkbox"/> Black <input type="checkbox"/> Beige	Length:	Quantity:

<input type="checkbox"/> Left Arm	Size:	Item #:
Color: <input type="checkbox"/> Black <input type="checkbox"/> Beige	Length:	Quantity:

Companion Products (sold separately)

<input type="checkbox"/> Dorsal Pocket Glove	Size:	Item #:
--	-------	---------

<input type="checkbox"/> Right Hand	Quantity:
-------------------------------------	-----------

<input type="checkbox"/> Dorsal Pocket Glove	Size:	Item #:
--	-------	---------

<input type="checkbox"/> Left Hand	Quantity:
------------------------------------	-----------

Dorsal Pocket Glove Size Chart (A–E) (Available in Black Only)

	Small	Medium	Large	X Large
E	15–19cm	19–22cm	22–26cm	26–29cm
A	14–20cm	16–22cm	18–24cm	20–26cm
Black				
Right	1201-DGR-R	1202-DGR-R	1203-DGR-R	1204-DGR-R
Left	1201-DGR-L	1202-DGR-L	1203-DGR-L	1204-DGR-L

Compreflex Standard Arm Size Chart (A–C)

	Small	Medium	Large	X Large
C	23–32cm	28–37cm	33–43cm	39–49cm
B1	22–30cm	26–34cm	30–39cm	35–44cm
B	20–27cm	24–31cm	28–35cm	32–39cm
A	14–19cm	15–20cm	16–21cm	18–23cm
Black				
Short	1401-ARS	1402-ARS	1403-ARS	1404-ARS
Regular	1401-ARR	1402-ARR	1403-ARR	1404-ARR
Long	1401-ARL	1402-ARL	1403-ARL	1404-ARL
Beige				
Short	1411-ARS	1412-ARS	1413-ARS	1414-ARS
Regular	1411-ARR	1412-ARR	1413-ARR	1414-ARR
Long	1411-ARL	1412-ARL	1413-ARL	1414-ARL



Coolflex Standard

ARM

Compression Level

20–40mmHg



Color



Black

Features

The Coolflex Standard Arm is ideal for individuals managing lymphedema following mastectomy or other surgical intervention. The garment is made of a lightweight, cooling material that transfers heat and releases moisture away from the skin, leaving it cool and dry. The Coolflex material has two layers which gives it a cushioned feel and helps to protect and provide support the limb.

The product features an easy to grip closure system with hooks for simple, independent donning. The unique pleated elbow design contains a stretch panel to encourage a full range of motion for comfortable all-day wear.

Pair the Coolflex Standard Arm with the Dorsal Pocket Glove for compression to the hand.

Ideal for


- Lymphedema
- Post-surgical edema
- Non-specific edema

Includes

- Coolflex Standard Arm
- 1 Pair of non-compressive Cotton Liners

Product Selection Tip

Select the Coolflex Standard Arm for post mastectomy patients

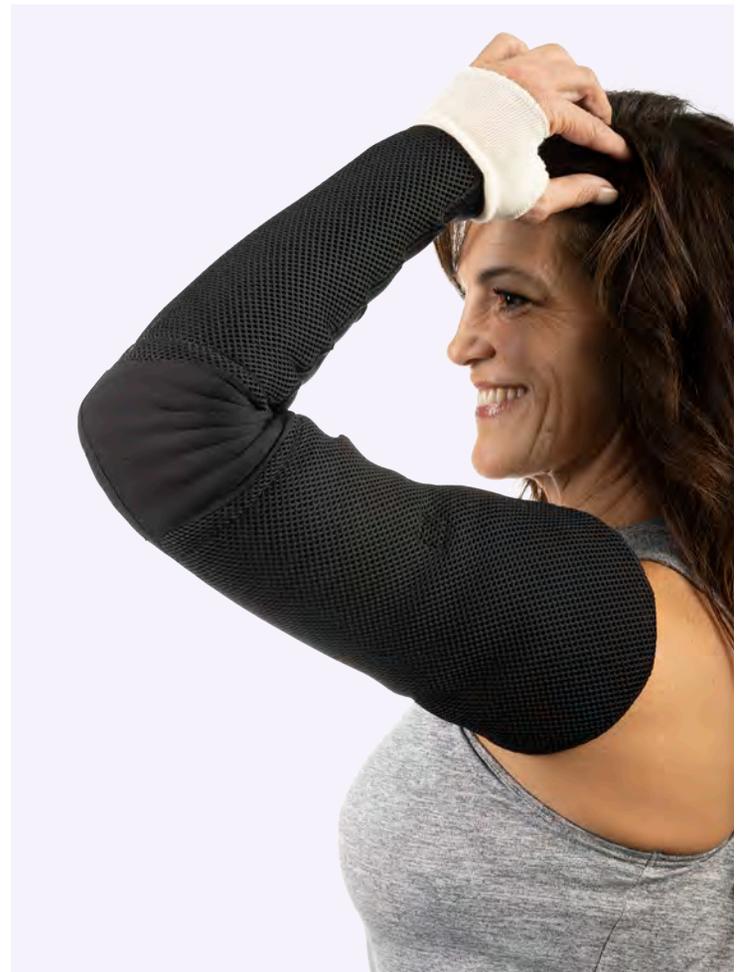

Note: Hand compression wear sold separately



Dorsal Pocket Glove



Secure Gauntlet

Coolflex Standard ARM

Patient Name:	Contact Name:
Account Name:	Account #:
Account Phone #:	P.O. #:
Ship Name & Address:	Date:

Product Information

Product includes one Coolflex Standard Arm and one pair of Cotton Liners.

<input type="checkbox"/> Right Arm	Size:	Item #:
Color: <input type="checkbox"/> Black <input type="checkbox"/> Beige	Length:	Quantity:
<input type="checkbox"/> Left Arm	Size:	Item #:
Color: <input type="checkbox"/> Black <input type="checkbox"/> Beige	Length:	Quantity:

Companion Products (sold separately)

<input type="checkbox"/> Dorsal Pocket Glove	Size:	Item #:
<input type="checkbox"/> Right Hand	Quantity:	
<input type="checkbox"/> Dorsal Pocket Glove	Size:	Item #:
<input type="checkbox"/> Left Hand	Quantity:	

Dorsal Pocket Glove Size Chart (A- E) (Available in Black Only)


	Small	Medium	Large	X Large
E	15-19cm	19-22cm	22-26cm	26-29cm
A	14-20cm	16-22cm	18-24cm	20-26cm
Black				
Right	1201-DGR-R	1202-DGR-R	1203-DGR-R	1204-DGR-R
Left	1201-DGR-L	1202-DGR-L	1203-DGR-L	1204-DGR-L

Circumference

	Left	Right
C	_____	_____
B1	_____	_____
B	_____	_____
A	_____	_____
E	_____	_____

D Length
 Short: 38-42.9cm
 Regular: 43-47.9cm
 Long: 48-53cm
 Measure length along dorsal aspect

D	_____
---	-------

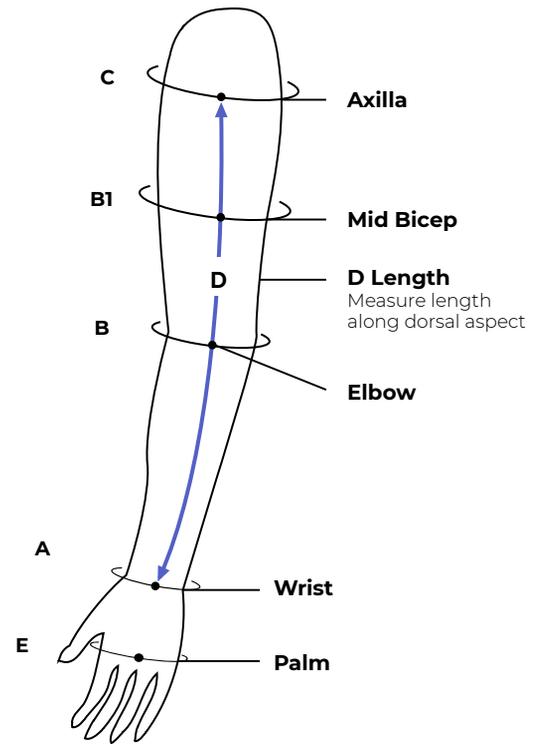

Coolflex Standard Arm Size Chart (A-C) (Available in Black Only)

	Small	Medium	Large	X Large
C	23-32cm	28-37cm	33-43cm	39-49cm
B1	22-30cm	26-34cm	30-39cm	35-44cm
B	20-27cm	24-31cm	28-35cm	32-39cm
A	14-19cm	15-20cm	16-21cm	18-23cm
Right				
Short	2101-ARS-R	2102-ARS-R	2103-ARS-R	2104-ARS-R
Regular	2101-ARR-R	2102-ARR-R	2103-ARR-R	2104-ARR-R
Long	2101-ARL-R	2102-ARL-R	2103-ARL-R	2104-ARL-R
Left				
Short	2101-ARS-L	2102-ARS-L	2103-ARS-L	2104-ARS-L
Regular	2101-ARR-L	2102-ARR-L	2103-ARR-L	2104-ARR-L
Long	2101-ARL-L	2102-ARL-L	2103-ARL-L	2104-ARL-L



sigvaris

Specialty **SECURE ARMSLEEVE, GLOVE, AND GAUNTLET 560**

The Secure Armsleeve combines compression and containment with comfort and durability that's optimal for managing lymphedema.

Compression Level

15–20mmHg

20–30mmHg

Features and Benefits

- More graduated compression and containment than traditional garments
- Full-fingered Glove and Gauntlet available
- Provides UPF 50+ sun protection

Fiber Content

- 78% Nylon, 22% Spandex (Armsleeve)
- 60% Nylon, 40% Spandex (Glove, Gauntlet)
- Non-latex

Best Suited for

- Chronic edema
- Lymphedema that is well-managed
- Post-operative edema
- Post-mastectomy care



Silicone band grip-top



Dusty Rose



Beige



Cocoa



Black



562GM78: Gauntlet Medium in Beige



562FM78: Glove Medium in Beige

Specialty SECURE ARMSLEEVE, GLOVE, AND GAUNTLET 560

Sizing Chart Arm


	Small				Medium				Large				X Large				XX Large			
Wrist	5.5"–6.5" (14–16cm)				6.5"–7" (16–18cm)				7"–8" (18–20cm)				8"–9" (20–23cm)				9"–10" (23–26cm)			
Elbow	Average		Full		Average		Full		Average		Full		Average		Full		Average		Full	
		8"–10" 20–26cm		10"–13.5" 26–34cm		9"–11.5" 23–29cm		11.5"–14.5" 29–37cm		10"–12.5" 26–32cm		12.5"–15.5" 32–40cm		11.5"–13.5" 29–34cm		13.5"–16.5" 34–42cm		12.5"–14.5" 32–37cm		14.5"–17.5" 37–45cm
Upper Arm	Average		Full		Average		Full		Average		Full		Average		Full		Average		Full	
	10"–13.5" 26–34cm		13.5"–15.5" 34–40cm		11.5"–14.5" 29–37cm		14.5"–17" 37–43cm		12.5"–15.5" 32–40cm		15.5"–18" 40–46cm		13.5"–16.5" 34–42cm		16.5"–19" 42–48cm		14.5"–17.5" 37–45cm		17.5"–20.5" 45–52cm	
Arm Length	Short	Long	Short	Long	Short	Long	Short	Long	Short	Long	Short	Long	Short	Long	Short	Long	Short	Long	Short	Long
	≤17" (43cm)	>17" (43cm)	≤17" (43cm)	>17" (43cm)	≤17" (43cm)	>17" (43cm)	≤17" (43cm)	>17" (43cm)	≤17" (43cm)	>17" (43cm)	≤17" (43cm)	>17" (43cm)	≤17" (43cm)	>17" (43cm)	≤17" (43cm)	>17" (43cm)	≤17" (43cm)	>17" (43cm)	≤17" (43cm)	>17" (43cm)
Size	S1	S2	S3	S4	M1	M2	M3	M4	L1	L2	L3	L4	X1	X2	X3	X4	G1	G2	G3	G4

Sizing Chart Hand

	X Small		Small		Medium		Large		X Large		XX Large	
Palm Circumference	6"–6.5" (16–17cm)		7"–7.5" (18–19cm)		8"–8.5" (20–21cm)		8.5"–9" (22–23cm)		9.5"–10" (24–25cm)		10.5"–11" (26–28cm)	
Size	XS		S		M		L		XL		G	

Styles

Sleeve, Glove, and Gauntlet = ●	Dusty Rose 65	Beige 78	Cocoa 84	Black 99
Armsleeve (A)	●	●	●	●
Glove (F)	●	●	●	●
Gauntlet (G)	●	●	●	●

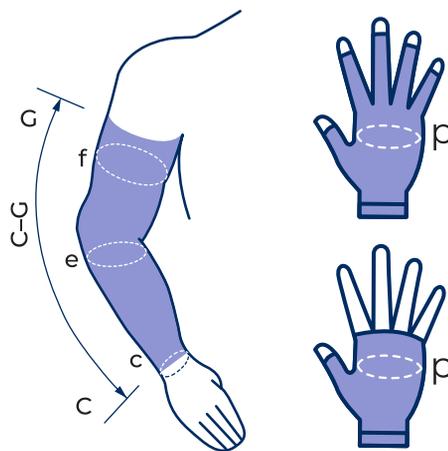

Armsleeve

Product Code Example: **561AL265**

Series	mmHg Level	Style	Size	Color
56	1	A	L2	65

Glove

Product Code Example: **562FM78**

Series	mmHg Level	Style	Size	Color
56	2	F	M	78

Gauntlet

Product Code Example: **562GM78**

Series	mmHg Level	Style	Size	Color
56	2	G	M	78

us_orders@sigvaris.com

P: 800-322-7744 • F: 800-481-5488

Medahand

Variable compression



Color



Black

Features

- Open-cell foam chips to gently massage tissue
- Conforming design to ensure proper fit and gradient compression
- Breathable, moisture-wicking, anti-odor, and soil-release fabric
- Effective as a compression garment or bandage liner where fingers are minimally or not affected by edema



Medaglove

Variable compression



Color



Black

Features

- Conforming design with finger separation ensures proper fit
- Open-cell foam chips gently massage tissue
- Gradient compression from fingers to wrist
- Effective as a compression garment or bandage liner where fingers are affected by edema
- Breathable, moisture-wicking, anti-odor, and soil-release fabric



Dorsal Pocket Glove

Variable compression



Color



Black

Features

- Ultra-thin profile for full hand function during use
- Full finger compression with doffing loops for easy application and removal
- Dorsal pocket holds Chip Pad Dorsum (sold separately) or other therapeutic insert
- Palm pad provides soft, mild compression




Hand Compression Wear

Patient Name:	Contact Name:
Account Name:	Account #:
Account Phone #:	P.O. #:
Ship Name & Address:	Date:

Circumference

	Left	Right
A	_____	_____
E	_____	_____

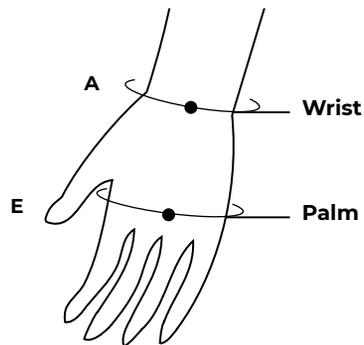

Product Information			
<input type="checkbox"/> Medahand Glove	Color: Black	Item #:	
<input type="checkbox"/> Right <input type="checkbox"/> Left	Size:	Quantity	
<input type="checkbox"/> Medaglove Glove	Color: Black	Item #:	
<input type="checkbox"/> Right <input type="checkbox"/> Left	Size:	Quantity	
<input type="checkbox"/> Dorsal Pocket Glove	Color: Black	Item #:	
<input type="checkbox"/> Right <input type="checkbox"/> Left	Size:	Quantity	

Hand Size Chart (A- E) (Available in Black Only)				
	Small	Medium	Large	X Large
E	15-19cm	19-22cm	22-26cm	26-29cm
A	14-20cm	16-22cm	18-24cm	20-26cm
Medahand				
Right	1201-HDR-R	1202-HDR-R	1203-HDR-R	1204-HDR-R
Left	1201-HDR-L	1202-HDR-L	1203-HDR-L	1204-HDR-L
Medaglove				
Right	1201-GLR-R	1202-GLR-R	1203-GLR-R	1204-GLR-R
Left	1201-GLR-L	1202-GLR-L	1203-GLR-L	1204-GLR-L
Dorsal Pocket Glove				
Right	1201-DGR-R	1202-DGR-R	1203-DGR-R	1204-DGR-R
Left	1201-DGR-L	1202-DGR-L	1203-DGR-L	1204-DGR-L

Measuring Instructions

Hand Size:

Measure circumference at palm and record on line E
Measure circumference at wrist and record on line A

Note: The products on this page have a manufacturing and shipping lead time of up to 10 days.







sigvaris

Lymphedema Wear

LOWER EXTREMITY

The Sigvaris portfolio of lymphedema garments includes products that are designed to support the patient from the decongestive phase to long term management of lymphedema. The product portfolio includes custom garments for individuals that are outside of the standard size range, as well as ready-to-wear products sized to fit most limbs.

Compreshort

Compression Level

10-15mmHg



Color



Black

Features

The Compreshort offers mild, supportive compression for the treatment of abdominal or truncal and upper thigh edema. Available in short and capri lengths, patients with lymphedema and lipedema benefit from the therapeutic compression and support the garment provides.

Custom option available page 106.

Ideal for


- Lymphedema
- Lipedema

Includes

- Compreshort

Product Selection Tip

May be helpful for patients with lipedema, lymphedema, or lipolymphedema in the trunk and thigh

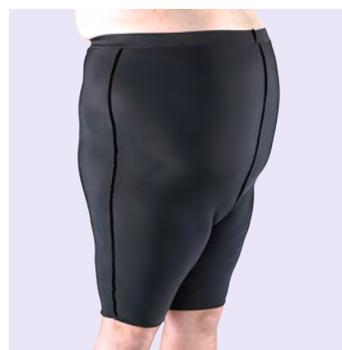

Compreshort

Product available in **short, capri, and custom lengths.**

Compreshort

Patient Name:	Contact Name:
Account Name:	Account #:
Account Phone #:	P.O. #:
Ship Name & Address:	Date:


Circumference

Waist _____

Hip _____

Midthigh _____

Calf _____

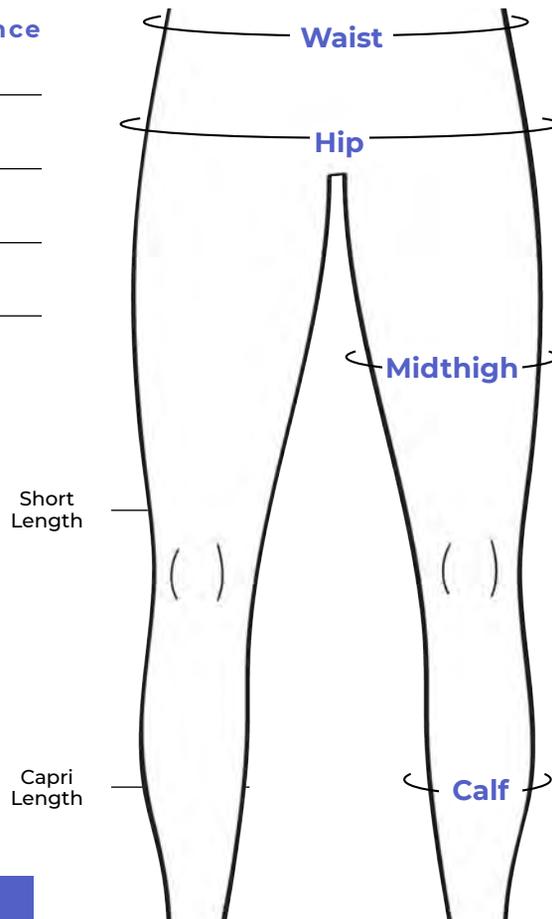

Product Information

Product includes one Compreshort.

Short Length Size: Item #:

Color: Black Quantity:

Capri Length Size: Item #:

Color: Black Quantity:

Caution: Elastic band on size 2XL and larger contains natural rubber latex.

Compreshort Sizing Chart (Waist, Hip, Midthigh)

(Available in Black Only)

	Waist	Hip	Midthigh	Calf	Product SKU (Short)	Product SKU (Capri)
Small	69-76cm	84-91cm	41-46cm	28-32cm	1301-CS	1301-CP
Small Max	76-86cm	91-102cm	41-46cm	28-32cm	1301-CM	1301-CPM
Medium	76-86cm	91-102cm	46-53cm	34-40cm	1302-CS	1302-CP
Medium Max	86-97cm	102-112cm	46-53cm	34-40cm	1302-CM	1302-CPM
Large	86-97cm	102-112cm	53-61cm	40-46cm	1303-CS	1303-CP
Large Max	97-112cm	112-127cm	53-61cm	40-46cm	1303-CM	1303-CPM
X Large	97-112cm	112-127cm	61-71cm	46-53cm	1304-CS	1304-CP
X Large Max	112-132cm	127-147cm	61-71cm	46-53cm	1304-CM	1304-CPM
2X Large	112-132cm	127-147cm	71-81cm	53-61cm	1305-CS	1305-CP
2X Large Max	132-152cm	147-168cm	71-81cm	53-61cm	1305-CM	1305-CPM
3X Large	132-152cm	147-168cm	81-91cm	61-71cm	1306-CS	1306-CP
3X Large Max	152-172cm	168-188cm	81-91cm	61-71cm	1306-CM	1306-CPM
4X Large	152-172cm	168-188cm	91-101cm	71-81cm	1307-CS	1307-CP
4X Large Max	172-193cm	188-208cm	91-101cm	71-81cm	1307-CM	1307-CPM
5X Large	172-193cm	188-208cm	101-112cm	81-91cm	1308-CS	1308-CP

Note: Custom sizes available.

us_orders@sigvaris.com

P: 800-322-7744 • F: 800-481-5488



For a PDF of the product ordering form, scan the QR code or visit

sigvaris.group/mceforms

Compreshape

Compression Level

10-15mmHg



Color



Black

Features

Using premium material that is designed to be worn as outerwear, the Compreshape offers mild, supportive compression for the treatment of abdominal or truncal and upper thigh edema. The stylish waistband provides abdominal coverage and the smooth seams give an attractive appearance. Available in short and capri lengths, patients with lymphedema, edema, and lipedema benefit from the therapeutic compression and support the garment provides.

Ideal for

- Lymphedema
- Edema
- Lipedema

Includes

- Compreshape

Product Selection Tip

May be helpful for patients with lipedema, edema, lymphedema, or lipolymphedema in the trunk and thigh



Smooth seams create stylish finish for outerwear design



Compreshape

Product available in **short** and **capri** lengths.

Compreshape

Patient Name:	Contact Name:
Account Name:	Account #:
Account Phone #:	P.O. #:
Ship Name & Address:	Date:


Circumference

Waist _____

Hip _____

Midthigh _____

Calf _____

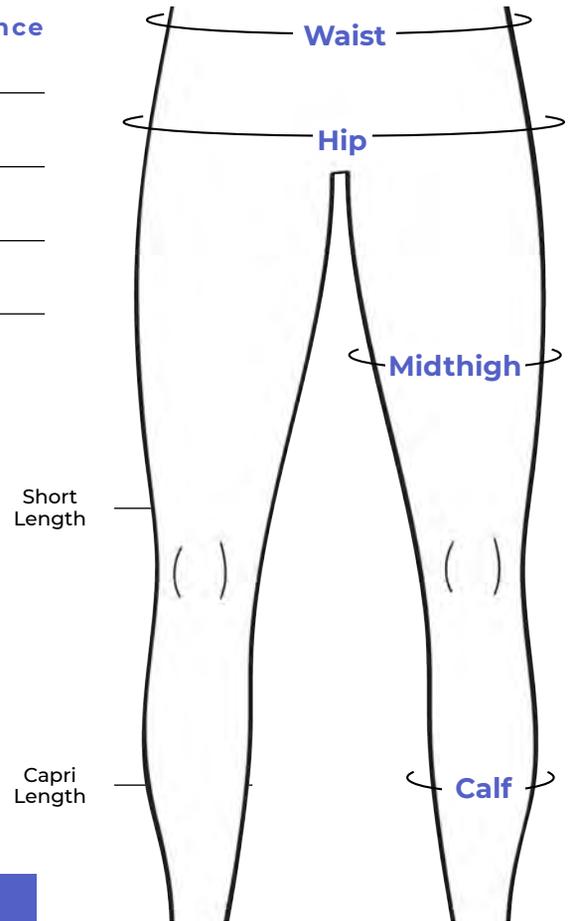

Product Information

Product includes one Compreshape.

Short Length Size: Item #:

Color: Black Quantity:

Capri Length Size: Item #:

Color: Black Quantity:

Caution: Elastic band on size 2XL and larger contains natural rubber latex.

Compreshape Sizing Chart (Waist, Hip, Midthigh)

(Available in Black Only)

	Waist	Hip	Midthigh	Calf	Product SKU (Short)	Product SKU (Capri)
Small	69-76cm	84-91cm	41-46cm	28-32cm	2301-CS	2301-CP
Small Max	76-86cm	91-102cm	41-46cm	28-32cm	2301-CM	2301-CPM
Medium	76-86cm	91-102cm	46-53cm	34-40cm	2302-CS	2302-CP
Medium Max	86-97cm	102-112cm	46-53cm	34-40cm	2302-CM	2302-CPM
Large	86-97cm	102-112cm	53-61cm	40-46cm	2303-CS	2303-CP
Large Max	97-112cm	112-127cm	53-61cm	40-46cm	2303-CM	2303-CPM
X Large	97-112cm	112-127cm	61-71cm	46-53cm	2304-CS	2304-CP
X Large Max	112-132cm	127-147cm	61-71cm	46-53cm	2304-CM	2304-CPM
2X Large	112-132cm	127-147cm	71-81cm	53-61cm	2305-CS	2305-CP
2X Large Max	132-152cm	147-168cm	71-81cm	53-61cm	2305-CM	2305-CPM
3X Large	132-152cm	147-168cm	81-91cm	61-71cm	2306-CS	2306-CP
3X Large Max	152-172cm	168-188cm	81-91cm	61-71cm	2306-CM	2306-CPM
4X Large	152-172cm	168-188cm	91-101cm	71-81cm	2307-CS	2307-CP
4X Large Max	172-193cm	188-208cm	91-101cm	71-81cm	2307-CM	2307-CPM
5X Large	172-193cm	188-208cm	101-112cm	81-91cm	2308-CS	2308-CP

us_orders@sigvaris.com

P: 800-322-7744 • F: 800-481-5488



For a PDF of the product ordering form, scan the QR code or visit

sigvaris.group/mceforms

Chipsleeve

Ideal for patients with lymphedema needing a comfortable garment that promotes lymph fluid movement and drainage.

Compression Level

20-25mmHg



Color




Black

Soft, comfortable, massaging nighttime wear

Sigvaris Chipsleeve garments are made for the upper and lower extremity. Filled with chip foam, each product is designed with directional channels, which work to efficiently direct lymph fluid and promote drainage. Patients with fibrotic limbs will enjoy the comfortable fit and massaging effect of the chip filled garment.

Each Chipsleeve garment is made of four-way stretch fabric. The fabric contains moisture wicking properties and is bacteriostatic to help control odor. Safe for day and nighttime use, the Chipsleeve garment may be used as an underliner for bandaging or as an independent compression garment.

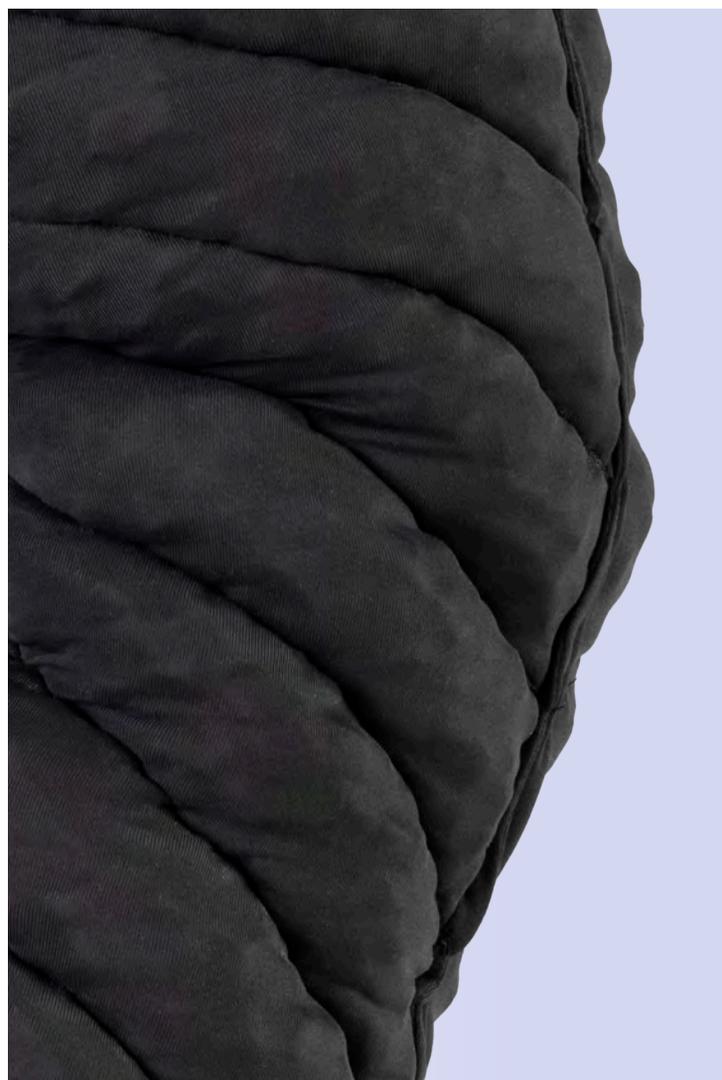

Ideal for

- Lymphedema

Lymphedema Wear

LOWER EXTREMITY



Channel formations

Channels designed to direct the flow of fluid

Chip foam

Open-cell foam chips massage tissue and stimulate fluid movement

Moisture-wicking bacteriostatic material

Each Chipsleeve garment is made with moisture-wicking material

Chipsleeve Standard

CALF & FOOT

Compression Level

20–25mmHg



Color



Black

Features

The Chipsleeve Standard Calf and Foot garment provides mild to moderate gradient compression from the foot to knee. The garment is a single piece that easily glides onto the limb.

Ideal for

- Lymphedema

Includes

- Chipsleeve Standard Calf and Foot
- Oversleeve (Black)
- 1 Pair of non-compressive Basic Liners

Product Selection Tip

May be right for patients with mild to moderate lymphedema in the calf and foot




Chipsleeve Standard

CALF & FOOT

Patient Name:	Contact Name:
Account Name:	Account #:
Account Phone #:	P.O. #:
Ship Name & Address:	Date:

Circumference

	Left	Right
C1	_____	_____
C	_____	_____
B	_____	_____
A	_____	_____

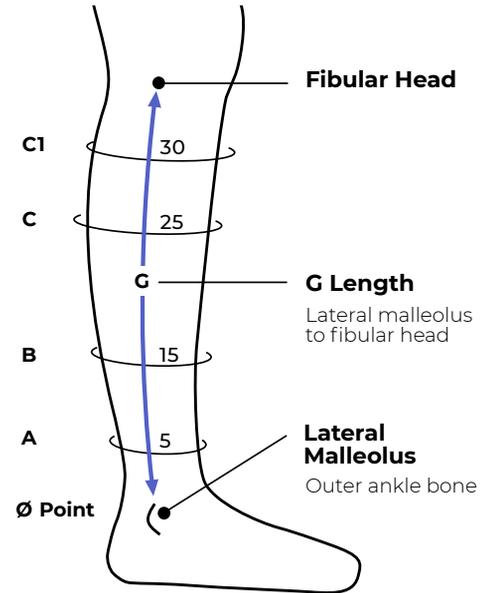

G Length

Regular: 30–35.9cm
Tall: 36–42cm

G	_____	_____
---	-------	-------

Product Information

Product includes one Chipsleeve Standard Calf & Foot, one pair of Basic Liners, and one black Oversleeve.

<input type="checkbox"/> Right Calf	Size:	Item #:
Color: Black	Length:	Quantity:
<input type="checkbox"/> Left Calf	Size:	Item #:
Color: Black	Length:	Quantity:

Oversleeve Calf & Foot Size Chart

Color	Small	Medium	Large	X Large	XX Large
Black	1801-BKR	1802-BKR	1803-BKR	1804-BKR	1805-BKR

Chipsleeve Standard Calf Size Chart (A–C1) (Available in Black Only)

	Black	Medium	Large	X Large	XX Large
C1	32–42cm	37–47cm	42–52cm	51–61cm	58–68cm
C	29–39cm	34–44cm	39–49cm	48–58cm	55–65cm
B	24–34cm	29–39cm	33–43cm	41–51cm	44–55cm
A	20–29cm	21–30cm	25–36cm	32–42cm	33–43cm

Regular	2631-BKR	2632-BKR	2633-BKR	2634-BKR	2635-BKR
Tall	2631-BKT	2632-BKT	2633-BKT	2634-BKT	2635-BKT



Chipsleeve FULL LEG

Compression Level

20–25mmHg



Color



Black

Features

The Chipsleeve Full Leg garment provides mild to moderate gradient compression from the foot to the gluteal fold.

Ideal for

- Lymphedema

Includes

- Chipsleeve Standard Full Leg
- Oversleeve (Black)
- 1 Pair of non-compressive Cotton Liners

Product Selection Tip

May be right for patients with mild to moderate lymphedema in upper and lower portion of the limb



Oversleeve
not shown.

Chipsleeve

FULL LEG

Patient Name:	Contact Name:
Account Name:	Account #:
Account Phone #:	P.O. #:
Ship Name & Address:	Date:


Product Information

Product includes one Chipsleeve Full Leg, one pair of Cotton Liners, and one black Oversleeve.

<input type="checkbox"/> Right Leg	Size:	Item #:
Color: Black	Length:	Quantity:
<input type="checkbox"/> Left Leg	Size:	Item #:
Color: Black	Length:	Quantity:

Circumference

	Left	Right
F	_____	_____
E	_____	_____
D	_____	_____
C1	_____	_____
C	_____	_____
B	_____	_____
A	_____	_____
H	_____	_____

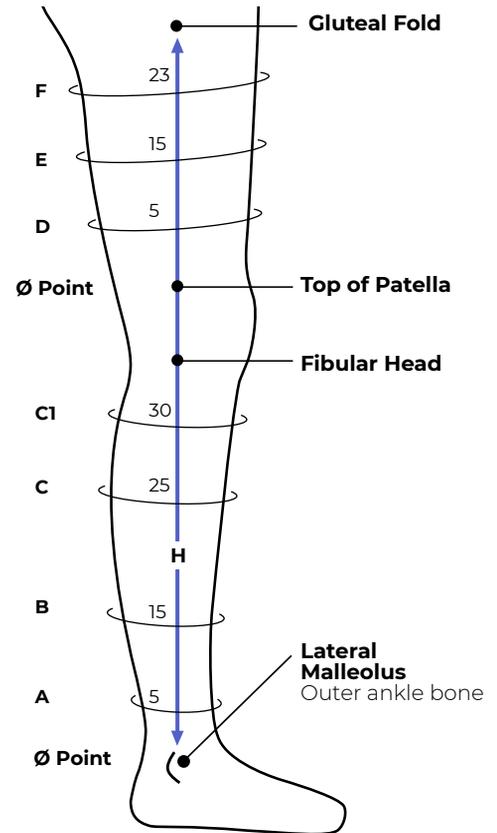

H Length
 Short: 61-70.9cm
 Regular: 71-80.9cm
 Tall: 81-91cm

Oversleeve Full Leg Size Chart

Color	Small	Medium	Large	X Large
Black	1801-THR	1802-THR	1803-THR	1804-THR

Chipsleeve Full Leg Size Chart (A-F) (Available in Black Only)

	Small	Medium	Large	X Large
F	48-58cm	56-66cm	64-74cm	74-84cm
E	43-53cm	51-61cm	58-68cm	68-78cm
D	38-48cm	46-56cm	53-63cm	63-73cm
C1	32-42cm	37-47cm	42-52cm	51-61cm
C	29-39cm	34-44cm	39-49cm	48-58cm
B	24-34cm	29-39cm	33-43cm	41-51cm
A	20-29cm	21-30cm	25-36cm	32-42cm

Note: This product has a manufacturing and shipping lead time of up to 10 days.

		Black			
Left	Short	2631-THS-L	2632-THS-L	2633-THS-L	2634-THS-L
	Regular	2631-THR-L	2632-THR-L	2633-THR-L	2634-THR-L
	Tall	2631-THT-L	2632-THT-L	2633-THT-L	2634-THT-L
Right	Short	2631-THS-R	2632-THS-R	2633-THS-R	2634-THS-R
	Regular	2631-THR-R	2632-THR-R	2633-THR-R	2634-THR-R
	Tall	2631-THT-R	2632-THT-R	2633-THT-R	2634-THT-R



Medaform

Ideal for patients who require high containment. Movable straps allow for adjustment during changes in fluid volume.

Compression Level

30–40mmHg



Color



Blue & Black

Durable, high-grade non-stretch compression

The Medaform garments are for the treatment and management of Stage II and III Lymphedema where high containment levels are desired. The garments are constructed using durable high grade fabric and have a high range of adjustability to accommodate fluctuations in limb fluid volume, common in acute lymphedema. The garments are engineered for durability and are a long-term solution for the treatment and management of lymphedema. Following continuous wear, a Medaform garment may be refreshed by ordering replacement straps or an interior Wavefoam liner.

Other key product features include: the removable Wavefoam liner which works to massage the superficial lymphatics and direct lymph fluid for removal, removable D-ring straps which are highly adjustable and can be positioned for the optimal application of compression, and the easy to grip donning loop which encourage independent application and removal of the product.

The Medaform is available for the lower leg, full leg, and arm, and offers complete limb coverage.



Ideal for

- Lymphedema
- Fibrotic Tissue

Lymphedema Wear

LOWER EXTREMITY



Highly adjustable

Compression system with straps that may be repositioned for a custom fit

Tissue stimulating

Inner Wavefoam liner creates variation in pressure promoting fluid movement

Premium quality

Fabricated using durable medical grade material for containment

Medaform Standard **CALF & FOOT**

Compression Level

30–40mmHg



Color



Blue &
Black

Features

The Medaform Calf and Foot garment provides firm gradient compression from the foot and ankle to the knee. The Medaform Standard Calf garment is paired with the Compreboot Plus Foot for foot compression.

Ideal for

- Lymphedema
- Fibrotic Tissue

Includes

- Medaform Standard Calf
- Compreboot Plus Foot
- 1 Pair of non-compressive Basic Liners

Product Selection Tip

May be right for patients with Stage II or III Lymphedema in calf & foot



Medaform Standard

CALF & FOOT

Patient Name:	Contact Name:
Account Name:	Account #:
Account Phone #:	P.O. #:
Ship Name & Address:	Date:

Product Information

Product includes one Medaform Standard Calf, one Compreboot Plus Foot wrap, and one pair of Basic Liners. Indicate below the required Compreboot Plus Foot size.

<input type="checkbox"/> Right Calf	Size:	Item #:
Color: Blue & Black	Length:	Quantity:
<input type="checkbox"/> Right Foot	Size:	Item #:
Color: Black	Length:	Quantity:
<input type="checkbox"/> Left Calf	Size:	Item #:
Color: Blue & Black	Length:	Quantity:
<input type="checkbox"/> Left Foot	Size:	Item #:
Color: Black	Length:	Quantity:

Note: Liner and additional Compreboot options available, each sold separately.


Compreboot Plus Foot Size Chart (I-K)

(Available in Black Only)

	Small		Medium / Large		X Large / XX Large	
	Regular	Long	Regular	Long	Regular	Long
I	28-36cm	28-34cm	≤39cm	≤39cm	≤44cm	≤44cm
J	22-26cm	22-25cm	≤30cm	≤28cm	≤33cm	≤33cm
K	14-18cm	19-23cm	18-20cm	23-25cm	20-22cm	25-27cm
Black	1801-BPR	1801-BPL	1803-BPR	1803-BPL	1805-BPR	1805-BPL

Circumference

	Left	Right
C1	_____	_____
C	_____	_____
B	_____	_____
A	_____	_____
I	_____	_____
J	_____	_____

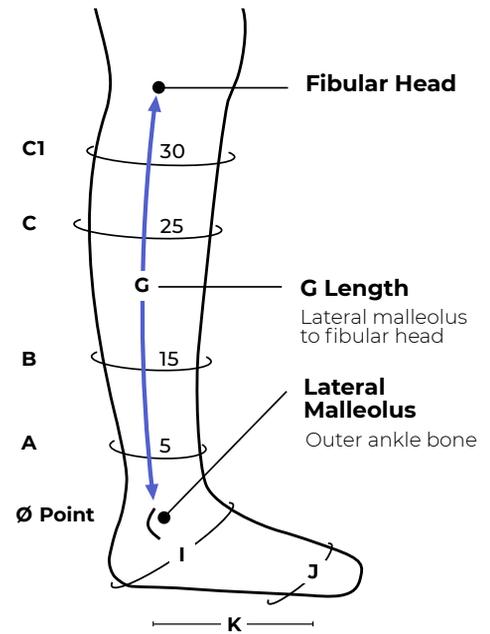

G Length
 Short: 25-29.9cm
 Regular: 30-35.9cm
 Tall: 36-42cm

G	_____	_____
K	_____	_____

K: Measure **medial** length from heel to 1st metatarsal head

Medaform Standard Calf Size Chart (A-C1)

(Available in Blue & Black Only)

	Small	Medium	Large	X Large	XX Large
C1	32-42cm	38-48cm	42-52cm	51-61cm	58-68cm
C	29-39cm	34-44cm	39-49cm	48-58cm	55-65cm
B	24-34cm	29-39cm	33-43cm	41-51cm	44-55cm
A	20-29cm	21-30cm	25-36cm	32-42cm	33-43cm

		Blue & Black				
Left	Short	1201-BKS-L	1202-BKS-L	1203-BKS-L	1204-BKS-L	1205-BKS-L
	Regular	1201-BKR-L	1202-BKR-L	1203-BKR-L	1204-BKR-L	1205-BKR-L
	Tall	1201-BKT-L	1202-BKT-L	1203-BKT-L	1204-BKT-L	1205-BKT-L
Right	Short	1201-BKS-R	1202-BKS-R	1203-BKS-R	1204-BKS-R	1205-BKS-R
	Regular	1201-BKR-R	1202-BKR-R	1203-BKR-R	1204-BKR-R	1205-BKR-R
	Tall	1201-BKT-R	1202-BKT-R	1203-BKT-R	1204-BKT-R	1205-BKT-R



Medaform

FULL LEG

Compression Level

30-40mmHg



Color



Blue &
Black

Features

The Medaform Full Leg garment provides firm gradient compression from the foot to the gluteal fold. The Medaform Full Leg is paired with the Compreboot Plus Foot for compression in the foot.

Ideal for

- Lymphedema
- Fibrotic Tissue

Includes

- Medaform Full Leg
- Compreboot Plus Foot
- 1 Pair of non-compressive Cotton Liners

Product Selection Tip

May be right for patients with Stage II or III Lymphedema requiring a highly adjustable garment for the full leg



Medaform

FULL LEG

Patient Name:	Contact Name:
Account Name:	Account #:
Account Phone #:	P.O. #:
Ship Name & Address:	Date:

Product Information

Product includes one Medaform Full Leg, one Compreboot Plus Foot, and one pair of Cotton Liners. Indicate below the required Compreboot Plus Foot size.

<input type="checkbox"/> Right Leg	Size:	Item #:
Color: Blue & Black	Length:	Quantity:
<input type="checkbox"/> Right Foot	Size:	Item #:
Color: Black	Length:	Quantity:
<input type="checkbox"/> Left Leg	Size:	Item #:
Color: Blue & Black	Length:	Quantity:
<input type="checkbox"/> Left Foot	Size:	Item #:
Color: Black	Length:	Quantity:

Note: Liner and additional Compreboot options available, each sold separately.

Compreboot Plus Foot Size Chart (I-K) (Available in Black Only)

	Small		Medium / Large		X Large / XX Large	
	Regular	Long	Regular	Long	Regular	Long
I	28-36cm	28-34cm	≤39cm	≤39cm	≤44cm	≤44cm
J	22-26cm	22-25cm	≤30cm	≤28cm	≤33cm	≤33cm
K	14-18cm	19-23cm	18-20cm	23-25cm	20-22cm	25-27cm
Black	1801-BPR	1801-BPL	1803-BPR	1803-BPL	1805-BPR	1805-BPL


Note: This product has a manufacturing and shipping lead time of up to 10 days.

Circumference

Left Right

F	_____	_____
E	_____	_____
D	_____	_____
C1	_____	_____
C	_____	_____
B	_____	_____
A	_____	_____
I	_____	_____
J	_____	_____
H	_____	_____
K	_____	_____

H Length
Short: 61-70.9cm
Regular: 71-80.9cm
Tall: 81-91cm

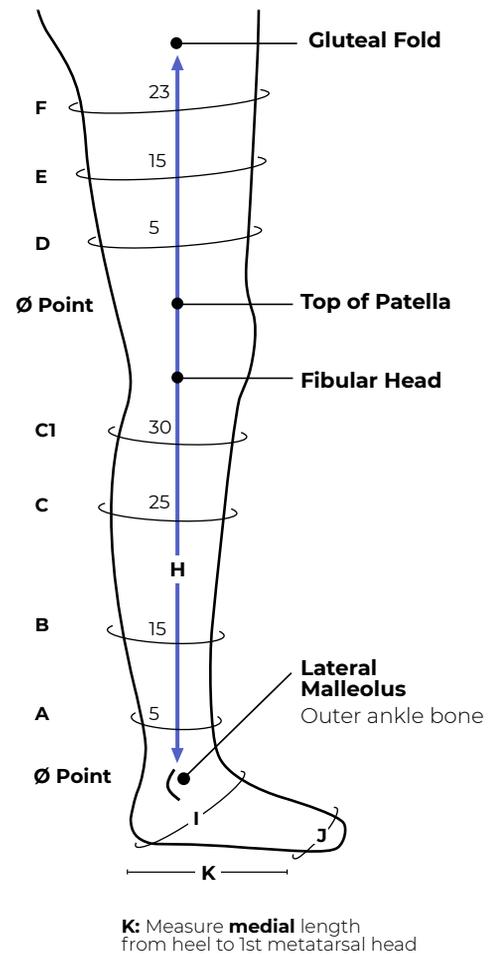

K: Measure **medial** length from heel to 1st metatarsal head

Medaform Full Leg Size Chart (A-F) (Available in Blue & Black Only)

	Small	Small Max	Medium	Medium Max	Large	Large Max	X Large
F	48-58cm	56-66cm	56-66cm	64-74cm	64-74cm	74-84cm	74-84cm
E	43-53cm	51-61cm	51-61cm	58-68cm	58-68cm	68-78cm	68-78cm
D	38-48cm	46-56cm	46-56cm	53-63cm	53-63cm	63-73cm	63-73cm
C1	32-42cm		37-47cm		42-52cm		51-61cm
C	29-39cm		34-44cm		39-49cm		48-58cm
B	24-34cm		29-39cm		33-43cm		41-51cm
A	20-29cm		21-30cm		25-36cm		32-42cm

		Blue & Black						
Left	Short	1201-THS-L	1201-TMS-L	1202-THS-L	1202-TMS-L	1203-THS-L	1203-TMS-L	1204-THS-L
	Regular	1201-THR-L	1201-TMR-L	1202-THR-L	1202-TMR-L	1203-THR-L	1203-TMR-L	1204-THR-L
	Tall	1201-THT-L	1201-TMT-L	1202-THT-L	1202-TMT-L	1203-THT-L	1203-TMT-L	1204-THT-L
Right	Short	1201-THS-R	1201-TMS-R	1202-THS-R	1202-TMS-R	1203-THS-R	1203-TMS-R	1204-THS-R
	Regular	1201-THR-R	1201-TMR-R	1202-THR-R	1202-TMR-R	1203-THR-R	1203-TMR-R	1204-THR-R
	Tall	1201-THT-R	1201-TMT-R	1202-THT-R	1202-TMT-R	1203-THT-R	1203-TMT-R	1204-THT-R



sigvaris

Compreflex Reduce

Ideal garment with easy adjustment for patients in active therapy.

Compression Level

Variable compression



Colors



Black



Beige

Complete Decongestive Therapy (CDT) Aid

The Compreflex Reduce is an innovative compression system used during the decongestive phase of lymphedema therapy. A replacement for tedious short-stretch bandage wrapping, the Compreflex Reduce is easy to size, customizable at the site of care, simple to assemble, and conveniently adjustable for daily wear.

Each garment includes features which promote patient independence and encourages adherence to a planned treatment protocol. These features include an interior stretch panel which holds the garment in place before the straps are closed and visual indicators called Tension Tabs that indicate where straps should be pulled for appropriate levels of tension. To promote gradience and maximize lymphatic drainage, each garment may be combined with foam or other padding materials.

This product is patented.



Ideal for

- Lymphedema
- Patients in the decongestive phase of lymphedema treatment

Tension Tab™



For accurate and consistent application of tension, the Compreflex Reduce garments include Tension Tabs. Tension Tabs are indicators that visually mark the appropriate tension levels as prescribed by the therapist. The Tension Tabs allow for easy readjustment of the garment for effective compression all day long.

Lymphedema Wear

LOWER EXTREMITY

Premium material

High density, durable short-stretch material that compresses and contains the limb

Easy to use

Three-step process to customize the garment during each patient visit

Worry-free compression wear

Visual indicators known as Tension Tabs for setting the correct level of tension with each application

Patient-focused design

Rounded straps for easy-to-grip closure and a secure fit



Compreflex Reduce CALF

Compression Level

Variable compression



Colors



Black

Beige

Features

The Compreflex Reduce Calf provides firm gradient compression to the calf. The garment can be easily modified to create a custom fit and transitions with the patient from active therapy to maintenance.

Pair the Compreflex Reduce Calf with the Compreflex Reduce Foot for complete lower leg compression.

Ideal for

- Lymphedema
- Patients in initial phases of decongestive phase

Includes

- Compreflex Reduce Calf
- 1 Pair of non-compressive Cotton Liners
- Tension Tab Devices
- Marking Kit

Product Selection Tip

May be right for patients in decongestive phase who present with lymphedema in the calf with exclusion of lymphedema in the foot



Lobe Strap

- Pairs with Legassist product to support lobule and hold in place
- Length: 80cm, Width: 8cm

	1 Set
Black	2400-LS
Beige	2410-LS

Compreflex Reduce

CALF


Patient Name:	Contact Name:
Account Name:	Account #:
Account Phone #:	P.O. #:
Ship Name & Address:	Date:

Circumference

	Left	Right
C	_____	_____
A	_____	_____
I	_____	_____
J	_____	_____

G Length
 Regular: 30–33.9cm
 Tall: 34–42cm

G	_____	_____
K	_____	_____

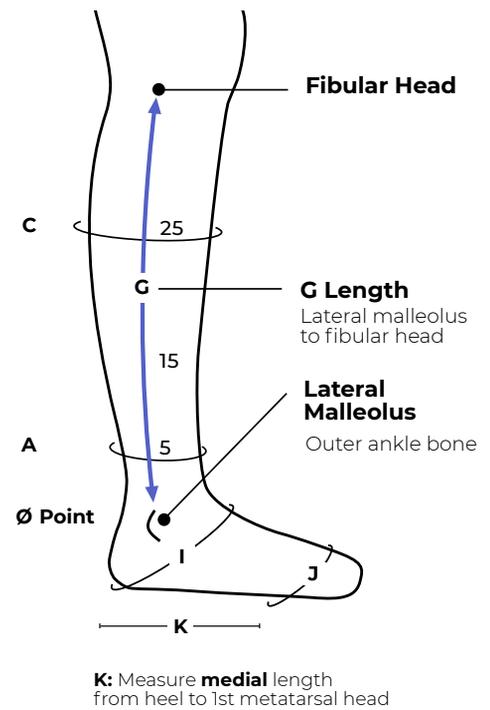

Product Information

Product includes one Compreflex Reduce Calf and one pair of Cotton Liners.

<input type="checkbox"/> Right Calf	Size:	Item #:
Color: <input type="checkbox"/> Black <input type="checkbox"/> Beige	Length:	Quantity:
<input type="checkbox"/> Left Calf	Size:	Item #:
Color: <input type="checkbox"/> Black <input type="checkbox"/> Beige	Length:	Quantity:

Companion Products (sold separately)

<input type="checkbox"/> Compreflex Reduce Foot	Color: <input type="checkbox"/> Black <input type="checkbox"/> Beige
Size:	Item #: Quantity:
Size:	Item #: Quantity:

Note: Compreboot options available, each sold separately.

Measuring Instructions

Compreflex Reduce Calf:

Measure circumferences where indicated, beginning 5cm above the lateral malleolus (outer ankle bone), continuing to 25cm. Record in the corresponding lines A and C (Left/Right). Measure length from the lateral malleolus (outer ankle bone) to lateral aspect of the fibular head, and record in box G.

Compreflex Reduce Foot:

Measure circumference of ankle bend and heel and record on line I. Measure circumference across metatarsal heads and record on line J. Measure length on medial side of foot from heel to 1st metatarsal head and record on line K.

Compreflex Reduce Calf Size Chart (A-C)

	Range 1	Range 3
C	<65cm	<90cm
A	<55cm	<75cm

Black		
Regular	2401-BKR	2403-BKR
Tall	2401-BKT	2403-BKT
Beige		
Regular	2411-BKR	2413-BKR
Tall	2411-BKT	2413-BKT

Compreflex Reduce Foot Size Chart (I-K)

	Range 1	Range 3
I	30-44cm	39-54cm
J	24-33cm	30-43cm
K	16-25cm	20-27cm
Black		
	2801-BTL	2803-BTL
Beige		
	2811-BTL	2813-BTL

us_orders@sigvaris.com

P: 800-322-7744 • F: 800-481-5488



For a PDF of the product ordering form, scan the QR code or visit

sigvaris.group/mceforms

Compreflex Reduce **THIGH**

Compression Level

Variable compression



Colors



Black

Beige

Features

The Compreflex Reduce Thigh provides firm gradient compression to the thigh. The garment is easily modified to create a custom fit and transitions with the patient from decongestive phase to maintenance. For added support and comfort, the product includes a Belt and Hip Attachment used to secure the product in place.

Pair the Compreflex Reduce Thigh with the Compreflex Reduce Knee for upper leg compression.

Ideal for

- Lymphedema
- Patients in all stages of decongestive phase

Includes

- Compreflex Reduce Thigh
- 1 Pair of non-compressive Cotton Liners
- Belt and Hip Attachment
- Tension Tab Devices
- Marking Kit

Product Selection Tip

May be right for patients in decongestive phase who present with lymphedema in the thigh



Lobe Strap

- Pairs with Legassist product to support lobule and hold in place
- Length: 80cm, Width: 8cm

	1 Set
Black	2400-LS
Beige	2410-LS

Compreflex Reduce THIGH

Patient Name:	Contact Name:
Account Name:	Account #:
Account Phone #:	P.O. #:
Ship Name & Address:	Date:


Circumference

Left Right

F _____

D _____

C1 _____

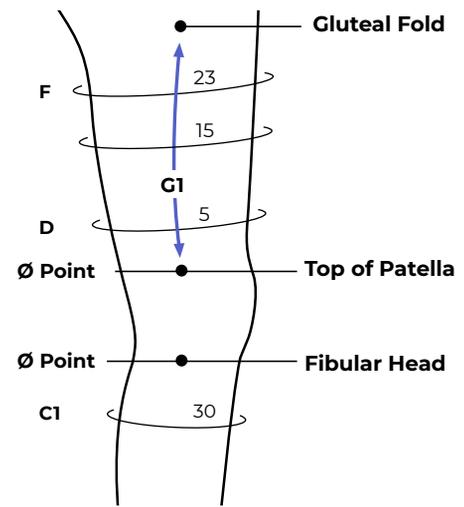

G1 Length
 Regular: 20–24.9cm
 Tall: 25–30cm

G1 _____

Product Information

Product includes one Compreflex Reduce Thigh wrap, one pair of Cotton Liners, and one Belt and Hip attachment.

Right Thigh Size: _____ Item #: _____
 Color: Black Beige Length: _____ Quantity: _____

Left Thigh Size: _____ Item #: _____
 Color: Black Beige Length: _____ Quantity: _____

Companion Products (sold separately)

Compreflex Reduce Knee Color: Black Beige

Size: _____ Item #: _____ Quantity: _____

Size: _____ Item #: _____ Quantity: _____

Measuring Instructions

Compreflex Reduce Thigh:

Measure circumferences where indicated on the measurement chart, beginning 5cm above the top of patella and continue 23cm up the thigh, just below the gluteal fold. Record in the corresponding lines D and F. Measure length from the lateral aspect top of patella to the gluteal fold, and record in box G1.

Compreflex Reduce Knee:

Measure circumferences where indicated, beginning 5cm below the lateral aspect fibular head. Continue 5cm above the lateral aspect top of patella. Record in the corresponding lines C1 and D.

Compreflex Reduce Thigh Size Chart

	Range 1	Range 3
F	<80cm	<125cm
D	<65cm	<110cm

Black

	Black	
	Regular	Tall
Left	2401-TCR-L	2403-TCR-L
Right	2401-TCR-R	2403-TCR-R

Beige

	Beige	
	Regular	Tall
Left	2411-TCR-L	2413-TCR-L
Right	2411-TCR-R	2413-TCR-R

Compreflex Reduce Knee Size Chart

	Range 1	Range 3
D	<75cm	<115cm
C1	<65cm	<90cm

Black

Black	2401-KP	2403-KP
-------	---------	---------

Beige

Beige	2411-KP	2413-KP
-------	---------	---------



Compreflex Reduce

KNEE

Compression Level

Variable compression



Colors



Black

Beige

Features

The Compreflex Reduce Knee provides firm gradient compression to the knee. The garment is easily modified to create a custom fit and transitions with the patient from decongestive phase to maintenance.

Pair the Compreflex Reduce Knee with the Compreflex Reduce Foot, Calf, and Thigh for total leg compression.

Ideal for

- Lymphedema
- Patients in initial phases of decongestive phase

Includes

- Compreflex Reduce Knee
- Butterfly Strap
- Tension Tab Devices
- Marking Kit

Product Selection Tip

May be right for patients in decongestive phase who present with lymphedema in the knee



Lobe Strap

- Pairs with Legassist product to support lobule and hold in place
- Length: 80cm, Width: 8cm

	1 Set
Black	2400-LS
Beige	2410-LS

Compreflex Reduce

KNEE

Patient Name:	Contact Name:
Account Name:	Account #:
Account Phone #:	P.O. #:
Ship Name & Address:	Date:

Product Information

Product includes one Compreflex Reduce Knee wrap and one Butterfly Strap

Right Knee Size: Item #:

Color: Black Beige Length: Quantity:

Left Knee Size: Item #:

Color: Black Beige Length: Quantity:

Companion Products (sold separately)

Compreflex Reduce Thigh Color: Black Beige

Size: Item #: Quantity:

Size: Item #: Quantity:

Measuring Instructions

Compreflex Reduce Knee:

Measure circumferences where indicated, beginning 5cm below the lateral aspect fibular head. Continue 5cm above the lateral aspect top of patella. Record in the corresponding lines C1 and D.

Compreflex Reduce Thigh:

Measure circumferences where indicated on the measurement chart, beginning 5cm above the top of patella and continue 23cm up the thigh, just below the gluteal fold. Record in the corresponding lines D and F. Measure length from the lateral aspect top of patella to the gluteal fold, and record in box G1.

Circumference

Left Right

F _____

D _____

C1 _____


G1 Length

Regular: 20–24.9cm

Tall: 25–30cm

*Top of patella to
gluteal fold*

G1 _____

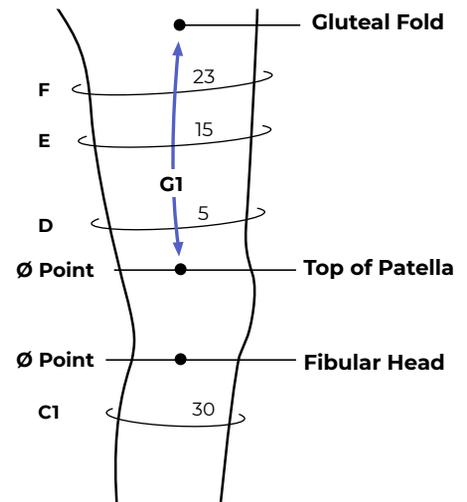

Compreflex Reduce Knee Size Chart

	Range 1	Range 3
D	<75cm	<115cm
C1	<65cm	<90cm
Black	2401-KP	2403-KP
Beige	2411-KP	2413-KP

Compreflex Reduce Thigh Size Chart

	Range 1	Range 3
F	<80cm	<125cm
D	<65cm	<110cm
Black		
Left	Regular	2401-TCR-L 2403-TCR-L
	Tall	2401-TCT-L 2403-TCT-L
Right	Regular	2401-TCR-R 2403-TCR-R
	Tall	2401-TCT-R 2403-TCT-R
Beige		
Left	Regular	2411-TCR-L 2413-TCR-L
	Tall	2411-TCT-L 2413-TCT-L
Right	Regular	2411-TCR-R 2413-TCR-R
	Tall	2411-TCT-R 2413-TCT-R



Compreflex Reduce **FOOT**

Compression Level

Variable compression



Colors



Black

Beige

Features

The Compreflex Reduce Foot provides firm gradient compression for the foot. The garment is easily modified to create a custom fit and transitions with the patient from decongestive phase to maintenance.

Ideal for

- Lymphedema
- Patients in initial phases of decongestive phase

Includes

- Compreflex Reduce Foot
- Butterfly Strap
- Marking Kit

Product Selection Tip

May be right for patients in decongestive phase who present with lymphedema in the foot



Compreflex Reduce

FOOT


Patient Name:	Contact Name:
Account Name:	Account #:
Account Phone #:	P.O. #:
Ship Name & Address:	Date:

Circumference

Left **Right**

I _____ _____ Ø Point

J _____ _____

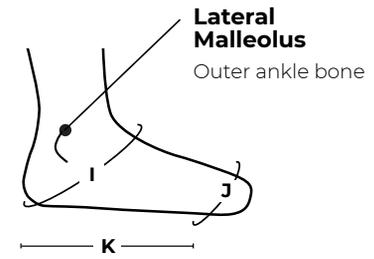

Length

K _____ _____

K: Measure **medial** length from heel to 1st metatarsal head

Product Information

Product includes one Compreflex Reduce Foot.

Right Foot Size: Item #:

Color: Black Beige Quantity:

Left Foot Size: Item #:

Color: Black Beige Quantity:

Note: Compreflex Reduce Calf options available, each sold separately.

Compreflex Reduce Foot Size Chart (I-K)

	Range 1	Range 3
I	30-44cm	39-54cm
J	24-33cm	30-43cm
K	16-25cm	20-27cm
Black	2801-BTL	2803-BTL
Beige	2811-BTL	2813-BTL

Measuring Instructions

Compreflex Reduce Foot:

Measure circumference of ankle bend and heel and record on line I. Measure circumference across metatarsal heads and record on line J. Measure length on medial side of foot from heel to 1st metatarsal head and record on line K.







sigvaris

Bandages and Therapy Supplies

SIGVARIS GROUP offers a full range of bandages, therapy supplies, and product accessories that support the patient and clinician. We offer a range of short-stretch bandages used during CDT as well as foam and other materials.

Our product accessories are designed to improve the function of the inelastic wraps and lymphedema garments. Our garments range from items that work to extend the life of the garment to those that improve comfort and maximize the potential for patient adherence.

sigvaris

Short Stretch Bandages

Biaform

Case holds 20 rolls

Premium short stretch bandage. Product offers superior stretch and conforming abilities and is durable.



Size	Qty	Product Code
6cm x 5m	1 Case	5052-C
8cm x 5m	1 Case	5053-C
10cm x 5m	1 Case	5054-C
12cm x 5m	1 Case	5055-C
10cm x 10m	1 Case	5056-C
12cm x 10m	1 Case	5057-C

Stockinette & Dressing Retainers

Cotton Stockinette

Ideal for use as an underliner with nighttime bandaging or compression garments. Cotton Stockinette is available in a range of sizes to fit almost any application. Product is non-compressive, 100% Cotton.



Size	Qty	Product Code
3" x 25yds	1 Roll	5403
4" x 25yds	1 Roll	5404
5" x 25yds	1 Roll	5405
6" x 25yds	1 Roll	5406
8" x 25yds	1 Roll	5408

Biagrip

Compressive stockinette, which may be washed and reused to augment compression garments. Product may be cut to size.



Size	Qty	Product Code
D 3" x 11yds	1 Roll	5433
E 3.5" x 11yds	1 Roll	5434
F 4" x 11yds	1 Roll	5435
G 5" x 11yds	1 Roll	5436
J 7" x 11yds	1 Roll	5437
K 8" x 11yds	1 Roll	5438

Gauze & Tape

Elastosoft

Soft cotton elastic gauze. Ideal for wrapping fingers & toes.



Size	Qty	Product Code
1" x 4.1yds	1 Bag / 24 Rolls	5601
2" x 4.1yds	1 Bag / 12 Rolls	5602
3" x 4.1yds	1 Bag / 12 Rolls	5603
4" x 4.1yds	1 Bag / 12 Rolls	5604

Durapore®

3M Silk Surgical Tape

Exceptionally durable, easy to cut to size surgical tape. Optimal for securing bandages.



Size	Qty	Product Code
1" x 10yds	1 Box / 12 Rolls	5691 - B
2" x 10yds	1 Box / 6 Rolls	5692 - B

Bandages and Therapy Supplies

Foam

Grayfoam

Moderately rigid open-cell foam helps to soften fibrotic areas. The product provides firm, stable support under Biaform short stretch bandages.



Size	Qty	Product Code
1/4"x18x72"	1 Sheet	5301
1/2"x18x72"	1 Sheet	5302

Wavefoam

One sided

High-pressure ridges gently massage thickened fibrotic areas, while low-pressure channels allow pathways for fluid evacuation.




Size	Qty	Product Code
1.2cm x 69cm x 90cm	1 Sheet	5202

Wavefoam II

Two sided

High-pressure ridges gently massage thickened fibrotic areas, while low-pressure channels allow pathways for fluid evacuation. Ridged or wave pattern is on both sides of the sheet.

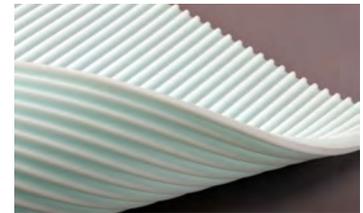

Size	Qty	Product Code
1.2cm x 69cm x 90cm	1 Sheet	5212

Padding

Biasoft

Case holds 10 rolls

Soft, durable, and seamless fleece padding for use under short stretch bandages. The product may be machine washed and dried for repeated use.



Size	Qty	Product Code
6cm x 3m	1 Roll	5251
6cm x 3m	1 Case	5251-C
10cm x 3m	1 Roll	5252
10cm x 3m	1 Case	5252-C
15cm x 3m	1 Roll	5253
15cm x 3m	1 Case	5253-C

Non-Woven

Soft Synthetic Padding Case holds 30 rolls

Optimal for use under short stretch bandages.



Size	Qty	Product Code
10cm x 3.5m (4" x 3.8yds)	1 Case	5661-C



sigvaris

Chip Pad Dorsum

Ideal for use on dorsum of hand or foot. Helps with softening scars from carpal tunnel surgery, or underneath compression bandages over pockets of lymphedema. Fits into Dorsal Pocket Glove.



Size	Single
Small	2231-DO (8cm x 9cm)
Large	2233-DO (9cm x 14cm)

Chip Pad Axilla

Use under arm for treatment of lymphedema, especially over areas of post-operative scarring.



Size	Single
Small	2231-AX (10cm x 19cm)
Large	2233-AX (18cm x 27cm)

Chip Pad Radial Half


Radial pathways create low pressure channels to promote fluid drainage. Use over affected areas or pockets of edema.



Size	Single
One size fits all	2230-RH (15cm x 25cm)

Chip Pad Radial Full

Versatile, effective design is ideal for post-surgical scarring and/or pockets of lymphedema.

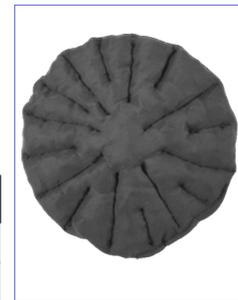

Size	Single
Small	2231-RF (18cm diameter)
Large	2233-RF (23cm diameter)

Chip Pad Bra


Ideal for post mastectomy/lumpectomy edema. Pair with Comprevest.

Size	Single
Black	
Half	2230-BRH
Full	2230-BRF
White	
Half	2280-BRH
Full	2280-BRF



Chip Pad Bra Half

Cup: 20cmX20cm
Extension: 15cm
Total Length: 35cm
Width: 14cm

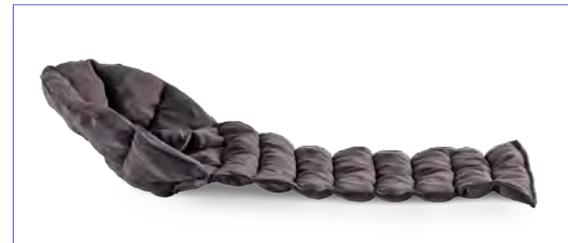

Chip Pad Bra Full

Cup: 20cmX20cm
Extension: 30cm
Total Length: 50cm
Width: 14cm


Note: Due to the nature of this product, dimensions may vary.

Bandages and Therapy Supplies

Chip Pads Genifit™

The Genifit series offers gender-specific Chip Pads designed for the treatment of genital edema. Each Chip Pad includes open-cell foam chips sewn in channel formations to assist in massaging tissue and facilitate drainage. The Chip Pads are made with moisture-wicking fabric and may be machine washed and dried for repeated use.

Male

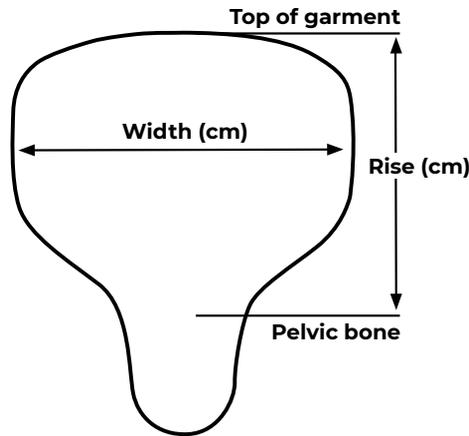

Tall


The rise of the tall male chip pad ranges in length from 32–44cm. See the measurement guide for reference.



Short

The rise of the short male chip pad ranges in length from 14–22cm. See the measurement guide for reference.

Female

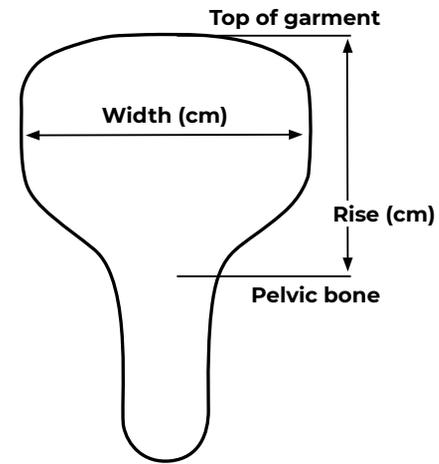

Tall

The rise of the tall Genifit female chip pad ranges in length from 32–41cm. See the measurement guide for reference.



Short

The rise of the short female chip pad ranges in length from 14–20cm. See the measurement guide for reference.

Male Genifit

	Small	Medium	Large	X Large	XX Large
Width	28cm	32cm	36cm	40cm	44cm
Scrotal Circ.	26cm	32cm	38cm	44cm	51cm
Rise – Short	14cm	16cm	18cm	20cm	22cm
Rise – Tall	32cm	35cm	38cm	41cm	44cm
Short	2231-GMS	2232-GMS	2233-GMS	2234-GMS	2235-GMS
Tall	2231-GMT	2232-GMT	2233-GMT	2234-GMT	2235-GMT

Female Genifit

	Small	Medium	Large	X Large
Width	28cm	32cm	36cm	40cm
Rise – Short	14cm	16cm	18cm	20cm
Rise – Tall	32cm	35cm	38cm	41cm
Short	2231-GFS	2232-GFS	2233-GFS	2234-GFS
Tall	2231-GFT	2232-GFT	2233-GFT	2234-GFT

Note: Due to the nature of this product, dimensions may vary.



sigvaris

Complete Liner

10–15mmHg in foot, ankle, and calf

- Comfortable knee-high liners
- Drirelease® moisture wicking and odor control
- Extra thick sole from toe to heel for added comfort
- Silver ionized treatment for bacteriostatic properties



Compreflex Complete Calf Size Chart (A–C)

SIZE	Ankle Circ. (A)	Calf Circ. (C)	Women's Shoe	Men's Shoe	Pair (Black)	Pair (Beige)
Small	19–26cm	29–39cm	<7½	<6½	1601-CK-FI	1611-CK-FI
Medium	21–30cm	34–44cm	7.5–9.5	6.5–8.5	1602-CK-FI	1612-CK-FI
Large	26–36cm	40–50cm	10–12	9–11	1603-CK-FI	1613-CK-FI
X Large	31–41cm	46–56cm	12.5–15	11.5–14	1604-CK-FI	1614-CK-FI
XX Large	36–46cm	54–64cm	12.5–15	11.5–14	1605-CK-FI	1615-CK-FI

Transition Liner

10–15mmHg in foot and ankle

- Soft, flat seam knee-high liners
- Drirelease moisture wicking and odor control
- Extra thick sole from toe to heel for added comfort
- Silver ionized treatment for bacteriostatic properties



Compreflex Transition Calf Size Chart (A–C)

SIZE	Ankle Circ. (A)	Calf Circ. (C)	Women's Shoe	Men's Shoe	Pair (Black)	Pair (Beige)
Small	19–26cm	29–39cm	<7½	<6½	1601-CK	1611-CK
Medium	21–30cm	34–44cm	7.5–9.5	6.5–8.5	1602-CK	1612-CK
Large	26–36cm	40–50cm	10–12	9–11	1603-CK	1613-CK
X Large	31–41cm	46–56cm	12.5–15	11.5–14	1604-CK	1614-CK
XX Large	36–46cm	54–64cm	12.5–15	11.5–14	1605-CK	1615-CK

Basic Liner

Non-compressive

- Non-compressive, protective, moisture-wicking liner
- Extra thick sole from toe to heel for added comfort
- Silver ionized treatment for bacteriostatic properties

Size: One size fits all

1 Pair	10 Pair Pack
1700-NK1	1700-NKC



Cotton Socks & Liner

Non-compressive

Soft, seamless Cotton Liners are non-compressive and can be worn under all compression products



Style	Size	1 Pack	3 Pack	10 Pack
Sock	S–XXL	1795-BK1 (1 Pair)	1795-BK3 (3 Pairs)	1795-BKB (10 Pairs)
Thigh High	S / M / L	1793-TH (One liner per pack)		
Thigh High	XL	1794-TH (One liner per pack)		
Arm	S–XL	1794-AR (One liner per pack)		

Foam Liner

- High density, open-cell foam with breathable, moisture-wicking mesh lining
- May be paired with inelastic wraps or worn separately

Size	Regular	Tall
Small	1101-FLR	1101-FLT
Medium	1102-FLR	1102-FLT
Large	1103-FLR	1103-FLT
X Large	1104-FLR	1104-FLT
XX Large	1105-FLR	1105-FLT



Liner pictured, not included. Use chart of Comprefit Standard Calf & Foot on page 22.

Bandages and Therapy Supplies

Ankle Pad

- Wave foam ankle padding
- Includes stretch band to secure padding in place
- Constructed using moisture-wicking fabric



Size	Each
S / M / L	2233-AP

Unistrap



- Neoprene wrap for added hand compression
- Length: 80cm, Width: 6.5cm

Each
1902

Comprefit Strap Extenders

- Adds 10cm of adjustable range to a single strap
- 1 set extends 1 strap
- Strap extender works only with the Comprefit brand of products



	1 Set	10 Set Pack
Black	1100-SE1	1100-SEB
Beige	1110-SE1	1110-SEB

Compreflex Strap Extenders

- Adds 10cm of adjustable range to a single strap
- 1 set extends 1 strap
- Strap extender works only with the Compreflex brand of products



	1 Set	10 Unit Pack
Black	1401-SE1	1401-SEB
Beige	1411-SE1	1411-SEB

Gripper Patches

- Adds traction and grip to foot compression garments
- Hook and loop bottom makes patches easy to apply
- 1 set includes two patches



1 Set
1800-GP

Lobe Strap

- Pairs with Legassist product to support lobule and hold in place
- Length: 80cm, Width: 8cm



	1 Set
Black	2400-LS
Beige	2410-LS

Accutab Devices

Compression measuring device design for use with the Compreflex line of inelastic garments for the lower extremity. See page 6 for more information about the Accutab device.



One set includes five Accutab devices

Size	1 Set
Small	1491-TB1
Medium	1492-TB1
Large	1493-TB1
X Large	1494-TB1
XX Large	1495-TB1







sigvaris

Custom Garments

Sigvaris offers a line of custom compression garments for patients who measure outside of the standard size range of the ready-to-wear garments or those with a complex limb shape. Custom garments may be the best solution for patients with advanced or complex lymphedema for patients with Stage III Lymphedema.

Manufactured with high grade textiles, our compression garments are the industry standard for custom lymphedema compression wear for day or night. Each custom garment is manufactured in our facility in Holland, Michigan.

To order Sigvaris custom garments, request a custom ordering packet by emailing

sigvariscares@sigvaris.com.

Chipsleeve Custom ARM

Compression Level

15-25mmHg



Color



Black

Features

The Chipsleeve Custom Arm provides mild to moderate gradient compression from the hand to the axilla.

The Chipsleeve Custom Arm features chip foam sewn in directional channels to stimulate tissue and encourage the movement of lymph fluid. The garment is made using a moisture-wicking bacteriostatic material to keep the limb dry to encourage good hygiene. Patients who measure outside of the Chipsleeve Standard Arm size range may select the Chipsleeve Custom Arm.

To review the standard product order form, please see page 44.

Ideal for

- Lymphedema
- Fibrotic tissue

Includes

- Chipsleeve Arm
- Oversleeve (Black)
- 1 Pair of non-compressive Cotton Liners

Product Selection Tip

May be right for patients requiring a nighttime garment for the management of lymphedema



Custom

Chipsleeve Custom ARM

Patient Name:		Contact Name:
Account Name:		Account #:
Account Phone #:		P.O. #:
Item #: 2739-AR	Height:	Weight:
Ship Name & Address:		Date:


Product Information

Product includes one Chipsleeve Custom Arm, one pair of Cotton Liners, and one black Oversleeve.

- | | |
|------------------------------------|--------------|
| <input type="checkbox"/> Right Arm | Color: Black |
| <input type="checkbox"/> Left Arm | Color: Black |

Circumference

	Left	Right
C	_____	_____
B1	_____	_____
B	_____	_____
A1	_____	_____
A	_____	_____
D Length		
<i>Measure length along dorsal aspect</i>		
D	_____	_____

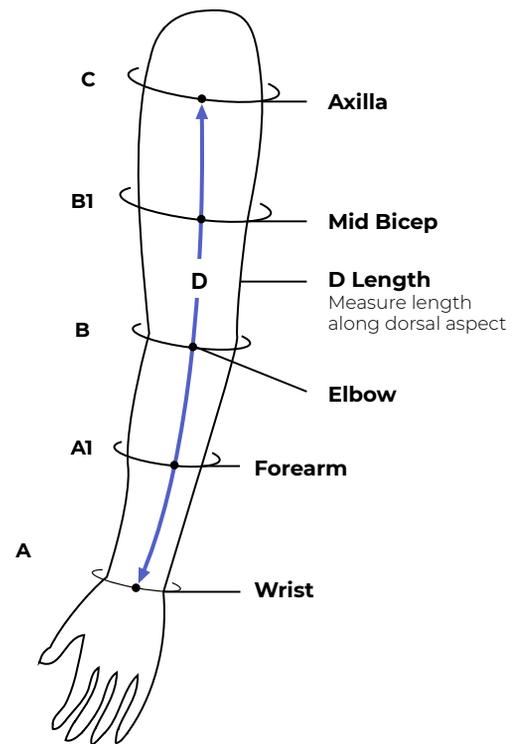

Important

Exact measurements are critical for this garment to ensure a proper fit. If you'd like to learn more about measuring and fitting, attend our Medically Complex Edema Basic Fitter Training. To sign up, visit sigvariseducation.com or contact your local territory manager for more information.

Alternatively, call Customer Care Solution Center at **800-322-7744**, or e-mail us_orders@sigvaris.com, to receive a remote consultation/training.

Supplies Needed

- Cell phone with camera. Photos of the arm with measurement markings must be emailed to: us_orders@sigvaris.com
- Measuring instructions and forms.
- SIGVARIS GROUP Measuring tape and body pen (or eyeliner pencil).
- Signed Custom Order Terms & Conditions Form (include with order).

MEASURING INSTRUCTIONS

With patient seated, place the arm extended and elbow bent slightly, with the palm down, on a flat surface.

CIRCUMFERENCES

- Measure circumference at palm and record on line E.
- Measure circumference at wrist, and record on line A.
- Mark dorsal aspect at distal edge of tape.
- Measure circumference at forearm, and record on line A1.
- Measure circumference at elbow, and record on line B.
- Measure circumference at mid bicep, and record on line B1.
- Measure circumference at axilla, and record on line C.
- Mark dorsal aspect at proximal edge of tape.

LENGTH

Measure length at dorsal aspect from mark at Point A to mark at Point C and record in box D.



Chipvest

Compression Level

Variable compression



Color



Black

Features

The Chipvest is a custom upper body compression garment that offers mild to moderate compression for individuals with lymphedema.

The product contains open-cell foam chips that help to soften fibrotic areas. The chip foam is formed into channels to help direct lymph fluid and facilitate drainage. To improve fit and comfort, the garment includes a back drawstring and corset design as well as hook and loop shoulder straps that are easily adjusted by the patient. Clinicians may further customize the product by selecting the location of the foam (unilateral or bilateral).

Ideal for


- Lymphedema
- Post-surgical or non-specific edema

Includes

- Chipvest

Product Selection Tip

Select this product for lymphedema and post-mastectomy patients

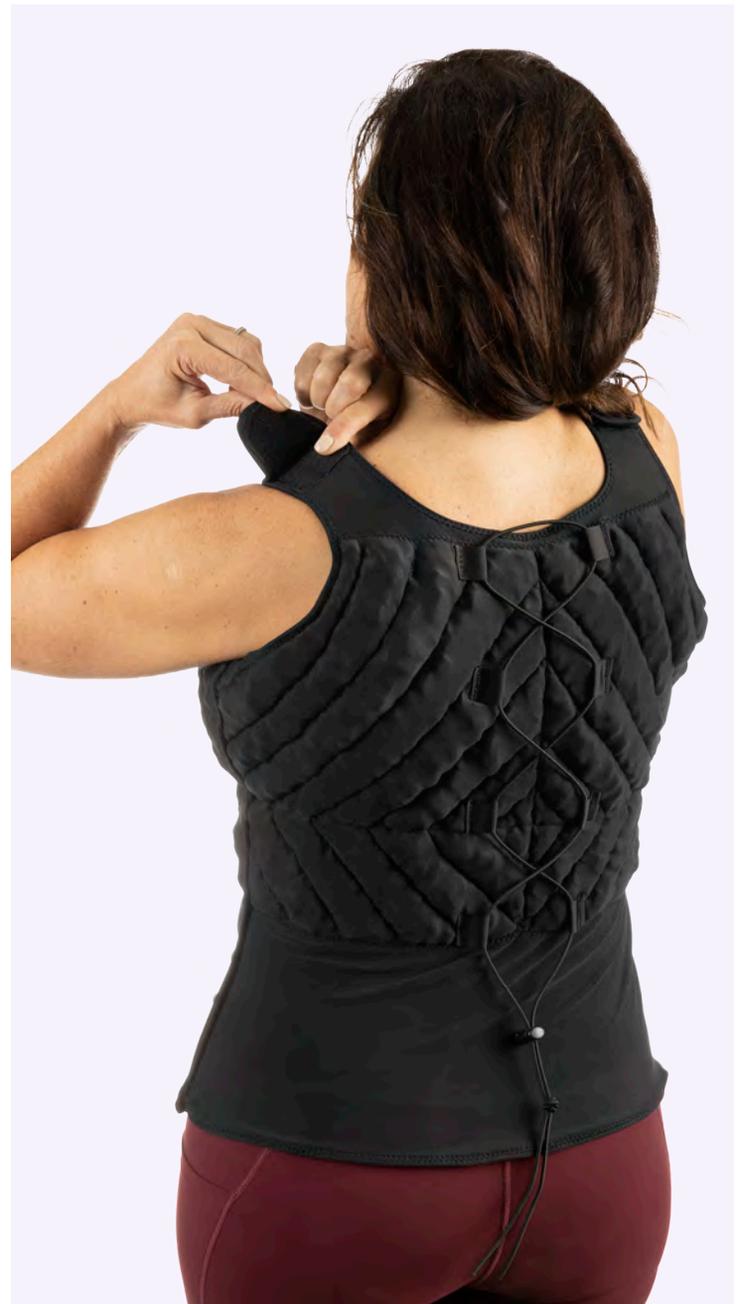

Custom

Chipvest

Patient Name:	Contact Name:
Account Name:	Account #:
Account Phone #:	P.O. #:
Height:	Weight:
Ship Name & Address:	Date:

Product Information

Product includes one Chipvest Custom.

Custom Chipvest Item #: **2239-VS**

- Full (Bilateral)
- Right Side (Unilateral)
- Left Side (Unilateral)

Measuring Instructions

Step One:

Measure the length in centimeters from the suprasternal notch to the waist.


Step Two:

Measure circumferences in centimeters.

Custom Size Measurements

Length	Centimeters
L	
Circumference	Centimeters
Axilla	
Largest Chest	
Xyphoid Process	
Waist	
Mid-Hip	

Caution: Drawstring contains natural rubber latex.

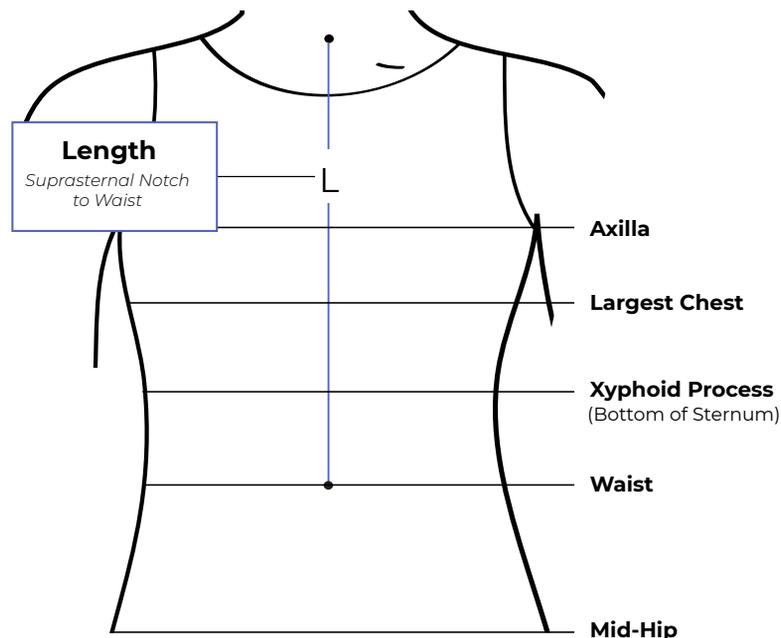

Important

Exact measurements are critical for this garment to ensure a proper fit. If you'd like to learn more about measuring and fitting, attend our Medically Complex Edema Basic Fitter Training. To sign up, visit sigvariseducation.com or contact your local territory manager for more information.

Alternatively, call Customer Care Solution Center at **800-322-7744**, or e-mail us_orders@sigvaris.com, to receive a remote consultation/training.

Supplies Needed

- Cell phone with camera. Photos of the torso with measurement markings must be emailed to: us_orders@sigvaris.com
- Measuring instructions and forms.
- SIGVARIS GROUP Measuring tape and body pen (or eyeliner pencil).
- Signed Custom Order Terms & Conditions Form (include with order).

us_orders@sigvaris.com
P: 800-322-7744 · F: 800-481-5488



For a PDF of the product ordering form, scan the QR code or visit sigvaris.group/mceforms

Compreshort

CUSTOM

Compression Level

10-15mmHg



Color



Black

Features

The Compreshort offers mild, supportive compression for the treatment of abdominal or truncal edema. Available in short and capri lengths, patients with lymphedema and lipedema benefit from the therapeutic compression and support the garment provides.

To review the standard product order form, please see page 56.

Ideal for


- Lymphedema
- Lipedema

Includes

- Compreshort

Product Selection Tip

May be helpful for patients with lipedema, lymphedema, or lipolymphedema in the trunk and thigh

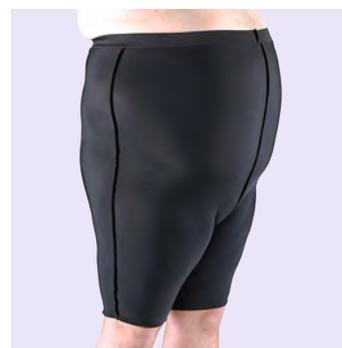

Compreshort

Product available in **short, capri, and custom lengths.**

Custom

Compresshort CUSTOM

Patient Name:	Contact Name:
Account Name:	Account #:
Account Phone #:	P.O. #:
Height:	Weight:
Ship Name & Address:	Date:

Product Information

Custom Compresshorts are ideal for patients outside of the standard size range. Product includes one Compresshort Custom.

Short Length

Color: Black

Quantity:

Capri Length

Color: Black

Quantity:

Caution: Elastic Band in sizes 2XL and larger contains natural rubber latex.
Note: Ready to wear sizes are available in short and capri length.


Circumference

Waist _____

Hip _____

Midthigh _____

Calf _____

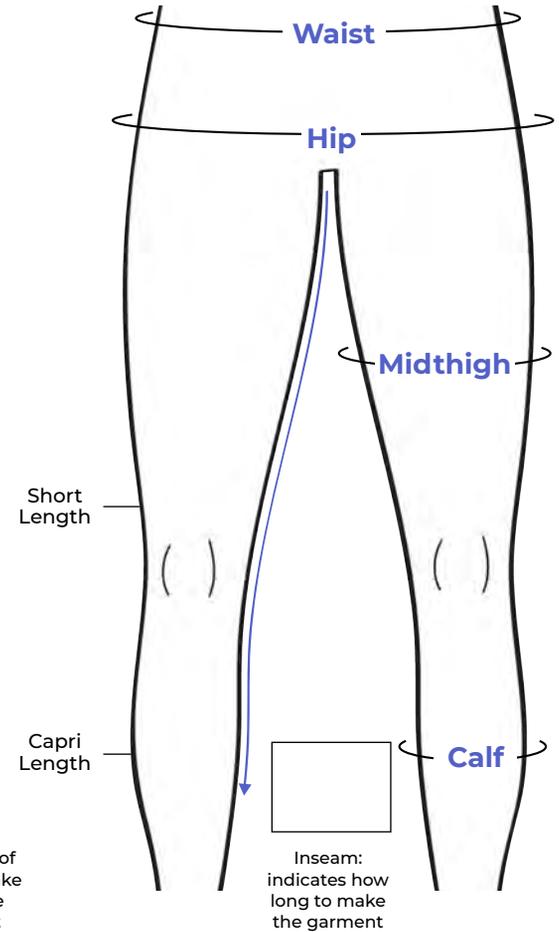

Measuring Instructions

Step One:

Record the waist, hip, and midthigh measurements

Step Two:

Measure from the groin down and indicate the garment inseam length. Record the inseam length

Step Three:

Mark the end of the inseam and take a circumference measurement at that point. Record the inseam circumference

Step Four:

After recording these measurements, please call Customer Care Solution Center: **800-322-7744**

Important

Exact measurements are critical for this garment to ensure a proper fit. If you'd like to learn more about measuring and fitting, attend our Medically Complex Edema Basic Fitter Training. To sign up, visit sigvariseducation.com or contact your local territory manager for more information.

Alternatively, call Customer Care Solution Center at **800-322-7744**, or e-mail us_orders@sigvaris.com, to receive a remote consultation/training.

Supplies Needed

- Cell phone with camera. Photos of relevant limbs with measurement markings must be emailed to: us_orders@sigvaris.com
- Measuring instructions and forms.
- SIGVARIS GROUP Measuring tape and body pen (or eyeliner pencil).
- Signed Custom Order Terms & Conditions Form (include with order).

us_orders@sigvaris.com

P: 800-322-7744 · F: 800-481-5488



For a PDF of the product ordering form, scan the QR code or visit

sigvaris.group/mceforms

Chipflow

Ideal for patients in need of an easy-to-don garment that pairs the massaging effect of foam with an inelastic wrap for a high level of containment.

Compression Level

Variable compression



Color



Black

Gradient compression, high containment

The Chipflow pairs an inelastic wrap with a chip foam liner to form a garment that can be worn day or night. This garment is customized to fit individuals with a small ankle and large calf, commonly seen in patients with lymphatic disorders.

Each Chipflow garment contains a channeled foam liner which creates a pressure gradient that promotes lymph fluid movement and healthy skin. It also contains a unique split-ankle zone that allows for independent application while maintaining compression. The inelastic wrap is applied over the liner to allow for proper distribution of compression and a high level of containment.

Each Chipflow garment includes the Compreboot with gripper patches to allow for foot compression without slippage. Additionally, each garment includes non-compressive liners for protection of the skin.

This product is patent pending.



Ideal for

- Lymphedema
- Post surgical or non-specified edema
- Lipedema

Custom Garments

LOWER EXTREMITY



Chipflow

CALF & FOOT

Compression Level

Variable compression



Color



Black

Features

The Chipflow Calf & Foot is an easy to apply garment made with comfortable, low profile, chip foam and an inelastic wrap. This combination provides gradient compression and high containment for the calf and foot.

Ideal for

- Lymphedema
- Lipedema
- Post surgical or non-specified edema

Includes

- 1 Chipflow Calf
- 1 Compreboot Foot
- 1 Chip Foam Liner
- 1 Pair of non-compressive Basic Liners

Product Selection Tip

May be right for patients who present Lymphedema in the calf and foot that requires compression paired with massaging mechanism of chip foam.



Custom

Chipflow CALF & FOOT

Patient Name:	Contact Name:
Account Name:	Account #:
Account Phone #:	P.O. #:
Ship Name & Address:	Date:

Product Information

Product includes one custom Chipflow Calf, one Compreboot Standard Foot wrap, one Chip Foam Liner, and one pair of non-compressive Basic Liners*.

<input type="checkbox"/> Right Calf	Color: Black	Quantity:
<input type="checkbox"/> Right Foot	Color: Black:	Quantity:
<input type="checkbox"/> Left Calf	Color: Black	Quantity:
<input type="checkbox"/> Left Foot	Color: Black:	Quantity:

*For the largest size, a pair of non-compressive Cotton Liners is included instead of the Basic Liners.

Important

Exact measurements are critical for this garment to ensure a proper fit. If you'd like to learn more about measuring and fitting, attend our Medically Complex Edema Basic Fitter Training. To sign up, visit sigvariseducation.com or contact your local territory manager for more information.


Alternatively, call Customer Care Solution Center at **800-322-7744**, or e-mail us_orders@sigvaris.com, to receive a remote consultation/training.

Supplies Needed

- Measuring instructions and forms.
- SIGVARIS GROUP Measuring tape and body pen (or eyeliner pencil).

Circumference

	Left	Right
C1	_____	_____
C	_____	_____
B	_____	_____
A	_____	_____
G	_____	_____
I	_____	_____
J	_____	_____
K	_____	_____

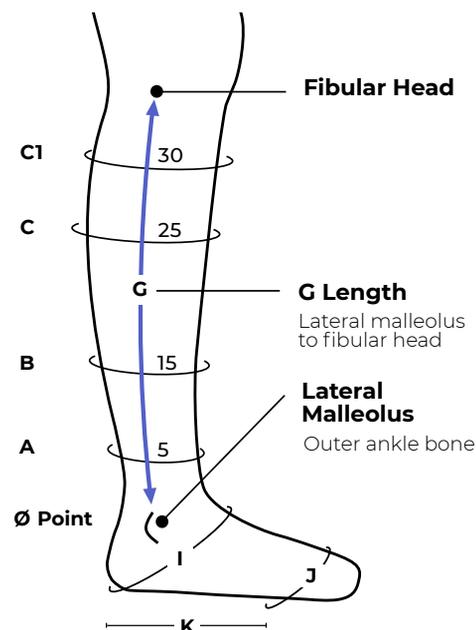

K: Measure **medial** length from heel to 1st metatarsal head

Chipflow Calf Size Chart (A-C1) (Available in Black Only)

Note: Product contains Measurement Limitations with Min & Max sizing for ankle and calf

	Min	Max
C1	31cm	69cm
C	31cm	69cm
B	26cm	59cm
A	19cm	51cm
Black		
G Length	27cm	42cm

Measurement Disclaimer

Max calf < 2x ankle



Chipflow

FULL LEG

Compression Level

Variable compression



Color



Black

Features

The Chipflow Full Leg is an easy to apply garment made with comfortable, low profile, chip foam and an inelastic wrap. This combination provides gradient compression and high containment throughout the full leg.

Ideal for

- Lymphedema
- Lipedema
- Post surgical or non-specified edema

Includes

- 1 Chipflow Full Leg
- 1 Compreboot Foot
- 1 Chip Foam Liner
- 1 Pair of non-compressive Cotton Liners

Product Selection Tip

May be right for patients who present Lymphedema in the full leg and require compression paired with massaging mechanism of foam.



Patient Name:	Contact Name:
Account Name:	Account #:
Account Phone #:	P.O. #:
Ship Name & Address:	Date:

Product Information

Product includes one custom Chipflow Full leg, one Compreboot Standard Foot wrap, one Chip Foam Liner, and one pair of non-compressive Cotton Liners.

- | | | |
|-------------------------------------|---------------|-----------|
| <input type="checkbox"/> Right Leg | Color: Black | Quantity: |
| <input type="checkbox"/> Right Foot | Color: Black: | Quantity: |
| <input type="checkbox"/> Left Leg | Color: Black | Quantity: |
| <input type="checkbox"/> Left Foot | Color: Black: | Quantity: |

Important

Exact measurements are critical for this garment to ensure a proper fit. If you'd like to learn more about measuring and fitting, attend our Medically Complex Edema Basic Fitter Training. To sign up, visit sigvariseducation.com or contact your local territory manager for more information.


Alternatively, call Customer Care Solution Center at **800-322-7744**, or e-mail us_orders@sigvaris.com, to receive a remote consultation/training.

Supplies Needed

- Measuring instructions and forms.
- SIGVARIS GROUP Measuring tape and body pen (or eyeliner pencil).

Circumference

	Left	Right
F	_____	_____
E	_____	_____
D	_____	_____
C1	_____	_____
C	_____	_____
B	_____	_____
A	_____	_____
H	_____	_____
I	_____	_____
J	_____	_____
K	_____	_____

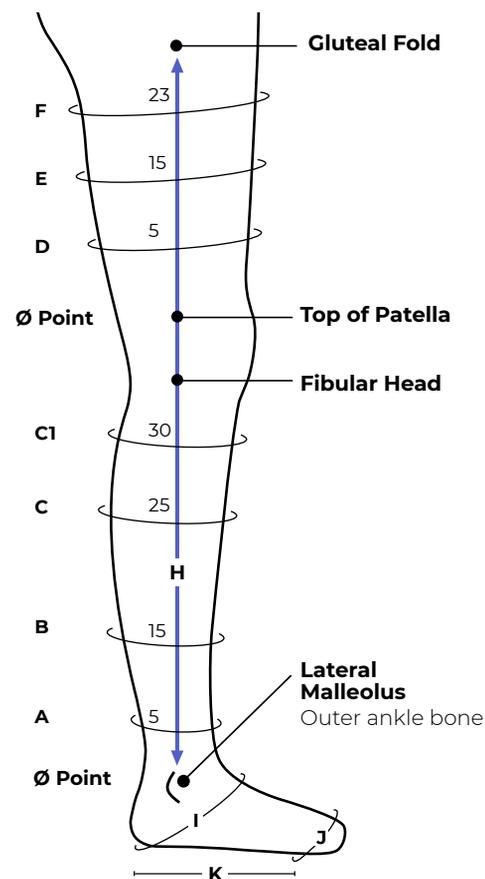

K: Measure **medial** length from heel to 1st metatarsal head

Chipflow Full Leg Size Chart (A-C1) (Available in Black Only)

Note: Product contains Measurement Limitations with Min & Max sizing for ankle, calf, and thigh

	Min	Max
F	48cm	90cm
E	43cm	84cm
D	38cm	81cm
C1	31cm	69cm
C	31cm	69cm
B	26cm	59cm
A	19cm	51cm
Black		
H Length	66cm	91cm

Measurement Disclaimer

Max thigh < 2x calf	Max calf < 2x ankle
---------------------	---------------------



Armassist

Ideal for patients who need a custom compression garment with high containment.

Compression Level

30–40mmHg



Color



Black & Green

Firm, durable compression with high containment

The Armassist offers 30–40mmHg of firm, gradient compression. Ideal for patients with moderate to severe lymphedema, the Armassist is constructed with a durable, high grade fabric and includes adjustable hook and loop straps. The adjustable straps may be positioned to properly shape and contain the limb. The strap positioning also ensures the gradient application of compression.

To support increased ease of application of the product, a zipper along the forearm is included. The zipper minimizes the need for arm strength and dexterity and encourages independent wear and removal. The Armassist also includes a hand piece made with palm and dorsal foam for hand and wrist compression. To increase compression at a specific location on the limb, apply the included Unistrap.

Clinicians may customize this garment to their patient's needs by selecting the preferred foam.



Ideal for

- Lymphedema

Includes

- Armassist
- 1 Unistrap
- 1 Pair of non-compressive Cotton Liners

Product Selection Tip

May be right for patients with moderate to severe lymphedema

Measuring and Ordering

CUSTOM GARMENTS

Measuring instructions and order forms for custom Armassist compression wear are found in the following pages, and can also be found online at the link on those pages.

Our Customer Care team is available to assist you through the process of ordering a custom Armassist. From determining if your patient is a good candidate for a custom garment to providing measuring guidance, we are here to help!

To assist in making the custom garment, we require photos of the patient as the measurements are taken, and the measurement marks must be visible in the photos. Please email the photos to us_orders@sigvaris.com.

All custom garments in the catalog are patterned and custom made according to the provided measurements at our manufacturing facility in Holland, Michigan.

Contact us:

1.800.322.7744, ext. 1

us_orders@sigvaris.com



Durability for long-term wear and comfort

Durable high grade textiles give the product unsurpassed quality and ensures long-term wear

Increased Adjustability

Straps extended for increased adjustability and improved application of compression

Seamless Compression

Palm and dorsum foam offer compression to the hand with the same garment

CUSTOM Measuring Instructions**Important**

Exact measurements are critical for this garment to ensure a proper fit. If you'd like to learn more about measuring and fitting, attend our Medically Complex Edema Basic Fitter Training. To sign up, visit sigvariseducation.com or contact your local territory manager for more information.

Alternatively, call Customer Care Solution Center at **800-322-7744**, or e-mail us_orders@sigvaris.com, to receive a remote consultation/training.

Supplies Needed

- Cell phone with camera.
- Photos of the arm with measurement markings must be emailed to: us_orders@sigvaris.com
- Measuring instructions and forms.
- SIGVARIS GROUP Measuring tape and body pen (or eyeliner pencil).
- Signed Custom Order Terms & Conditions Form (include with order).

INSTRUCTIONS

All measurements should be recorded in centimeters. Apply slight tension to hold the tape measure in place.

STEP 1: Locate Landmarks

1. Encircle arm with tape measure located as high as possible at axilla. Mark on the proximal edge of tape at the anterior, dorsal, and posterior aspects of arm, and record the circumference on the line labeled posterior axilla.
2. Encircle elbow with the distal edge of a tape measure aligned at the elbow crease and olecranon process (elbow). Mark the dorsal aspect of the elbow along this edge, and record the circumference on the line labeled elbow (Ø point).
3. Encircle wrist with the distal edge of a tape measure aligned with the ulnar styloid process. Mark the dorsal aspect of the wrist along this edge, and record the circumference on the line labeled ulnar styloid.

STEP 2: Locate Circumference Measurement Positions

1. Starting with "0" of the tape measure at your mark on the dorsal aspect of the elbow, mark the dorsal aspect of the upper arm in 5cm increments up to the mark made from the axilla on the dorsal aspect of the arm.
2. Starting with the "0" of the tape measure at your mark on the dorsal aspect of the elbow, mark the dorsal aspect of the forearm in 5cm increments down to the mark made for the ulnar styloid on the dorsal aspect of the arm.

STEP 3: Measure Lengths

1. Measure upper arm from anterior axilla (inside armpit) to the elbow crease and record length on box A.
2. Measure upper arm from posterior axilla (outside armpit) to the olecranon process (elbow) and record length on box B.
3. Measure the forearm along dorsal aspect from elbow (as marked at the mid-point between olecranon process and elbow crease) to the ulnar styloid as indicated on diagram and record length on box C.
4. Measure hand from your mark on dorsal aspect of wrist to the third metacarpal head and record length on box D.
5. Measure width of hand across dorsal metacarpal heads and record width on box E.

STEP 4: Measure Circumferences

1. Encircle arm with a tape measure at each mark located between the elbow and axilla, and record each circumference on the line corresponding to each position.
2. Encircle arm with a tape measure at each mark located between the elbow and ulnar styloid, and record each circumference on the line corresponding to each position.

TAKE DIGITAL PHOTOS (REQUIRED) AFTER MEASUREMENT MARKINGS ARE MADE.

EMAIL TO: us_orders@sigvaris.com

NOTE: Armassist is available in: 7099-AR FlatFoam, 7199-AR Wavefoam, 7299-AR Super FlatFoam, 7399-AR Super WaveFoam.

Patient Name: _____ Contact Name: _____ Height: _____ Weight: _____ Date: _____

Ship Name & Address: _____ Account Name: _____ Account #: _____ P.O #: _____


Product Information

Product includes one Armassist and one pair of Cotton Liners.

Right Arm Left Arm Foam
 Flat
 Wavefoam

Circumference Measurements

● = Locations measured along **Dorsal** aspect

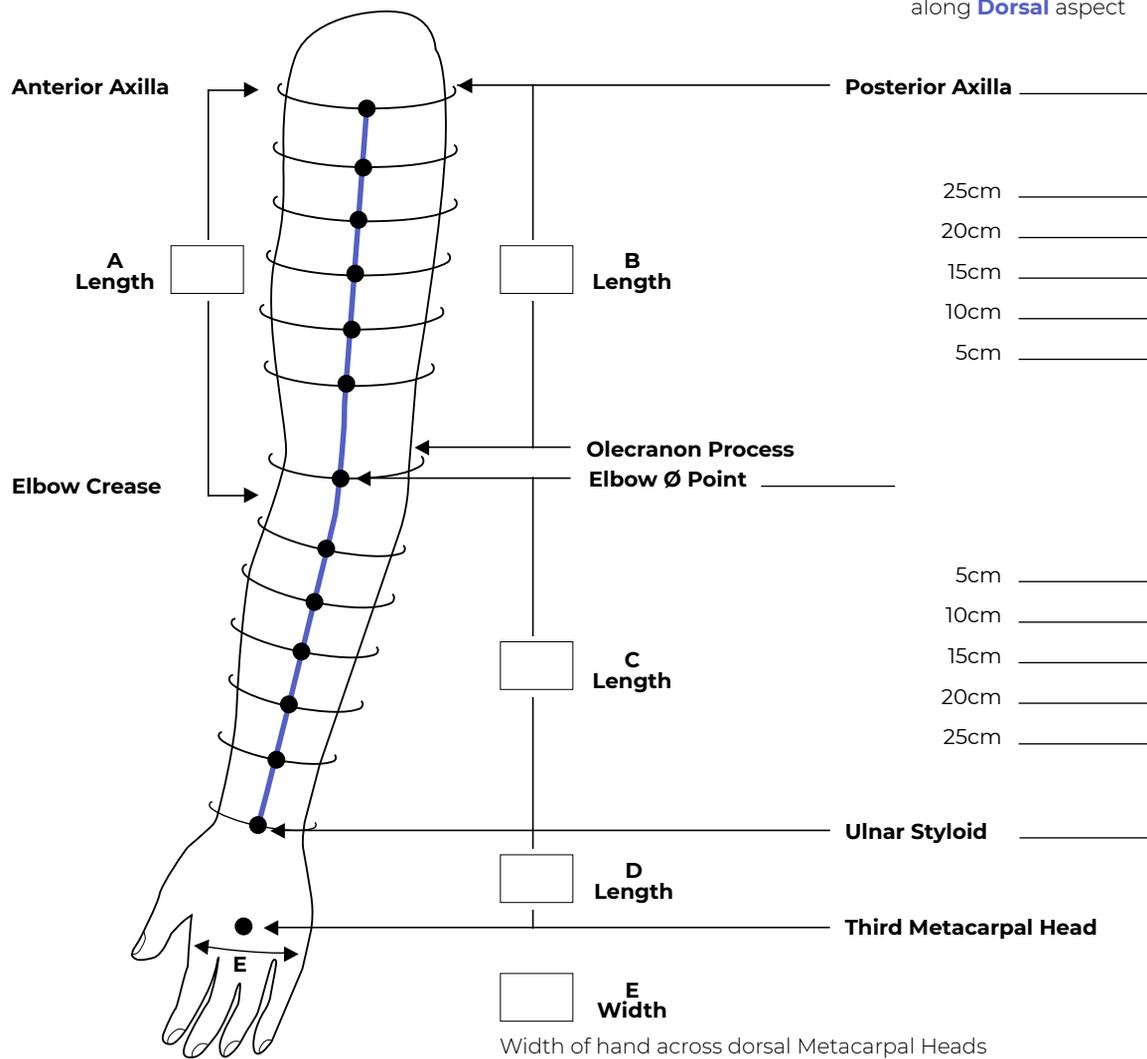

Legassist

Ideal for patients who need a custom compression garment with high containment.

Compression Level

30–40mmHg



Color



Black & Green

Firm, durable compression with high containment

Sigvaris Legassist garments are ideal for patients with moderate to severe lymphedema. Manufactured with medical grade textiles, the Legassist line of custom compression garments are the industry standard for custom lymphedema compression wear for day or night. Constructed with a non-stretch medical grade fabric that is non-static and non-absorbent, the Legassist garments are manufactured with precision and a high standard of quality to the custom specifications provided by the clinician.

The Legassist product is made with Flat Foam or Wavefoam (selected by clinician and patient) which encourage the movement and drainage of lymph fluid and contains a zipper panel for easy application and removal. Clinicians and patients may further customize the garment by adding strapping over the knee or by selecting the optional Belt and Hip Attachment.

The Legassist line of products is clinically proven to be effective for the containment and reduction of edema and have helped thousands of lymphedema patients worldwide experience an improved quality of life.



Ideal for

- Lymphedema

Product Selection Tip

Select this custom-sized garment for patients with moderate to severe lymphedema

Custom Garments

LOWER EXTREMITY



Medical Grade Fabric

Durable fabric is non-stretch, non-static, and non-absorbent

Superior Product Construction

Product is manufactured using premium materials and intelligent processes

sigvaris

Legassist

CALF & FOOT

Compression Level

30-40mmHg



Color



Black & Green

Includes

- Legassist Calf and Foot
- Compreboot Plus Foot
- 1 Pair of non-compressive Cotton Liners

Features

The Legassist Calf and Foot is a low profile and lightweight custom garment offering gradient compression from ankle to knee bend.

See page 114 for measuring chart and instructions



Legassist

CALF & FOOT WITH SHELF

Compression Level

30-40mmHg



Color



Black & Green

Includes

- Legassist Calf and Foot
- Compreboot Plus Foot
- 1 Pair of non-compressive Cotton Liners

Features

The Legassist Calf and Foot is a low profile and lightweight custom garment offering gradient compression from ankle to knee bend.

See page 116 for measuring chart and instructions

Legassist Lobule Compression System

LCS

Compression Level

30-40mmHg



Color



Black & Green

Includes

- Legassist Lobule Compression System
- 1 Pair of non-compressive Cotton Liners

Features

The Legassist LCS is a custom garment that lifts enlarged areas. The garment may be paired with the Calf and Foot or Belt and Hip Attachment for additional support.

See page 120 for measuring chart and instructions



Legassist

FULL LEG

Compression Level

30-40mmHg



Color



Black & Green

Includes

- Legassist Full Leg
- Compreboot Plus Foot
- 1 Pair of non-compressive Cotton Liners

Features

The Legassist Full Leg is a low profile and lightweight custom garment offering gradient compression from the foot to gluteal fold.

See page 124 for measuring chart and instructions



Measuring and Ordering

CUSTOM GARMENTS

Measuring instructions and order forms for custom Legassist Calf & Foot compression wear are found in the following pages, and can also be found online at the link below.

Our Customer Care team is available to assist you through the process of ordering a custom Legassist Calf & Foot. From determining if your patient is a good candidate for a custom garment to providing measuring guidance, we are here to help!

To assist in making the custom garment, we require photos of the patient as the measurements are taken, and the measurement marks must be visible in the photos. Please email the photos to us_orders@sigvaris.com.

All custom garments in the catalog are patterned and custom made according to the provided measurements at our manufacturing facility in Holland, Michigan.

Contact us:

1.800.322.7744, ext. 1

us_orders@sigvaris.com



For a PDF of the product ordering form, scan the QR code or visit

sigvaris.group/mceforms


Patient Name:	Contact Name:	Height:	Weight:	Date:
Ship Name & Address:	Account Name:	Account #:	P.O #: :	

Product Information

Product includes one Legassist Calf, one Compreboot Plus Foot, and one pair of Cotton Liners.

<input type="checkbox"/> Right Leg Size: <input type="checkbox"/> Regular <input type="checkbox"/> Super	<input type="checkbox"/> Left Leg Size: <input type="checkbox"/> Regular <input type="checkbox"/> Super	Foam <input type="checkbox"/> Flat <input type="checkbox"/> Wavefoam	Additional Foot Option <input type="checkbox"/> Custom Medaboot (additional charge)
<input type="checkbox"/> Right Leg Size: <input type="checkbox"/> Regular <input type="checkbox"/> Super	<input type="checkbox"/> Left Leg Size: <input type="checkbox"/> Regular <input type="checkbox"/> Super	Foam <input type="checkbox"/> Flat <input type="checkbox"/> Wavefoam	Additional Foot Option <input type="checkbox"/> Custom Medaboot (additional charge)

Note: If the greatest circumference measurement is >60cm, order a Super.

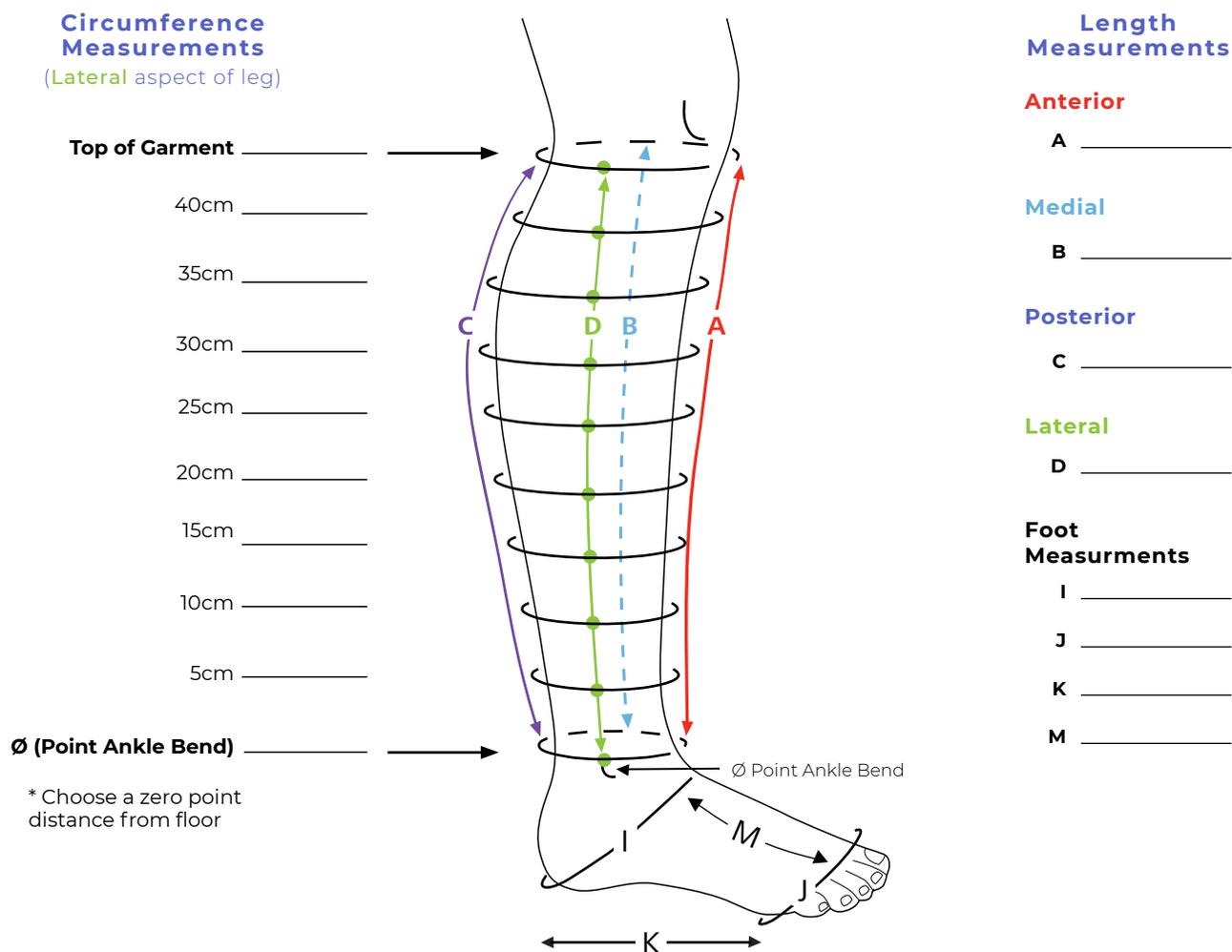

CUSTOM Measuring Instructions

Important

Exact measurements are critical for this garment to ensure a proper fit. If you'd like to learn more about measuring and fitting, attend our Medically Complex Edema Basic Fitter Training. To sign up, visit sigvariseducation.com or contact your local territory manager for more information.

Alternatively, call Customer Care Solution Center at **800-322-7744**, or e-mail us_orders@sigvaris.com, to receive a remote consultation/training.

Supplies Needed

- Cell phone with camera. Photos of the leg with measurement markings must be emailed to: us_orders@sigvaris.com
- Measuring instructions and forms.
- SIGVARIS GROUP Measuring tape and body pen (or eyeliner pencil).
- Signed Custom Order Terms & Conditions Form (include with order).

INSTRUCTIONS

All measurements should be recorded in centimeters. Apply slight tension to hold the tape measure in place.

STEP 1: Locate Landmarks

1. The ideal placement for the top edge of the garment is two finger widths below the popliteal crease, and across the fibular head and tibial tuberosity. However, the position of the top of garment can be customized, if needed, to overlap better with a Legassist LCS thigh piece. Encircle leg with top edge of tape measure at the desired location. Using the SIGVARIS GROUP body pen, trace a line around the leg along the top edge of the tape and measure and record the circumference on top of garment line.
2. Encircle leg with bottom edge of tape measure aligned with the ankle bend. Mark leg on this edge of tape at the anterior, medial, posterior and lateral aspects, and record the circumference on the line labeled ankle bend (ø point).

STEP 2: Locate Circumference Measurement Positions

1. Starting with "0" of the tape measure at your mark on the lateral aspect of the ankle bend, mark the leg in 5cm increments up to your mark on the lateral aspect of the top of garment.

STEP 3: Measure Lengths

1. Measure leg from your mark on the anterior aspect of the ankle bend to the anterior aspect of top of garment and record length on line marked A.
2. Measure leg from your mark on the medial aspect of the ankle bend to the medial aspect of top of garment and record length on line marked B.
3. Measure leg from your mark on the posterior aspect of the ankle bend to the posterior aspect of top of garment and record length on line marked C.
4. Measure leg from your mark on the lateral aspect of the ankle bend to the lateral aspect of top of garment and record length on line marked D.

STEP 4: Measure Circumferences

1. Encircle leg with a tape measure at each mark located between the ankle bend and top of garment, and record each circumference on the line corresponding to each position.

STEP 5: Measure Foot

1. Measure foot from heel to 1st metatarsal head (or desired boot length) and record length on line K.
2. Measure top of foot from the 3rd metatarsal head to ankle bend and record length on line M.
3. Encircle the instep of foot and apex of heel with a tape measure and record the circumference on line I.
4. Encircle the foot across the metatarsal heads and record the circumference on Line J.

TAKE DIGITAL PHOTOS (REQUIRED) AFTER MEASUREMENT MARKINGS ARE MADE. EMAIL TO: us_orders@sigvaris.com

NOTE: Legassist Calf & Foot is available in:

7099-BK Flat Foam, 7199-BK WaveFoam, 7299-BK Super, Flat Foam, 7399-BK Super, WaveFoam.

us_orders@sigvaris.com

P: 800-322-7744 · F: 800-481-5488

Measuring and Ordering

CUSTOM GARMENTS

Measuring instructions and order forms for custom Legassist Calf & Foot with Shelf compression wear are found in the following pages, and can also be found online at the link below.

Our Customer Care team is available to assist you through the process of ordering a custom Legassist Calf & Foot with Shelf. From determining if your patient is a good candidate for a custom garment to providing measuring guidance, we are here to help!

To assist in making the custom garment, we require photos of the patient as the measurements are taken, and the measurement marks must be visible in the photos. Please email the photos to us_orders@sigvaris.com.

All custom garments in the catalog are patterned and custom made according to the provided measurements at our manufacturing facility in Holland, Michigan.

Contact us:
1.800.322.7744, ext. 1
us_orders@sigvaris.com



For a PDF of the product ordering form, scan the QR code or visit sigvaris.group/mceforms

Patient Name: _____ Contact Name: _____ Height: _____ Weight: _____ Date: _____


Ship Name & Address: _____ Account Name: _____ Account #: _____ P.O #: _____

Product Information

Product includes one Legassist Calf, one Compreboot Plus Foot, and one pair of Cotton Liners.

Right Leg Foam Additional Foot Option
 Left Leg Flat Wavefoam Custom Medaboot (additional charge)

Note: If the greatest circumference measurement is >60cm, order a Super.

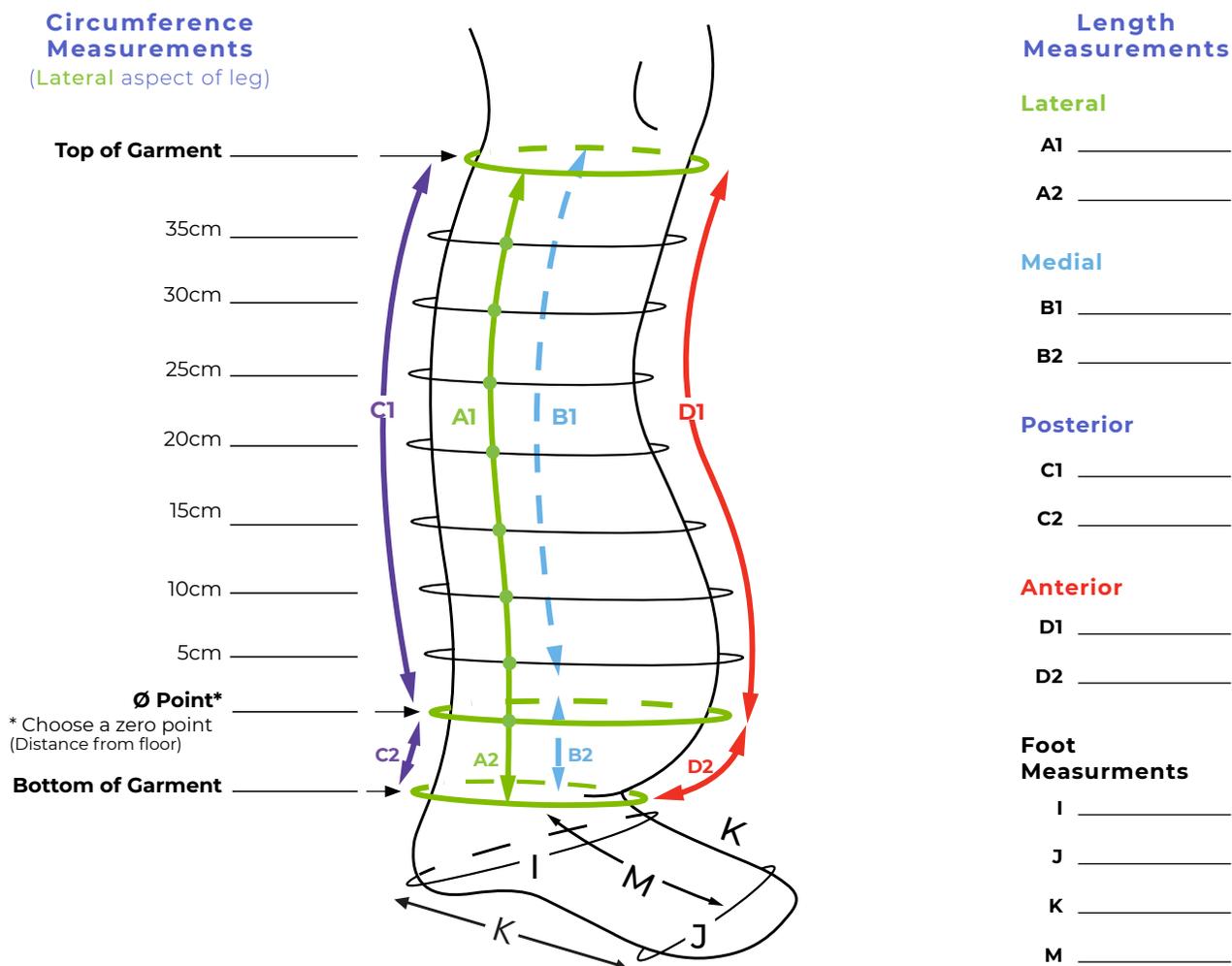

Scan the QR code to view a Measuring Instruction video



For a PDF of the product ordering form, scan the QR code or visit sigvaris.group/mceforms

Important

Exact measurements are critical for this garment to ensure a proper fit. If you'd like to learn more about measuring and fitting, attend our Medically Complex Edema Basic Fitter Training. To sign up, visit sigvariseducation.com or contact your local territory manager for more information.

Alternatively, call Customer Care Solution Center at **800-322-7744**, or e-mail us_orders@sigvaris.com, to receive a remote consultation/training.

Supplies Needed

- Cell phone with camera. Photos of the leg with measurement markings must be emailed to: **us_orders@sigvaris.com**
- Measuring instructions and forms.
- SIGVARIS GROUP Measuring tape and body pen (or eyeliner pencil).
- Signed Custom Order Terms & Conditions Form (include with order).

INSTRUCTIONS

All measurements should be recorded in centimeters. Apply slight tension to hold the tape measure in place.

STEP 1: Choose The Zero Point (Ø point)

The zero point is unique to the Legassist Calf & Foot with Shelf garment. It is the reference point from which all other measurements are taken, in order to determine the lengths needed for the custom garment to fit the patient's individual ankle shape. Encircle leg with bottom edge of tape measure immediately above the shelf, at the lowest aspect of the curve of the ankle lobule where a tape measure can be applied without sliding off. Trace a line around the leg along the bottom edge of the tape using the SIGVARIS GROUP body pen, and record the circumference on Ø point line. Document your zero point location by measuring the distance from the floor.

STEP 2: Choose Where The Garment Will End Distally - "bottom of garment"

Encircle leg with bottom edge of tape measure where the ankle lobule touches the ankle and/or foot. Do not pull tape under shelf. Mark leg on this edge of tape at the anterior, medial, posterior and lateral aspects, and record the circumference on the line labeled bottom of garment. The smallest ankle measurement is not relevant for this product. The product is designed to hug and lift the ankle lobule.

STEP 3: Choose Where The Garment Will End Proximally - "top of garment"

The ideal placement for the top edge of the garment is two finger widths below the popliteal crease, and across the fibular head and tibial tuberosity. However, the position of the top of garment can be customized, if needed, to overlap better with a Legassist LCS thigh piece. Encircle leg with top edge of tape measure at the desired location. Using the SIGVARIS GROUP body pen, trace a line around the leg along the top edge of the tape and measure and record the circumference on top of garment line.

STEP 4: Measure And Record Lengths

Measure each of the lengths below from the Zero Point, following the contour of the leg and ankle, but not extending the tape measure into folds. Lengths measured distal to the Zero Point should end where the ankle lobule touches the ankle or foot at the bottom of garment versus extending the tape measure up under the lobule. Record lengths as follows:

Lateral (A) Lengths:

A1: Measure & record lateral length from Ø point to top of garment

A2: Measure & record lateral length from Ø point to bottom of garment

Medial (B) Lengths

B1: Measure & record medial length from Ø point to top of garment

B2: Measure & record medial length from Ø point to bottom of garment

Posterior (C) Lengths

C1: Measure & record posterior length from Ø point to top of garment

C2: Measure & record posterior length from Ø point to bottom of garment

Anterior (D) Lengths

D1: Measure & record anterior length from Ø point to top of garment

D2: Measure & record anterior length from Ø point to bottom of garment

STEP 5: Measure And Record Circumferences

On the lateral aspect of the leg only, mark every 5cm moving proximally from the zero point. Next, record the circumferences at each point, making sure that the measuring tape is parallel to the floor.

STEP 6: Measure Foot

1. Measure the circumference from where the shelf touches instep of the foot to the apex of the heel and record on line I.
2. Measure the foot circumference at the metatarsal heads and record on line J.
3. Measure the length from the 1st metatarsal head to the posterior aspect of the heel, and record on line K.
4. Measure the length from the dorsal aspect of the third metatarsal head to where the shelf touches the instep of the foot and record on line M.

STEP 7: Pictures

Now that you have completed your measurements, take pictures of the leg (front, back, side) that show all of the measurement markings, as well as the overall shape of the leg. This will assist our expert design team to further customize the garment according to your patient's individual needs. Email the pictures with your completed forms to: us_orders@sigvaris.com. Don't forget to record the patient's age, height, and weight, as this also helps with customizing the garment.

STEP 8: Select Additional Options

You are ready to complete your forms and place the order.

Choose right or left leg

Choose Size: If the greatest leg circumference is >60cm, order a Super.

Choose Flat or Wavefoam (recommended)

Choose Compreboot Plus (included) or Custom Medaboot (additional charge)

TAKE DIGITAL PHOTOS (REQUIRED) AFTER MEASUREMENT MARKINGS ARE MADE. EMAIL TO: us_orders@sigvaris.com

NOTE: Legassist Calf & Foot With Shelf is available in:

7099-BK Flat Foam, 7199-BK WaveFoam, 7299-BK Super, Flat Foam, 7399-BK Super, WaveFoam.



Measuring and Ordering **CUSTOM GARMENTS**

Measuring instructions and order forms for custom Legassist Lobule Compression System (LCS) compression wear are found in the following pages, and can also be found online at the link below.

Our Customer Care team is available to assist you through the process of ordering a custom Legassist Lobule Compression System (LCS). From determining if your patient is a good candidate for a custom garment to providing measuring guidance, we are here to help!

To assist in making the custom garment, we require photos of the patient as the measurements are taken, and the measurement marks must be visible in the photos. Please email the photos to **us_orders@sigvaris.com**.

All custom garments in the catalog are patterned and custom made according to the provided measurements at our manufacturing facility in Holland, Michigan.

Contact us:
1.800.322.7744, ext. 1
us_orders@sigvaris.com



For a PDF of the product ordering form, scan the QR code or visit **sigvaris.group/mceforms**


Patient Name:	Contact Name:	Height:	Weight:	Date:
Ship Name & Address:	Account Name:	Account #:	P.O #: :	

Product Information

Product includes one Legassist Lobule Compression System (LCS) and one pair of Cotton Liners.

<input type="checkbox"/> Right Leg	<input type="checkbox"/> Left Leg	Foam	Optional
Size: <input type="checkbox"/> Regular <input type="checkbox"/> Super	Size: <input type="checkbox"/> Regular <input type="checkbox"/> Super	<input type="checkbox"/> Flat	<input type="checkbox"/> Hip Attachment
		<input type="checkbox"/> Wavefoam	

Note: If the greatest circumference measurement is >90cm, order a Super.

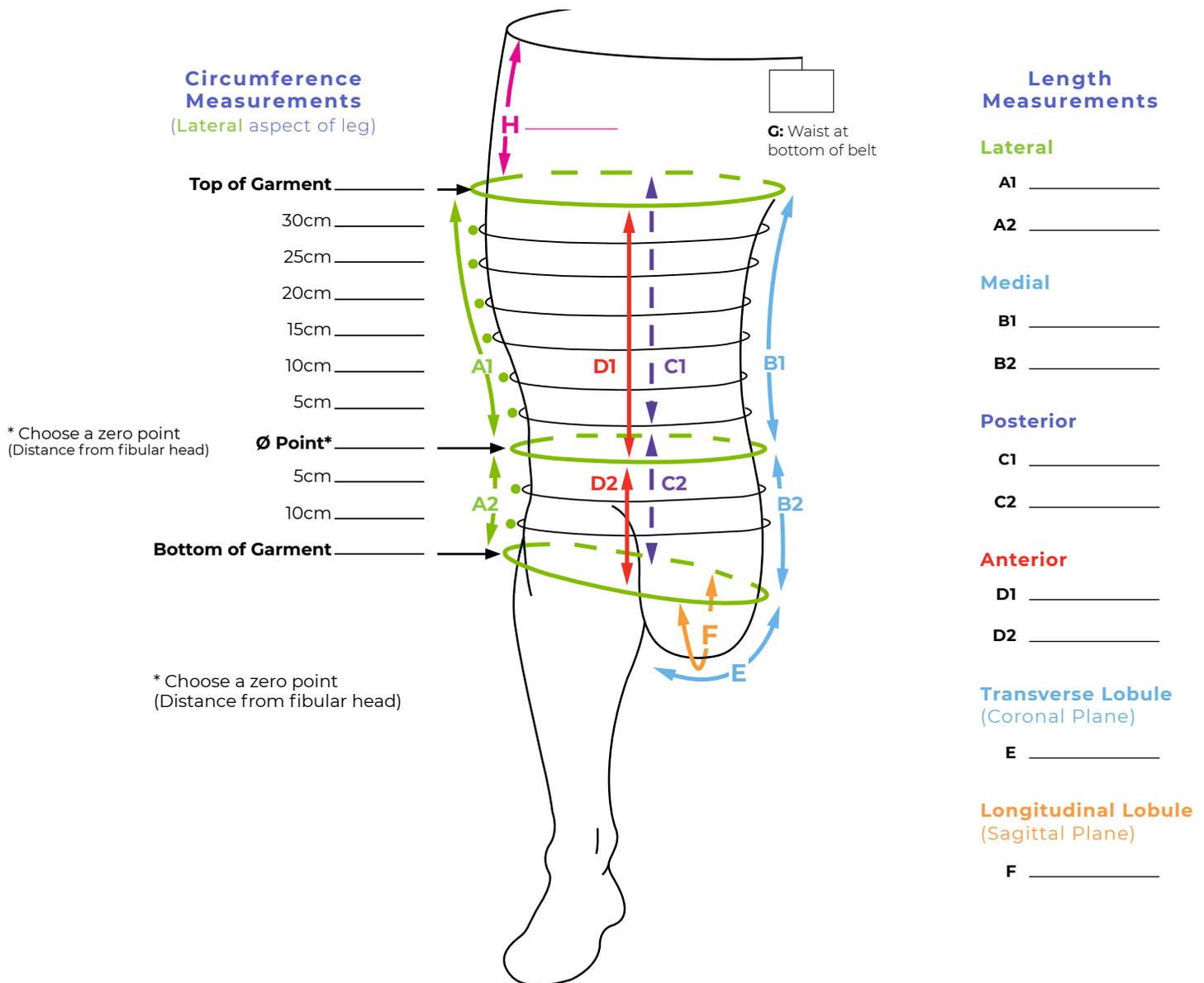

Important

Exact measurements are critical for this garment to ensure a proper fit. If you'd like to learn more about measuring and fitting, attend our Medically Complex Edema Basic Fitter Training. To sign up, visit sigvariseducation.com or contact your local territory manager for more information.

Alternatively, call Customer Care Solution Center at **800-322-7744**, or e-mail us_orders@sigvaris.com, to receive a remote consultation/training.

Supplies Needed

- Cell phone with camera. Photos of the leg with measurement markings must be emailed to: us_orders@sigvaris.com
- Measuring instructions and forms.
- SIGVARIS GROUP Measuring tape and body pen (or eyeliner pencil).
- Signed Custom Order Terms & Conditions Form (include with order).

INSTRUCTIONS

All measurements should be recorded in centimeters. Apply slight tension to hold the tape measure in place.

STEP 1: Choose The Zero Point (Ø point)

The zero point is unique to the Legassist LCS garment. It is the reference point from which all other measurements are taken, in order to determine the lengths needed for the custom garment to fit the patient's individual leg shape. Choose a location proximal to the knee joint where the leg shape is the most cylindrical and the tape measure can be positioned circumferentially around the leg and exactly parallel to the floor. Trace a line around the leg at this Zero Point using the SIGVARIS body pen and record the circumference on Ø point line. If possible, document the location of the zero point by measuring the distance from the head of the fibula to the zero point.

STEP 2: Choose Where The Garment Will End Distally - "bottom of garment"

The length of the Legassist LCS will depend on how low the medial knee lobule hangs down. Start by marking the lowest aspect of the curve of the lobule where a tape measure could be applied without sliding off (lobule apex). Next, choose and mark a point on the lateral aspect of the leg where the garment will overlap the lower leg piece by at least 3 centimeters and create no more than a 30degree downward angle or no more than 10 degree upward angle to the lobule apex. Encircle the leg from this point to the lobule apex, with the tape measure positioned just above your marking points. Trace this circumference with your body pen along the bottom edge of the tape and record the circumference measurement on bottom of garment line.

Problem Solving

If the bottom of garment is proximal to the knee, due to a lobule positioned very high on the thigh, you will need to measure for a customized Legassist TH system that overlaps with the Legassist LCS.

If the bottom of garment extends down to the mid-lower leg or ankle, due to a lobule positioned very low on the leg, you will need to measure for a customized Legassist Calf & Foot that is shortened to best overlap with the Legassist LCS.

STEP 3: Choose Where The Garment Will End Proximally - "top of garment"

You may choose to end the garment circumferentially at the gluteal crease, or extend the garment proximally on the lateral aspect of the thigh to the greater trochanter. Mark your preferred position for the top of the garment along the medial thigh, at least 4cm distal to the groin. Next, mark your preferred position for the top of the garment along the lateral thigh, anywhere from the level of the gluteal fold to the lateral greater trochanter. Encircle the leg with the proximal edge of the tape measure positioned just below your markings, trace a line around the leg along the top edge of the tape with your body pen, and record the circumference on top of garment line.

STEP 4: Mark the Midline of Garment and Record Lengths

Step back to visualize the leg and your marking points to determine where the midline of the garment would be, medial to lateral, encompassing the lobules. Mark the midline of the garment anteriorly (D) and posteriorly (C) at the top of garment, zero point, and bottom of garment lines. Measure each of the lengths from the zero point, following the contour of the leg but not extending the tape measure into folds. Record lengths as follows:

CUSTOM Measuring Instructions**Lateral (A) Lengths:**

A1: Measure & record lateral length from Ø point to top of garment

A2: Measure & record lateral length from Ø point to bottom of garment

Medial (B) Lengths

B1: Measure & record medial length from Ø point to top of garment

B2: Measure & record medial length from Ø point to bottom of garment

Posterior (C) Lengths

C1: Measure & record posterior length from Ø point to top of garment

C2: Measure & record posterior length from Ø point to bottom of garment

Anterior (D) Lengths

D1: Measure & record anterior length from Ø point to top of garment

D2: Measure & record anterior length from Ø point to bottom of garment

STEP 5: Measure And Record Circumferences

Starting from the zero point, on the lateral aspect of the leg only, mark every 5cm proximally from the zero point and mark every 5cm distally from the zero point. Record the circumferences at each point, making sure that the measuring tape is parallel to the floor (Note that some of the spaces to record circumferences may not be used, or you may need to add additional spaces, depending on the length of the limb).

STEP 6: Measure Foot

Measure the size of the lobule for the part of the LCS system that will lift and compress the lobule, called the "Lobule Sling".

Starting at the medial aspect of the lobule, measure the distance from the bottom of garment line to the medial aspect of the leg. Follow the curve of the lobule only to where it touches the leg. Do not extend the tape measure up into the crevices between the lobule and the leg. Record this length as E. At the midline of the lobule, measure from the anterior line of the bottom of garment, underneath the lobule and up to the posterior line of the bottom of garment. Record this length as F.

Hip Attachment (Optional)

1. Measure from lateral aspect of gluteal fold to lateral aspect of hip at the desired location for the bottom of the belt, and record length on line H.
2. Measure circumference of torso at the desired location for the bottom of the belt, and record circumference on line G.

STEP 8: Pictures

Now that you have completed your measurements, take pictures of the leg (front, back, side) that show all of the measurement markings, as well as the overall shape of the leg. This will assist our expert design team to further customize the garment according to your patient's individual needs. Email the pictures with your completed forms to: us_orders@sigvaris.com. Don't forget to record the patient's age, height, and weight, as this also helps with customizing the garment.

STEP 9: Select Additional Options

You are ready to complete your forms and place the order.

Choose right or left leg

Choose Size: If the greatest leg circumference is >90cm, order a super.

Choose Flat or Wavefoam (recommended)

Choose Hip Attachment, if desired (recommended)

Fitting Tips

The ropery, thicker seam should be positioned laterally, and the tag is positioned at the top of the garment straps pull lateral to medial.

TAKE DIGITAL PHOTOS (REQUIRED) AFTER MEASUREMENT MARKINGS ARE MADE. EMAIL TO: us_orders@sigvaris.com



Measuring and Ordering

CUSTOM GARMENTS

Measuring instructions and order forms for custom Legassist Full Leg compression wear are found in the following pages, and can also be found online at the link below.

Our Customer Care team is available to assist you through the process of ordering a custom Legassist Full Leg. From determining if your patient is a good candidate for a custom garment to providing measuring guidance, we are here to help!

To assist in making the custom garment, we require photos of the patient as the measurements are taken, and the measurement marks must be visible in the photos. Please email the photos to us_orders@sigvaris.com.

All custom garments in the catalog are patterned and custom made according to the provided measurements at our manufacturing facility in Holland, Michigan.

Contact us:

1.800.322.7744, ext. 1

us_orders@sigvaris.com



For a PDF of the product ordering form, scan the QR code or visit sigvaris.group/mceforms


Patient Name:	Contact Name:	Height:	Weight:	Date:
Ship Name & Address:	Account Name:	Account #:	P.O #: :	

Product Information

Product includes one Legassist Full Leg, one Compreboot Plus Foot, and one pair of Cotton Liners.

<input type="checkbox"/> Right Leg	<input type="checkbox"/> Left Leg	Foam	Additional Foot Option	Optional
Size: <input type="checkbox"/> Regular <input type="checkbox"/> Super	Size: <input type="checkbox"/> Regular <input type="checkbox"/> Super	<input type="checkbox"/> Flat <input type="checkbox"/> Wavefoam	<input type="checkbox"/> Custom Medaboot (additional charge)	<input type="checkbox"/> Hip Attachment

Note: If the greatest circumference measurement is >90cm, order a Super.

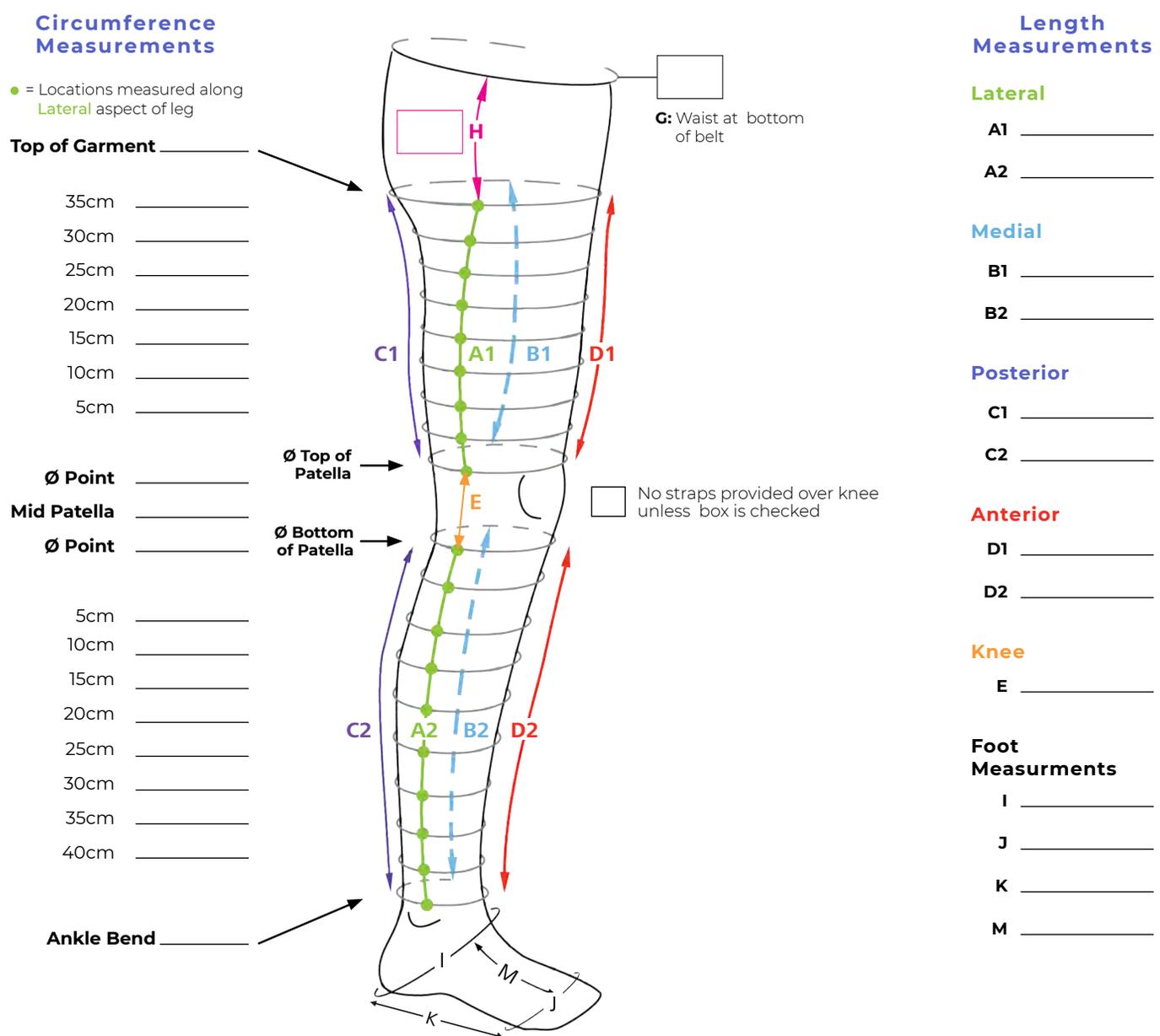

CUSTOM Measuring Instructions

Important

Exact measurements are critical for this garment to ensure a proper fit. If you'd like to learn more about measuring and fitting, attend our Medically Complex Edema Basic Fitter Training. To sign up, visit sigvariseducation.com or contact your local territory manager for more information.

Alternatively, call Customer Care Solution Center at **800-322-7744**, or e-mail us_orders@sigvaris.com, to receive a remote consultation/training.

Supplies Needed

- Cell phone with camera. Photos of the leg with measurement markings must be emailed to: us_orders@sigvaris.com
- Measuring instructions and forms.
- SIGVARIS GROUP Measuring tape and body pen (or eyeliner pencil).
- Signed Custom Order Terms & Conditions Form (include with order).

INSTRUCTIONS

All measurements should be recorded in centimeters. Apply slight tension to hold the tape measure in place.

STEP 1: Locate Landmarks

1. Encircle leg with top edge of tape measure located at the bottom of patella. Mark leg on this edge of tape at the anterior, lateral, posterior and medial aspects, recording the circumference on the line labeled bottom of patella (Ø point).
2. Encircle leg with the bottom edge of the tape measure aligned with the ankle bend. Mark leg on this edge of the tape at the anterior, lateral, posterior and medial aspects, recording the circumference on the line labeled ankle bend.
3. Encircle leg with bottom edge of the tape measure aligned with the top of patella. Mark leg on this edge of the tape at the anterior, lateral, posterior and medial aspects, recording the circumference on the line labeled top of patella.
4. You may choose to end the garment circumferentially at the gluteal crease, or extend the garment proximally on the lateral aspect of the thigh to the greater trochanter. Mark your preferred position for the top of the garment along the medial thigh, at least 4cm distal to the groin. Next, mark your preferred position for the top of the garment along the lateral thigh, anywhere from the level of the gluteal fold to the lateral greater trochanter. Encircle the leg with the proximal edge of the tape measure positioned just below your markings, trace a line around the leg along the top edge of the tape with your body pen, and record the circumference on top of garment line.

STEP 2: Locate Circumference Measurement Positions

1. Starting with "0" of the tape measure at your mark on the lateral aspect of the bottom of patella (Ø point), mark the leg in 5cm increments down to your mark on the lateral aspect of the ankle bend.
2. Repeat step 1, only this time from the top of the patella (Ø point) up to the lateral aspect of the top of garment.

STEP 3: Measure Lengths

Lateral (A1) Lengths:

Measure upper leg from your mark on the lateral aspect of top of patella to your mark on lateral aspect of top of garment. Record length on line A1.

Medial (B1) Lengths

Measure upper leg from your mark on medial aspect of top of patella to your mark on medial aspect of top of garment. Record length on line B1.

Posterior (C1) Lengths

Measure upper leg from your mark on posterior aspect of top of patella to your mark on posterior aspect of top of garment. Record length on line C1.

Anterior (D1) Lengths

Measure upper leg from your mark on anterior aspect of top of patella to your mark on anterior aspect of top of garment. Record length on line D2.

CUSTOM Measuring Instructions**Lateral (A2) Lengths:**

Measure lower leg from your mark on the lateral aspect of bottom of patella to your mark on lateral aspect of ankle bend. Record length on Line A2.

Medial (B2) Lengths

Measure lower leg from your mark on medial aspect of bottom of patella to your mark on medial aspect of ankle bend. Record length on Line B2.

Posterior (C2) Lengths

Measure lower leg from your mark on posterior aspect of bottom of patella to your mark on posterior aspect of ankle bend. Record length on line C2.

Anterior (D2) Lengths

Measure lower leg from your mark on anterior aspect of bottom of patella to your mark on anterior aspect of ankle bend. Record length on line D2.

Knee Space

Measure leg from your mark on the lateral aspect of top of patella to your mark on lateral aspect of bottom of patella. Record length on line E.

STEP 4: Measure And Record Circumferences

Upper: Encircle leg with tape measure at each mark located on lateral aspect between the top of patella and top of garment, and record each circumference on the corresponding line.

Knee: Encircle leg with tape measure at mid patella and record each circumference on the corresponding line.

Lower: Encircle leg with tape measure at each mark located on the lateral aspect between bottom of patella and ankle bend, and record each circumference on the corresponding line.

STEP 5: Measure Foot

1. Measure foot from heel to 1st metatarsal head (or desired boot length) and record length on line K.
2. Measure top of foot from the 3rd metatarsal head to ankle bend and record length on line M.
3. Encircle the instep of foot and apex of heel with a tape measure and record the circumference on line I.
4. Encircle the foot across the metatarsal heads and record the circumference on line J.

Hip Attachment (Optional)

1. Measure from lateral aspect of gluteal fold to lateral aspect of hip at the desired location for the bottom of the belt, and record length on line H.
2. Measure circumference of torso at the desired location for the bottom of the belt, and record circumference on line G.

TAKE DIGITAL PHOTOS (REQUIRED) AFTER MEASUREMENT MARKINGS ARE MADE. EMAIL TO: us_orders@sigvaris.com

NOTES: Legassist Full Leg is available in:

7099-TH Flat Foam, 7199-TH WaveFoam, 7299-TH Super Flat Foam, 7399-TH Super WaveFoam, TH with Hip Attachment is available in 7099-TA Flat Foam, 7199-TA Wavefoam, 7299-TA Super Flat Foam, 7399-TA Super WaveFoam (Includes MedaBoot Custom 7099-MB)



Medaboot Custom **FOOT**

Compression Level

20–50mmHg



Color



Black

Features

The Medaboot Custom Foot is a compression garment that contains and compresses the foot and ankle. Ideal for individuals with a non-standard foot size, or shape, the Medaboot Custom Foot contains added padding at the dorsum to stimulate and compress areas where edema is present.

Ideal for

- Lymphedema

Includes

- Medaboot Custom Foot

Product Selection Tip

May be right for lymphedema patients in need of containment with a non-standard foot size



Custom

Medaboot Custom FOOT

Patient Name:	Contact Name:
Account Name:	Account #:
Account Phone #:	P.O. #:
Height:	Weight:
Ship Name & Address:	Date:

Product Information

Product includes one Medaboot.

- | | | |
|-------------------------------------|--------------|-----------|
| <input type="checkbox"/> Right Foot | Color: Black | Quantity: |
| <input type="checkbox"/> Left Foot | Color: Black | Quantity: |

Important

Exact measurements are critical for this garment to ensure a proper fit. If you'd like to learn more about measuring and fitting, attend our Medically Complex Edema Basic Fitter Training. To sign up, visit sigvariseducation.com or contact your local territory manager for more information.

Alternatively, call Customer Care Solution Center at **800-322-7744**, or e-mail us_orders@sigvaris.com, to receive a remote consultation/training.

Supplies Needed


- Cell phone with camera. Photos of the foot with measurement markings must be emailed to: us_orders@sigvaris.com
- Measuring instructions and forms.
- SIGVARIS GROUP Measuring tape and body pen (or eyeliner pencil).
- Signed Custom Order Terms & Conditions Form (include with order).

Circumference

	Left	Right
I	_____	_____
J	_____	_____
L	_____	_____

Length

	Left	Right
K	_____	_____
M	_____	_____

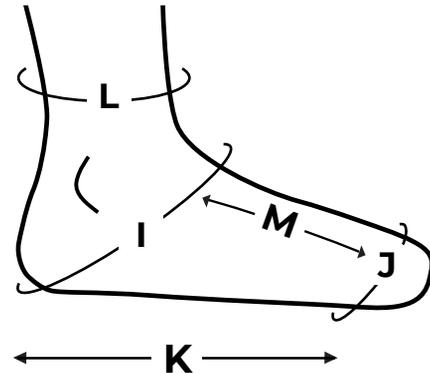

K: Measure **medial** length from heel to 1st metatarsal head

MEASURING INSTRUCTIONS

All measurements should be recorded in centimeters. Apply slight tension to hold the tape measure in place.

MEASURE FOOT

1. Measure foot medially from heel to 1st metatarsal head (or desired boot length) and record length on line K.
2. Measure top of foot from the 3rd metatarsal head to ankle bend and record length on line M.
3. Encircle the ankle bend and heel with a tape measure and record the circumference on line I.
4. Encircle the foot across the metatarsal heads and record the circumference on line J.
5. Measure the ankle circumference and record on line L.







sigvaris

Medical Compression Stockings

Sigvaris Medical Compression Stockings offer wearers four categories: Essential, Style, Motion, and Specialty. Each category includes unique styles, patterns, and compression levels that will fit any lifestyle.

Below, you will find a small sample of Sigvaris medical compression stockings. To see all of the product series, check out the latest Sigvaris medical compression stockings catalog.

sigvaris

Essential

COTTON 230

Essential Cotton combines a classic ribbed design and all-day cotton comfort that is ideal for sensitive skin.

Compression Level

20–30mmHg

30–40mmHg

Features and Benefits

- Supima® cotton for highest level of comfort and breathability
- Ideal for sensitive skin and allergies
- Double-covered inlay yarns provide comfort and allow product to easily glide on and off

Fiber Content

- 50% Nylon, 25% Spandex, 25% Supima Cotton (content may vary slightly by style)
- Non-latex

Best Suited for

- Containment of various types of edema including venous edema, phlebolympheidema, and early stages of lymphedema
- Ideal for dress or casual shoes
- Ideal for long-distance travel



White


Navy

Light Beige

Black



232NSLW66: 20–30mmHg, Closed Toe
Thigh Small Long in Light Beige

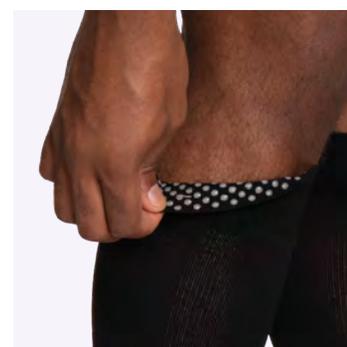

Silicone dotted grip-top available
for calf and thigh

Sizing Chart

	Small		Medium		Large		X Large	
Ankle Circumference	7"-8.5" (18-21.5cm)		8.5"-10" (21.5-25.5cm)		10"-11.5" (25.5-29cm)		11.5"-14" (29-36cm)	
Calf Circumference	11"-15.5" (28-39cm)		12.5"-17.5" (31.5-45cm)		14"-20" (35.5-51cm)		16"-24" (41-61cm)	
Calf Length	Short	Long	Short	Long	Short	Long	Short	Long
	≤16" (41cm)	>16" (41cm)	≤16" (41cm)	>16" (41cm)	≤16" (41cm)	>16" (41cm)	≤16" (41cm)	>16" (41cm)
Thigh Circumference	18"-26" (45.5-66cm)		19"-28.5" (48-73cm)		20.5"-33" (52-84cm)		NA	
Leg Length: Women & Open Toe	Short	Long	Short	Long	Short	Long		
	≤28.5" (72cm)	>28.5" (72cm)	≤28.5" (72cm)	>28.5" (72cm)	≤28.5" (72cm)	>28.5" (72cm)		
Leg Length: Men	Short	Long	Short	Long	Short	Long		
	≤30" (76cm)	>30" (76cm)	≤30" (76cm)	>30" (76cm)	≤30" (76cm)	>30" (76cm)		
Size	SS	SL	MS	ML	LS	LL	XS	XL

Women's Shoe Sizes: Short Leg 5.5-9.5, Long Leg 8-12 **Men's Shoe Sizes:** Short Leg 7-11, Long Leg 9.5-14

Styles

Colors

Closed Toe = ●
Open Toe = ○

White
00

Navy
10

**Light
Beige**
66

Black
99

	White 00	Navy 10	Light Beige 66	Black 99
Calf (C)	●	●	●○	●
Calf w/Grip-Top (C/S) Silicone band grip-top			●*	●**
Thigh w/Grip-Top (N) Silicone band grip-top			●○	●

* 20-30mmHg only

** 30-40mmHg calf with grip-top only available in Black, 30-40mmHg extra large calf only available in Black with grip-top

Product Code Example: 232CMLM99

Series	mmHg Level	Style	Size	Gender/ Open Toe	Color
23	2	C	ML	M	99

Product Code Example: 232CLSW66/S

Series	mmHg Level	Style	Size	Gender/ Open Toe	Color	Grip-top
23	2	C	LS	W	66	/S

us_orders@sigvaris.com

P: 800-322-7744 • F: 800-481-5488

Essential

OPAQUE 860

Essential Opaque offers durability you can trust in our most extensive size range, making it your go-to choice for everyday hosiery. *(For Men's Essential Opaque, see page 130)*

Compression Level

20–30mmHg

30–40mmHg

Features and Benefits

- Double-covered inlay yarns provide comfort and allow product to easily glide on and off
- Cotton-lined gusset for pantyhose styles
- The Pantyhose plus style is a great option during pregnancy

Fiber Content

- 64% Nylon, 36% Spandex
(content may vary slightly by style)
- Non-latex

Best Suited for

- Containment of various types of edema including venous edema, phlebolympheidema, and early stages of lymphedema
- Post-vein treatment including thermal and non-thermal techniques
- Patients looking for petite-sized and plus-sized options



Honey

Golden

Light Beige

Black



862CLLW33: 20–30mmHg, Closed Toe
Calf Large Long in Honey



Silicone dotted grip-top available
for calf and thigh

Women's Sizing Chart (for Men's Sizing Chart, see page 131)

	Small			Medium			Large		X Large	
Ankle Circumference	7"–8.5" (18–21.5cm)			8.5"–10" (21.5–25.5cm)			10"–11.5" (25.5–29cm)		11.5"–14" (29–36cm)	
Calf Circumference	11"–15.5" (28–39cm)			12.5"–17.5" (31.5–45cm)			14"–20" (35.5–51cm)		16"–24" (41–61cm)	
Calf Circumference: Calf Plus w/Grip-Top	NA			15.5"–21" (39.5–53cm)			18"–24" (45.5–61cm)		20"–28" (51–71cm)	
Calf Length	Petite	Short	Long	Petite	Short	Long	Short	Long	Short	Long
	≤14" (35.5cm)	≤16" (41cm)	>16" (41cm)	≤14" (35.5cm)	≤16" (41cm)	>16" (41cm)	≤16" (41cm)	>16" (41cm)	≤16" (41cm)	>16" (41cm)
Thigh Circumference	18"–26" (45.5–66cm)			19"–28.5" (48–73cm)			20.5"–33" (52–84cm)		25"–37.5" (64–95cm)	
Leg Length: Women & Open Toe	Petite	Short	Long	Petite	Short	Long	Short	Long	Short	Long
	≤25.5" (65cm)	≤28.5" (72cm)	>28.5" (72cm)	≤25.5" (65cm)	≤28.5" (72cm)	>28.5" (72cm)	≤28.5" (72cm)	>28.5" (72cm)	≤28.5" (72cm)	>28.5" (72cm)
Hip Circumference Guidelines	32.5"–44.5" (83–112.5cm)			34.5"–50" (87–127cm)			38.5"–59.5" (98–150.5cm)		NA	
Size	SP	SS	SL	MP	MS	ML	LS	LL	XS	XL
Calf Plus w/Grip-Top				MS/P	ML/P	LS/P	LL/P	XS/P	XL/P	

Women's Shoe Sizes: Petite Leg 4–8, Short Leg 5.5–9.5, Long Leg 8–12

Plus (Maternity) Hip Circumference Guidelines Small ≤56" (142cm), Medium ≤62" (157.5cm), Large ≤70" (178cm)

Women's Styles (for Men's, see page 131)

Colors

Closed Toe = ●
Open Toe = ○

Honey 33 **Golden** 36 **Light Beige** 66 **Black** 99

	Honey 33	Golden 36	Light Beige 66	Black 99
Women's Calf (C)	●	●	○	●●
Women's Calf Petite Sizes (C) Only available in 20–30mmHg, not available with silicone band grip-top	●	●		●
Women's Calf w/Grip-Top (C/S) Silicone band grip-top			●●	●●
Women's Calf Plus w/Grip-Top (C/P) Only available with silicone band grip-top			●●	●●
Women's Thigh w/Grip-Top (N) Silicone band grip-top*	●	●	○	●●
Women's Thigh w/Grip-Top Petite (N) Only available in 20–30mmHg	●	●		●
Women's Pantyhose (P) Not available in XS or XL sizes	●	●		●
Women's Pantyhose Petite (P) Only available in 20–30mmHg	●	●		●
Women's Pantyhose Plus (P/P)			●	●
Thigh w/Waist Attachment (W)** Indicate leg LEFT (/L) or RIGHT (/R)			○	

* Women's Thigh w/Grip-Top (N) and Open Toe Thigh w/Grip-Top (N) sizes XS and XL are only available in Light Beige and Black

** Thigh w/Waist Attachment not available in XL sizes

Product Code Example: 862NSSO99

Product Code Example: 862CLLW66/P

Series	mmHg Level	Style	Size	Gender/ Open Toe	Color
86	2	N	SS	O	99

Series	mmHg Level	Style	Size	Gender/ Open Toe	Color	Grip-top
86	2	C	LL	W	66	/P

us_orders@sigvaris.com

P: 800-322-7744 • F: 800-481-5488

Essential

OPAQUE 860

Essential Opaque offers durability you can trust in our most extensive size range, making it your go-to choice for everyday hosiery. *(For Women's Essential Opaque, see page 128)*

Compression Level

20–30mmHg

30–40mmHg

Features and Benefits

- Double-covered inlay yarns provide comfort and allow product to easily glide on and off
- Men's calf and thigh products feature ribbed fabric; open toe products are a smooth hosiery fabric
- Flat, low friction toe seam

Fiber Content

- 64% Nylon, 36% Spandex
(content may vary slightly by style)
- Non-latex

Best Suited for

- Containment of various types of edema including venous edema, phlebolymphe-
dema, and early stages of lymphedema
- Post-vein treatment including thermal and non-thermal techniques
- Patients looking for plus-sized options



**Ribbed
Light Beige**



**Ribbed
Black**



**Smooth
Light Beige**



**Smooth
Black**



862CML099: 20–30mmHg
Open Toe Calf Medium Long in Black



Silicone dotted grip-top available
for calf and thigh

Men's Sizing Chart (for Women's Sizing Chart, see page 129)

	Small		Medium		Large		X Large	
Ankle Circumference	7"-8.5" (18-21.5cm)		8.5"-10" (21.5-25.5cm)		10"-11.5" (25.5-29cm)		11.5"-14" (29-36cm)	
Calf Circumference	11"-15.5" (28-39cm)		12.5"-17.5" (31.5-45cm)		14"-20" (35.5-51cm)		16"-24" (41-61cm)	
Calf Circumference: Calf Plus w/Grip-Top	NA		15.5"-21" (39.5-53cm)		18"-24" (45.5-61cm)		20"-28" (51-71cm)	
Calf Length	Short	Long	Short	Long	Short	Long	Short	Long
	≤16" (41cm)	>16" (41cm)	≤16" (41cm)	>16" (41cm)	≤16" (41cm)	>16" (41cm)	≤16" (41cm)	>16" (41cm)
Thigh Circumference	18"-26" (45.5-66cm)		19"-28.5" (48-73cm)		20.5"-33" (52-84cm)		25"-37.5" (64-95cm)	
Leg Length: Men	Short	Long	Short	Long	Short	Long	Short	Long
	≤30" (76cm)	>30" (76cm)	≤30" (76cm)	>30" (76cm)	≤30" (76cm)	>30" (76cm)	≤30" (76cm)	>30" (76cm)
Leg Length: Open Toe	Short	Long	Short	Long	Short	Long	Short	Long
	≤28.5" (72cm)	>28.5" (72cm)	≤28.5" (72cm)	>28.5" (72cm)	≤28.5" (72cm)	>28.5" (72cm)	≤28.5" (72cm)	>28.5" (72cm)
Size	SS	SL	MS	ML	LS	LL	XS	XL
Calf Plus w/Grip-Top			MS/P	ML/P	LS/P	LL/P	XS/P	XL/P

Men's Shoe Sizes: Short Leg 7-11, Long Leg 9.5-14

Men's Styles (for Women's Styles and Colors, see page 129)

Colors

Closed Toe = ● closed toe men's products are ribbed hosiery
 Open Toe = ○ open toe men's products are smooth hosiery

Light
Beige
66

Black
99

	Light Beige 66	Black 99
Men's Calf (C)	●●	●●
Men's Calf w/Grip-Top (C/S) Silicone band grip-top	●●	●●
Men's Calf Plus w/Grip-Top (C/P) Only available with silicone band grip-top	●●	●●
Men's Thigh w/Grip-Top (N) Silicone band grip-top	●●	●●
Thigh w/Waist Attachment (W)* Indicate leg LEFT (/L) or RIGHT (/R)	○	

* Thigh w/Waist Attachment not available in XL sizes

Product Code Example: 863CMSM99

Series	mmHg Level	Style	Size	Gender/ Open Toe	Color
86	3	C	MS	M	99

Product Code Example: 862CLLO66/P

Series	mmHg Level	Style	Size	Gender/ Open Toe	Color	Grip-top
86	2	C	LL	O	66	/P

us_orders@sigvaris.com

P: 800-322-7744 • F: 800-481-5488

sigvaris

Specialty

SECURE 550

Secure combines compression and containment with comfort and durability that's optimal for comfortably managing advanced venous edema, lymphedema, and wound care.

Compression Level

20–30mmHg

30–40mmHg

40–50mmHg*

Features and Benefits

- More graduated compression and containment than other traditional garments
- Durable double-covered inlay yarns allow comfortable all-day wear

Fiber Content

- 78% Nylon, 22% Spandex (Calf style)
- 76% Nylon, 24% Spandex (Thigh style)
- Non-latex

Best Suited for

- Chronic Venous Insufficiency (CVI)¹
- Lymphedema that is well-managed
- Lipedema
- Varicosities
- Deep vein thrombosis (DVT)
- Superficial phlebitis
- Post-thrombotic syndrome
- Orthostatic hypotension
- People with sensitivity to latex



Beige



Black



554NM4W77: 40–50mmHg, Closed Toe Thigh Medium Full Long in Beige



Silicone dotted grip-top on calf and thigh

Sizing Chart

	Medium				Large				X Large				XX Large			
Ankle Circumference	8.5"-10" (21.5-25.5cm)				10"-11.5" (25.5-29cm)				11.5"-13" (29-33cm)				13"-14.5" (33-37cm)			
Calf Circumference	Average		Full		Average		Full		Average		Full		Average		Full	
	12"-17.5" (31-45cm)		17.5"-23" (45-59cm)		13.5"-19" (35-49cm)		19"-25" (49-63cm)		16"-22" (41-56cm)		22"-28" (56-71cm)		17.5"-23" (45-59cm)		23"-29" (59-74cm)	
Calf Length	Short	Long	Short	Long	Short	Long	Short	Long	Short	Long	Short	Long	Short	Long	Short	Long
	≤16" (41cm)	>16" (41cm)	≤16" (41cm)	>16" (41cm)	≤16" (41cm)	>16" (41cm)	≤16" (41cm)	>16" (41cm)	≤16" (41cm)	>16" (41cm)	≤16" (41cm)	>16" (41cm)	≤16" (41cm)	>16" (41cm)	≤16" (41cm)	>16" (41cm)
Thigh Circumference	Average		Full		Average		Full		Average		Full		Average		Full	
	20"-27" (51-69cm)		23.5"-31" (60-78cm)		24"-30" (61-77cm)		26.5"-34" (68-86cm)		27"-33.5" (69-85cm)		30"-37" (76-94cm)		29.5"-36.5" (75-93cm)		33"-40" (84-102cm)	
Leg Length	Short	Long	Short	Long	Short	Long	Short	Long	Short	Long	Short	Long	Short	Long	Short	Long
	≤28.5" (72cm)	>28.5" (72cm)	≤28.5" (72cm)	>28.5" (72cm)	≤28.5" (72cm)	>28.5" (72cm)	≤28.5" (72cm)	>28.5" (72cm)	≤28.5" (72cm)	>28.5" (72cm)	≤28.5" (72cm)	>28.5" (72cm)	≤28.5" (72cm)	>28.5" (72cm)	≤28.5" (72cm)	>28.5" (72cm)
Size	M1	M2	M3	M4	L1	L2	L3	L4	X1	X2	X3	X4	G1	G2	G3	G4

Women's Shoe Sizes: Short Leg 5.5-9.5, Long Leg 8-12 Men's Shoe Sizes: Short Leg 7-11, Long Leg 9.5-14

Styles

Colors

Closed Toe = ●
Open Toe = ○

Beige
77

Black
99

	Beige 77	Black 99
Calf w/Grip-Top (C/S) Silicone band grip-top	●●	●●
Thigh w/Grip-Top (N) Silicone band grip-top	●●	●●

Product Code Example: 552CM2W77/S

Series	mmHg Level	Style	Size	Gender/ Open Toe	Color	Grip-top
55	2	C	M2	W	77	/S

us_orders@sigvaris.com

P: 800-322-7744 • F: 800-481-5488

sigvaris

Specialty

DIABETIC COMPRESSION SOCKS 602

Diabetic Compression Socks are the only patented, clinically-proven graduated compression socks constructed for people with diabetes.

Compression Level

18–25mmHg

Also available in 8–15mmHg, see Well Being Lookbook

Features and Benefits

- Padded soles and invisible toe seam add comfort and protection
- Compression-free foot for optimum comfort
- Double-covered inlay yarns provide comfort and allow product to easily glide on and off
- *drirelease*® yarns add moisture wicking properties

Fiber Content

- 76% *drirelease* (65% Polyester, 11% Cotton), 19% Nylon, 5% Spandex
- Non-latex

Best Suited for

- Individuals with diabetes and symptoms related to Chronic Venous Disease (CVD)¹
- Individuals with peripheral neuropathy
- Patients with an ABI of >0.6



White



602CSSW00: 18–25mmHg, Closed Toe Calf Small Short in White



Padded soles and invisible toe seam add comfort and protection

Specialty

DIABETIC COMPRESSION SOCKS 602

Sizing Chart

	Small		Medium		Large		X Large	
Ankle Circumference	7"-8.5" (18-21.5cm)		8.5"-10" (21.5-25.5cm)		10"-11.5" (25.5-29cm)		11.5"-14" (29-36cm)	
Calf Circumference	11"-15.5" (28-39cm)		12.5"-17.5" (31.5-45cm)		14"-20" (35.5-51cm)		16"-24" (41-61cm)	
Calf Length	Short	Long	Short	Long	Short	Long	Short	Long
	≤16" (41cm)	>16" (41cm)	≤16" (41cm)	>16" (41cm)	≤16" (41cm)	>16" (41cm)	≤16" (41cm)	>16" (41cm)
Size	SS	SL	MS	ML	LS	LL	XS*	XL*

Women's Shoe Sizes: Short Leg 5.5-9.5, Long Leg 8-12 **Men's Shoe Sizes:** Short Leg 7-11, Long Leg 9.5-14

*XS & XL Only Available In Men's Sizing

Style

Color

Closed Toe = ●

White
00

Calf (C)	●
-----------------	---

Product Code Example: 602CMLM00

Series	mmHg Level	Style	Size	Gender/ Open Toe	Color
60	2	C	ML	M	00

© 2023 Copyright SIGVARIS, Inc.
Patent No. US 9,566,206 B2

us_orders@sigvaris.com

P: 800-322-7744 • F: 800-481-5488





sigvaris

How to Measure and Order

COMPRESSION GARMENTS

Ordering a Sigvaris compression garment for standard or non-standard patients is easy. To order, simply call our Customer Care agents. We are available Monday–Friday from 8am–7pm EST.

Online: sigvaris-online.com

Call: 1.800.322.7744

Fax: 1.800.481.5488

Sigvaris follows ISO Standard 13485:2016 and 9001:2015 for our products manufactured in our North American facilities.

Measuring and Ordering **COMPRESSION GARMENTS**

Precise measurements are essential to guarantee the indicated level of compression and to ensure a proper fit. We recommend that measurements be taken by a healthcare professional.

Detailed instructions for measuring for standard (non-custom) garments are listed below. The measurement chart and graphic for each product is also included on the corresponding order form. Follow the chart to identify the correct positioning and type of measurement needed.

Instructions for measuring for non-standard (custom) garments are included on each ordering form. If you need additional assistance or support with measuring for both standard and non-standard garments, our sales and customer care teams are available to assist you.

For guidance with measuring or fitting a patient for a garment, please contact Customer Care at: **1.800.322.7744**, ext. 1.

STANDARD Measuring Instructions

ARM

With patient seated, place the arm extended and **Elbow** bent slightly, with the palm down, on a flat surface.

Circumferences:

Measure circumference at **Palm** and record on line E (Medaform arm only).

Measure at **Wrist**, and record on line A. Mark **dorsal** aspect at distal edge of tape.

Measure at **Elbow**, and record on line B.

Measure at **Mid Bicep**, and record on line B1.

Measure at **Axilla**, and record on line C. Mark **dorsal** aspect at proximal edge of tape.

Length:

Measure length at **dorsal** aspect from mark at Point A to mark at Point C and record in box D.

CALF & FOOT

Circumferences:

Measure circumferences where indicated by the black dots, beginning with the dot 5cm above the **Lateral Malleolus** (outer ankle bone), continue to the 15cm, 25cm and 30cm dots. Record in the corresponding lines A, B, C, C1(left/right).

Length:

Measure length from the **Lateral Malleolus** (outer ankle bone) to lateral aspect of the **Fibular Head**, and record in box G.

Foot:

Measure circumference of ankle bend and heel and record on line I. Measure circumference across metatarsal heads and record on line J. Measure length on medial side of foot from heel to 1st metatarsal head and record on line K.

How to Measure and Order

COMPRESSION GARMENTS

THIGH COMPONENT

Circumferences:

Measure circumferences where indicated by the black dots, beginning with the dot 5cm above the lateral aspect **Top of Patella** and continue to the 15cm, and 23cm dots. Record in the corresponding lines D, E, F (left/right).

Length: (Compreflex only)

Measure length from the lateral aspect **Top of Patella** to lateral aspect **Gluteal Fold**, and record in box G1.

THIGH HIGH

Circumferences:

Measure circumferences where indicated by the black dots, beginning with the dot 5cm above the **Lateral Malleolus** (outer ankle bone) and continue to the 15cm, 25cm and 30cm dots. Record in the corresponding lines A, B, C, C1(Left/right). Measure circumferences at the 5cm dot from the lateral aspect top of patella and continue to the 15cm, and 23cm dots. Record in the corresponding lines D, E, F (left/right).

Length: (Compreflex only)

Measure length from the **Lateral Malleolus** (outer ankle bone) to lateral aspect of the **Gluteal Fold**, and record in box H.

KNEE

Circumferences:

Measure circumference where indicated on the measuring chart. Begin 30cm above the **Lateral Malleolus** (outer ankle bone), then continue to 5cm above the **Fibular Head**. Record the 30cm and 5cm measurements on the corresponding lines C1 and D.

Length: (Compreflex only)

Measure length from the **Lateral Malleolus** (outer ankle bone) to lateral aspect of the **Gluteal Fold**, and record in box H.

FOOT

Circumferences:

Measure circumference of ankle bend and heel and record on line I. Measure circumference across metatarsal heads and record on line J.

Length:

Measure length on medial side of foot from heel to 1st metatarsal head and record on line K.

For more information, please speak with your SIGVARIS GROUP Sales Representative or Customer Care Solution Center: 1.800.322.7744.





sigvaris

Accessories

Help yourself

Why choose **Accessories**?

Created for specific medical needs, including:

- Aid in putting on and removing stockings
- Stocking cleaning to maximize usage
- Adhesives and extra liners for added comfort

Gloves

Rubber Gloves – Ridged
The practical and protective way to apply and remove your compression stockings.



Package: 12 per Pack

Size	Product Code
Small	591R400S
Medium	591R400M
Large	591R400L
X-Large	591R400X

Package: Single Pair

Size	Product Code
Small	591RPRS
Medium	591RPRM
Large	591RPRL
X-Large	591RPRX

Latex-Free Gloves
Gripping donning gloves recommended for use in a non-latex environment.



Package: 12 per Pack

Size	Product Code
Small	592R400S
Medium	592R400M
Large	592R400L
X-Large	592R400X

Package: Single Pair

Size	Product Code
Small	592RPRS
Medium	592RPRM
Large	592RPRL
X-Large	592RPRX

Glove Sizing Chart

Measure Circumference of the Hand at the Palm

6.5"–7" (16.5–18cm)	7.5"–8" (19–20.5cm)	8.5"–9" (21.5–23cm)	9.5"–10" (24–25.5cm)
Small	Medium	Large	X-Large



For more Sigvaris videos, including donning tips and product demonstrations, visit our YouTube channel: [youtube.com/sigvarisgroupusa](https://www.youtube.com/sigvarisgroupusa)

Caution: Rubber Gloves contain natural rubber latex which may cause allergic reaction

Donning Aids



Sim-Slide



Video: How to Use the Sim-Slide
To watch the video go to: bit.ly/SIGssyt

Sim-Slide

Size	Product Code
Small	587SLIDESM
Medium	587SLIDEMD
Large	587SLIDELG
X-Large	587SLIDEXL

Sizing Chart

For use with Open toe styles

Women's Shoe Size	≤4–6	6.5–8	8.5–9.5	≥10
Men's Shoe Size	≤3–4.5	5–6.5	7–8	≥8.5
Size	S	M	L	XL



Magnide On/Off



Video: How to Use the Magnide On/Off
To watch the video go to: bit.ly/SIGmofyt

Magnide On/Off

Size	Product Code
Medium	587MAGMD
Large	587MAGLG
X-Large	587MAGXL

Sizing Chart

For use with Closed toe styles

Women's Shoe Size	≤8.5	9–12	≥12.5
Men's Shoe Size	≤7.5	8–11	≥11.5
Size	M	L	XL*

*XL suitable for extra-wide feet

Accessories

Donning Aid



Cotton Underliner

10–15mmHg stocking. Can be used as an underlayer to protect fragile skin, secure wound dressings, to add 10–15mmHg of compression, or create an easy glide surface to aid in putting on and taking off compression stockings. For sizing, please refer to the Secure Calf sizing on page 139.

Package: 2 per pack

Size	Product Code
Small / Short	588USS
Small / Long	588USL
Medium / Short	588UMS
Medium / Long	588UML
Large-XL / Short	588ULS
Large-XL / Long	588ULL

Washing Solution



Washing Solution

Gentle washing solution specifically formulated to help maintain the quality of SIGVARIS GROUP compression socks and stockings. Also great for washing delicate garments!

Package Options

Contains	Product Code
Six 16 oz. Bottles (473ml)	586W900
One 16 oz. Bottle (473ml)	586W950
Twelve 2 oz. Bottles (60ml)	586W975
One 2 oz. Bottle (60ml)	586W985

Body Adhesive



IT STAYS!

A gentle roll-on body adhesive that washes off with water and may be used to hold compression garments in place against the skin. IT STAYS! is pliable and moves with the body.

Package Options

Contains	Product Code
Twelve 2 oz. Bottles (60ml)	581I300
One 2 oz. Bottle (60ml)	581I350



Donning Aids



Foot Slip

Comfortable applicator for use with open toe stockings.

Package: 6 per pack

Size	Product Code
One Size	580F200



Sigvaris DOFF N' DONNER

A soft device that allows you to quickly guide a compression sock or stocking up and down your leg or arm, making it simple to put on and remove all styles of graduated compression. It is perfect for both self-donning and caregiver use. Non-latex.

Size	Product Code
One Size	583D100



Sigvaris DOFF N' DONNER CONE

This loading device makes it easy to roll stockings onto the Sigvaris DOFF N' DONNER. The bottom of the base firmly attaches to most clean flat surfaces via suction.



Video: How to Use the Sigvaris DOFF N' DONNER
To watch the video go to: bit.ly/sigdnd

Size	Product Code
One Size	583D200

Caution: this product contains natural rubber latex which may cause allergic reactions



MELANY

The donning and doffing device for compression stockings in case of limited mobility.

Sizing Chart				
Size	X-SMALL	SMALL	MEDIUM	LARGE
Calf Circumference	10"-12" (25-30cm)	12"-15" (30-38cm)	15"-18" (38-46cm)	18"-22" (46-56cm)
Product Code	579MXS	579MSM	579MMD	579MLG

Measure at the widest part of the calf



Video: How to Use the MELANY
To watch the video go to: goto.sigvaris.com/bbiMdz

WARNING: This product can expose you to chemicals, including Carbon black, which are known to the State of California to cause cancer. For more information go to www.P65Warnings.ca.gov.

Dealer Portal & Education

SIGVARIS GROUP Dealer Portal

The SIGVARIS GROUP Dealer Portal is your online location for securely ordering SIGVARIS GROUP products and managing your account. Place and track your orders, review order history, chat with Customer Care, and much more at sigvaris-online.com.


Dealer Portal benefits

- **24/7** access-order any time!
- Sizing Wizard
- Three ways to securely place an order online: Catalog, Quick Entry or Excel Upload
- Access to order history and tracking
- Shipping options including drop ship and easy access to an address table for quick shipment of reorders
- Check your account balance and print invoices
- Credit card payments
- New product launch information
- Current promotions
- Inelastic compression order forms

Email sigvariscares@sigvaris.com and you will receive your active account within one business day!

Can't find the product you're searching for?

Type the product SKU in the search function and select for ordering.

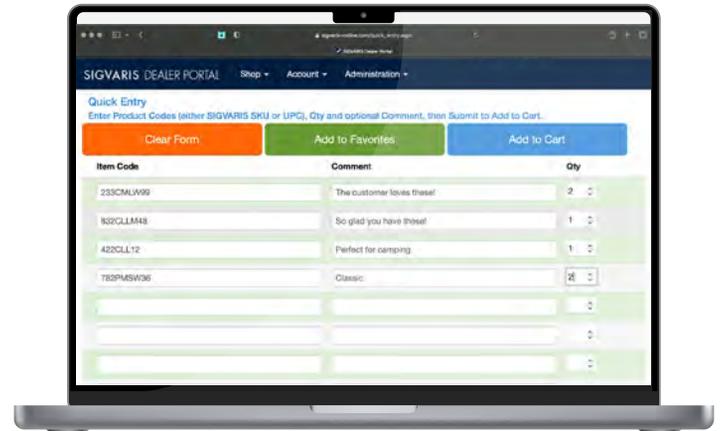

SIGVARIS GROUP Education & Training



SIGVARIS GROUP has several new and exciting educational opportunities for 2023 to help you get a leg up on leg health knowledge. These opportunities include online learning, webinars, and custom training sessions to educate fitters, clinicians, and therapists. Additionally, we offer a variety of courses with opportunities to earn continuing education credits.

Learn more at sigvariseducation.com.

Become a SIGVARIS Certified Fitter!

Contact your Territory Manager or email medical.education@sigvaris.com

SIGVARIS GROUP

SIGVARIS, INC.
1119 Highway 74 South
Peachtree City, GA 30269
Tel 800-322-7744
sigvaris.com

**SIGVARIS
GROUP**

Sigvaris, Inc.
1119 Highway 74 South
Peachtree City, GA 30269

ORDERING

Online: **sigvaris-online.com**

Call: **1-800-322-7744, ext. 1**

Fax: **1-800-481-5488**

Customer Care:
Monday–Friday, 8am–7pm (EST)

QUALITY STANDARD

Holland, MI:
Quality System ISO 9001
FDA Registered
Products are inspected based
on specifications and guidelines.

Peachtree City, GA:
Certified Quality System ISO 13485
FDA Registered

Please consult your doctor, or other qualified health care professional before using any product discussed within the literature. See product insert for complete information.

REFERENCES

1. Partsch, Hugo, et al. "Measurement of lower leg compression in vivo: recommendations for the performance of measurements of interface pressure and stiffness." *Dermatologic Surgery* 32.2 (2006): 224–233.
2. <http://www.tagungsmanagement.org/comp/index.php/icc-consensus>
3. Mosti, Giovanni, and Susanna Rossari. "The importance of measuring sub bandage pressure and presentation of new measuring device." *Acta Vulnol* 6 (2008): 31–36.
4. Mosti, Giovanni, Maria Letizia Iabichella, and Hugo Partsch. "Compression therapy in mixed ulcers increases venous output and arterial perfusion." *Journal of vascular surgery* 55.1 (2012): 122–128.
5. Mosti, Giovanni. "Compression in mixed ulcers: venous side." *Phlebology* 29.1_suppl (2014): 13–17.

Sigvaris, SIGVARIS GROUP, Accutab, Armassist, Chipvest, Compreflex, Comprefit, Compreboot, Compreshort, Comprevest, Doff N' Donner, Legassist, Medahand, Chipsleeve, Medaglove, Medaboot, and Medaform are registered trademarks of SIGVARIS AG, CH-9014 St.Gallen/Switzerland, in many countries worldwide.

Coolflex, Compreknee, Dorsal Pocket, Tension Tab, Genifit, EZ-Close, 90° Flexi-Zone, and Chipflow are trademarks of SIGVARIS AG, CH-9014 St.Gallen/Switzerland, in many countries worldwide.

direrelease is an Optimer Brands company.
Durapore is a registered trademark of 3M.
Supima is a registered trademark of the Supima Association of America Corporation.
PicoPress is a trademark of Microlab Elettronica.

ORDERING

Online
sigvaris-online.com
Call
1-800-322-7744, ext. 1
Fax
1-800-481-5488

CUSTOMER CARE

Monday–Friday
8am–7pm (EST)

Visit us on the web:

sigvaris.com



SIGVARIS GROUP USA